

MANUAL


Multifunctional negative pressure massage instrument



LW-123

Thank you for purchasing and using this product
If the instrument needs to be transferred to others, please transfer it together with the manual

Control panel function introduction

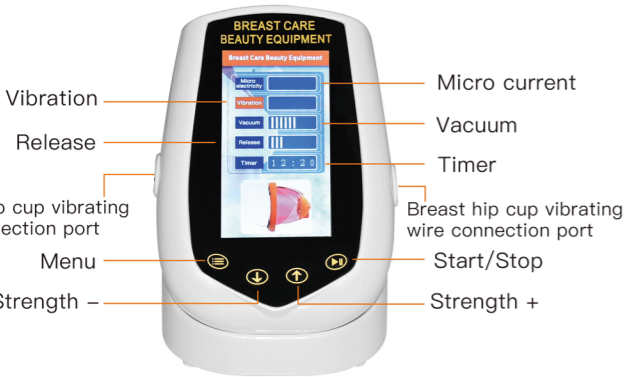


Precautions

1. Please read the manual carefully before using the instrument.
2. Please consult a doctor if you have the following conditions:
 - (1) The body has varicose veins and should not be used
 - (2) Not suitable for use during pregnancy or menstruation
 - (3) Epidermal damage and infection should not be used
 - (4) Consult a doctor for hypoglycemia and high blood pressure
 - (5) After a meal, the blood is concentrated in the digestive system and should not be used
 - (6) Gua Sha operation: Brush in one direction when using. Do not go up and down, left and right.
 - (7) It is recommended to start with a low degree of negative pressure to avoid discomfort caused by excessive suction.

Note: Please keep it out of reach of children

Wide range of applications

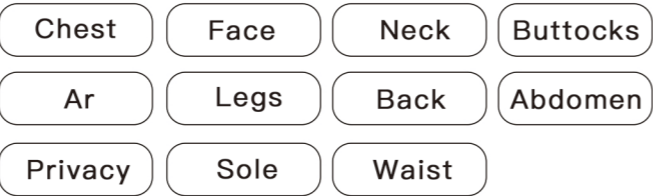

Product Principle

This product is a collection of breast care massage, scraping and cupping, and acupoint massage. It can scrape Gua Sha, remove acid and dampness, relax nerves, mobilize meridians, adjust functions, dispel cold and dehumidify, relieve fatigue, and tighten skin. It is a new generation of health equipment suitable for every family. Pulse energy therapy can quickly dredge the lymphatic vessels and meridians of the body, and then play a health and beauty effect.

Negative pressure massage:

- A. The dual purpose of cupping and scraping with high suction and negative pressure can remove dampness and detoxification.
- B. It can dredge the meridians, promote the blood circulation of the human body, open up the acupuncture points of the human body, and play a role in strengthening the body

Breast care massage:

- A. Micro-current: low-frequency bio-current can stimulate breast expansion, provide better nutrients and improve breast elastic fibers and glial tissue, thereby Make the chest soft and stretch.
- B. Vibration: Let the muscle layer and fat layer of the chest move to make The chest is fuller and more stylish.
- C. Automatic suction and release, physical massage on the chest using the principle of vacuum suction and release, and use of different frequency negative pressure circulation modes to dredge the breast meridians
- D. Red light: Red light can improve the local microcirculation of the mammary gland through photochemistry with abundant breast blood vessels and meridians. Warm and comfortable red light can relieve pain and eliminate lumps.



Product operation instructions

Instructions

1. After cleansing the skin, apply essential oil or gel
2. Connect the power supply to boot
3. Choose the cup that needs to work and connect it, and adjust the strength
4. The suction power can be controlled and adjusted according to the demand. It is recommended to start slowly from low
6. The oil outlet must be tightened to prevent air leakage
5. After the time countdown is over, it will automatically stop working or turn it off manually
6. After using the instrument, unplug the power and turn off the instrument. After wiping the instrument, place it in a dry place

Technical specifications

Power input voltage: 110-240V
Power output voltage: 12V
Power output current: 3A
Total instrument power: 36W
Time: 0-30 minutes
Inhalation: 10 levels adjustable
Deflation: 10 levels adjustable
Micro current: 10 levels adjustable
Vibration: 10 levels adjustable
Product material: ABS, PC, hardware

About maintenance

After each use, it is recommended to clean the surface of the machine. When cleaning, please unplug the power plug. When you encounter stubborn stains, you can wipe it with a sponge or cloth moistened with a diluted neutral detergent. Please do not use strong alkali, alcohol, volatile oil, delustering liquid, etc. to wipe. This instrument is not waterproof, and you cannot use it in water or wash the host and power supply with water.

Note, due to irregular product upgrades, the model of the manual may be different from the actual product, this manual is for reference only.

After-sales service

The product is guaranteed for one year from the date of purchase (excluding consumables) in the event of a malfunction during use from the date of sale.

1. Under the premise of following the relevant instructions of this product, in the event of malfunction, free maintenance can be carried out.
2. During the warranty period, if the following conditions are met, it is a paid maintenance.
 - ①Incorrect use, incorrect repair, modification or damage due to malfunction.
 - ②After purchasing this product, it is damaged due to man-made external forces.
 - ③Failure and damage due to abnormal voltage caused by natural disasters.
 - ④ Abnormal connection with other electrical appliances causes product failure.

Warranty Card

1. This product enjoys a free warranty within one year from the date of purchase. It has been used in normal operation and has not been disassembled and repaired. The warranty personnel have confirmed that it is a product quality problem and enjoy the warranty service. Repairs are provided under any normal use.
2. The warranty card should be filled out at the time of purchase. The warranty card is irrelevant. The company's warranty is invalid. Please bring the warranty card when repairing the product.

Customer name						
The way of buying						
Product number						
Purchase date	Year	month	day-	year	month	day
Customer name	Year	month	day-	year	month	day