


# User Manual

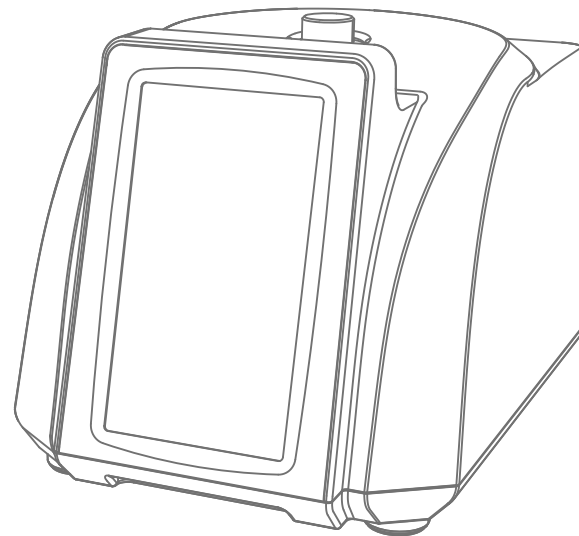

# Contents

Safety Precautions .....	1
Contraindications .....	2
Overview .....	2
Interface Introduction .....	2
Usage Instruction .....	4
Technical Specification .....	5
Use Precautions .....	6
Maintenance & Cleaning .....	6
FAQ .....	7
Ask & Answer .....	8
Operation Guide .....	9
Warranty & Support .....	9

# Safety Precautions

- 1.Please read this manual carefully before use and familiarize yourself with the operation of this instrument and related technology.
- 2.This instrument is a high-precision electronic device that requires proper storage. Please avoid storing it in humid conditions, direct sunlight, or near flammables and explosives to prevent damage.
- 3.Always use the designated accessory provided or recommended by the manufacturer to ensure performance throughout its service life.
- 4.Do not place the instrument in high-temperature conditions, as this may negatively impact its service life and performance.
- 5.Do not apply the instrument on the thyroid, parathyroid, testicles, pregnant women's abdomen, pacemakers, etc.
- 6.Patients should refrain from using this instrument unless they have obtained permission from their doctor.
- 7.Ensure that the instrument is applied on moist skin and avoid the operation on dry skin.
- 8.Power off the instrument when not in use to ensure its safety.
- 9.If you are currently taking anti-adiposity drugs, please discontinue their use for 1 - 2 months before using this instrument. Failure to do so may result in an extended course of treatment.
- 10.To maintain hygiene and extend the instrument's service life, sterilize it with saline.

## Contraindications


**Before using this instrument, it is essential to consult your doctor or seek advice from professionals. Details below:**

1. Pregnant or lactating woman.
2. Individuals with heart disease or a pacemaker.
3. Individuals with unhealed wounds or those in the recovery period from surgery.
4. Individuals suffering from epilepsy, severe diabetes, and hyperthyroidism.
5. Individuals suffering from malignant tumors, hemophilia, and hemorrhage.
6. Individuals with skin disease or infectious disease.
7. Individuals who have had a cesarean section within the last six months, a miscarriage within the last three months, or a natural birth within the last two months should refrain from using this instrument.
8. People who had plastic surgery, with implants or metal objects at the treatment area in 3 months.


## Overview

The Vacuum Butt Enhancer utilizes the latest and safest technology to help you attain your desired body shape. It provides the perfect solution for people seeking improvements in body figures without injections, drugs, or side effects. This external-use device delivers noticeable improvements day by day.


## Interface Introduction

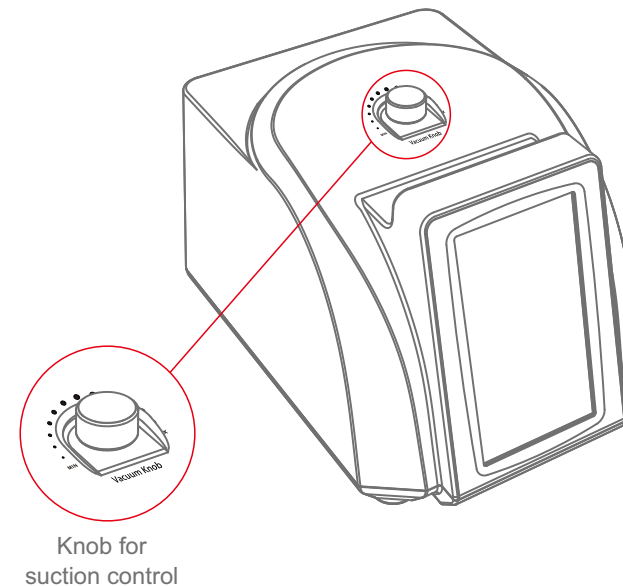

### Welcome Interface

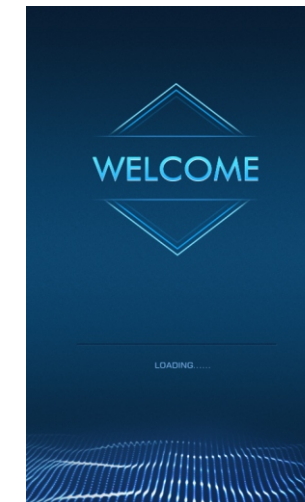

Main interface

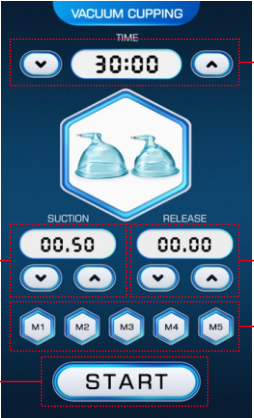

Time display and set buttons (0 - 60 minutes)

Suction time display and set buttons (range: 0.5 - 15)


Release time display and set buttons (range: 0 – 9)

Mode selection:  
M1: Suction only.  
M2: 3 seconds of suction followed by 1 second of release.  
M3: 5 seconds of suction followed by 1 second of release.  
M4: 7 seconds of suction followed by 1 second of release.  
M5: 9 seconds of suction followed by 1 second of release.

START/STOP button


Usage Instruction

How to replace the cotton filter

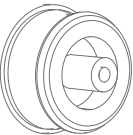


Open the inlet by twisting off the cap, insert a cotton filter, and then reattach and secure the cap by twisting it clockwise.

How to connect the hose




How to connect the Y-shaped hose:  
connect one end of the hose to the inlet on the main unit, and connect another two ends with two butt cups. (In case only one cup is needed for the treatment, tie a knot on the spare end. And untie the knot if 2 cups are needed for the treatment.)

Cup instruction

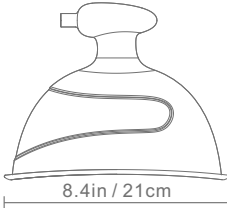


8.4in / 21cm

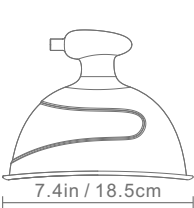

7.4in / 18.5cm

Select the suitable cup based on the size of the buttocks.

Technical Specification

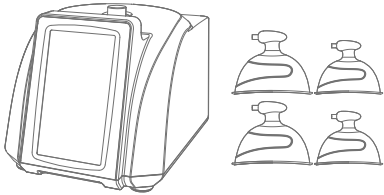	Power supply Input: AC110V-220v 50HZ-60Hz	Net Weight: 4.9KG
		Gross Weight: 6.1KG
Dimensions: 33.5*22*20.6cm		
Input Power: 20W		

	No-load Power: 19.8W	Airflow: 10LPM
	Load Power: 20W	Noise level: < 65db (30cm away)
	Vacuum: -73KPa	

## Use Precautions

- 1.Strictly follow the instructions in the manual to use the instrument and train operator.
- 2.Before using the instrument, use a plug with a grounding pin and ensure the power socket is grounded.
- 3.Make sure the voltage of the instrument is adapted. If the local power supply voltage is unstable, we recommend adding a power-matched regulated power supply between the mains and the instrument.
- 4.Remove all the metal objects inside the body to avoid causing accidents or affecting treatment effects.
- 5.Do not use the device on an empty abdomen, and take the treatment 1 hour later after the meal.
- 6.Always start with the lowest suction intensity.

## Maintenance & Cleaning

- 1.Disconnect the power cord when not in use and store the instrument in a dry, cool, and well-ventilated place.
- 2.To clean the cups, gently wipe them with a wet wipe or saline solution to sanitize, then thoroughly dry them with a clean towel.
- 3.After cleaning the instrument, cover it with a dustcloth or plastic wrap. Before use, ensure to sanitize the handle again for hygienic reasons.

4.Disassembling the instrument and attempting to clean the inside is strictly prohibited. This task should only be performed by professionals.

## FAQ

### 1.Are there any adverse effects linked to vacuum butt enlargement treatment?

A: No, there are no side effects at all. Vacuum butt enlargement therapy is entirely based on physical massage, suction, and release under negative pressure. With the integration of multiple advanced technologies, it provides users with immediate effects without any side effects or rebound effects.

### 2.How does negative pressure aid in detoxification during treatment?

A: Negative pressure can stimulate capillary congestion and boost cellular vitality. The alternation of negative pressure promotes the opening and closing of pores in the targeted area, improving skin breathing, increasing oxygen uptake, and aiding in the elimination of wastes. In summary, this treatment promotes Qi and circulation, dispels dampness and cold, clears meridians, detoxifies, and purges heat.

Ask & Answer

Ask		Cause		Answer
Fail to turn on the instrument.	→	The power cord is disconnected from the wall socket.	→	Check if the power cord is well connected to the wall socket and ensure that the socket is functional.
	→	The fuse is loose or burned.	→	Replace the fuse on the back of the instrument and restart it.
Absence of suction or weak suction from the cup.	→	The hose is not securely attached.	→	Ensure the hose is tightly secured to the main unit.
	→	The cap on the inlet is loose or the rubber ring inside is broken or absent.	→	Secure the cap on the inlet tightly or replace the rubber ring inside.
	→	Check if you have set the knob.	→	Set the knob.

Operation Guide

Vacuum butt enlargement treatment

(Step 1:   20 mins)



Parameter Settings  
Recommended working time: 20 minutes.  
Recommended suction time: 6.0  
Recommended release time: 2.5  
Adjust the suction intensity according to individual endurance.  
Product required: massage cream

Warranty & Support

We provide a one-year warranty and a three-month return and replacement policy from the date of purchase.

Product Name:

Vacuum Butt Enhancer

Product Model:

MS-13H2

Quality Inspection:

PASS