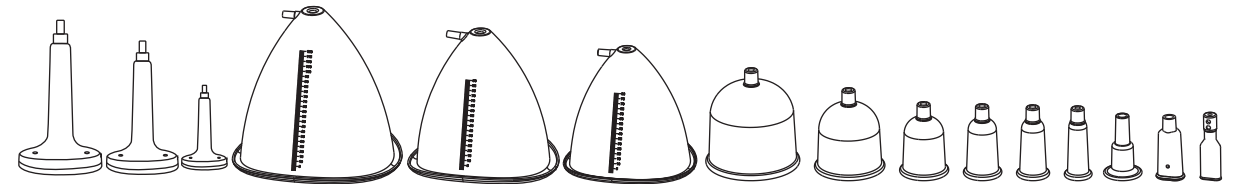
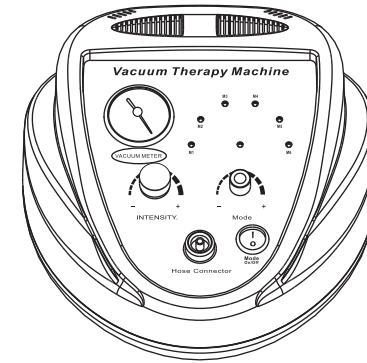
EN ES

User Manual

Vacuum Therapy Machine

Thank you for selecting our product. We trust you will have a delightful experience. For optimal results and a seamless experience, please thoroughly review this manual before use.

Safety Precautions

Contraindications ⓘ

This device is not suitable for everyone. Individuals with the following conditions should either refrain from using it or consult a doctor beforehand.

- 1.Pregnancy or lactation.
- 2.Individuals with heart disease or a pacemaker.
- 3.Unhealed wounds or during the recovery phase of surgery.
- 4.Epilepsy, severe diabetes, and hyperthyroidism.
- 5.Malignant tumors, hemophilia, and hemorrhage.
- 6.Skin disease or infectious disease.
- 7.Cesarean section within the last six months, a miscarriage within the last three months, or a natural birth within the last two months.
- 8.Had plastic surgery within 3 months or with implants or metal objects in the treatment area.

Precautions ⚠

- 1.Adhere strictly to the instructions in this manual when operating this device.
- 2.Ensure the target area is moist during use; avoid operation on dry skin.

- 3.If you are currently taking weight-loss medication, delay using this instrument for 1-2 months after discontinuing the medication, or an extended treatment course will be necessary.
- 4.Disassembling this device or replacing any components is prohibited.
- 5.Do not attempt to repair this instrument on your own.
- 6.Always place the instrument on a flat and stable surface.
- 7.Ensure the instrument is powered off before cleaning.
- 8.Avoid applying excessive essential oil or makeup during use.
- 9.Always set the Vacuum Intensity to the minimum level after use.
- 10.Limit the duration of a single treatment to within 30 minutes.
- 11.Protect the instrument from humidity, direct sunlight, inflammables, and explosives.
- 12.Do not place the instrument in high-temperature conditions, as this may negatively impact its service life and performance.
- 13.Avoid using this device on sensitive areas, including the eyes, thyroid, parathyroid, testes, the abdomen of a pregnant woman, or areas with a pacemaker installed.

Summary

Vacuum butt lift uses suction and mechanical massage to "lift" your buttocks and improve the overall appearance of your skin. Vacuum therapy is non-invasive, with minimal side effects. Treatment sessions typically last 30-45 minutes, once or twice a week. The best results usually require around six to 10 sessions. Monthly treatments are necessary to maintain results of a vacuum butt lift.

Vacuum Butt (Breast) Lift Therapy

During a vacuum butt(breast) lift, one suction cup is placed on each butt cheek(breast). Then turn on the machine and adjust the suction that is comfortable for you. The suction cups are clear to allow your beauty instructor to see the responding of your skin.

Therapy times generally last 20-35 mins. After therapy, you are advised to drink plenty of water to help flush toxins from your body. If you feel discomfort after therapy, rest until discomfort subsides and reduce exercise.

Initial treatments is advised 1 or 2 times per week for 6 to 10 sessions. After this point, monthly sessions are recommended to help maintain your results.

Benefits of The Vacuum Therapy Machine

Vacuum therapy was originally developed to treat scarring from burns, and it originated from cupping therapy—an ancient Chinese medicine treatment that uses vacuum cups to create suction on the skin to increase blood flow to a specific area of the body. During these treatments, it was also noticed that vacuum therapy reduced the appearance of cellulite.

- *Improved blood circulation, increasing lymphatic flow to remove toxins and water retention
- *exfoliating the skin, thereby making it smoother and more toned in appearance
- *stimulating the middle layers of skin for increased toning effects
- *decreasing the appearance of cellulite, smoothing skin
- *Increased skin elasticity (less "sagging")
- *Decreased muscle tension and soreness
- *Tightening, Lift and larger, rounder-looking butt(breast)

When performed for several sessions on a weekly basis, you may notice that your skin is tighter and more toned. This can provide the appearance of a butt lift.

How Does Vacuum Therapy Feel?

Vacuum therapy uses suction and mechanical massage. While this treatment is not painful, you might feel some stretching or "pulling" sensations during your session. Your provider can adjust the amount of suction to make you comfortable during treatment.

Vacuum Massage Devices For Cellulite

Cellulite is a term used to describe an area of the body where skin has a dimpled appearance. According to research, 80-90% of adult women have some cellulite and it's generally seen on the hips, buttocks, and thighs. It can also occur on the lower belly and upper arms.

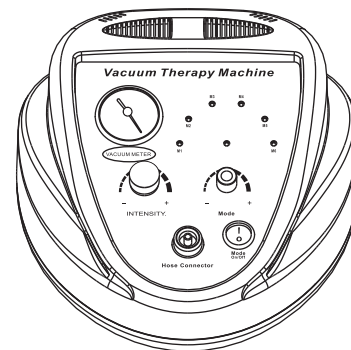
Vacuum massage may be able to improve the appearance of cellulite by:

- *Draining excess body fluid
- *Redistributing fat cells
- *Improving circulation
- *Plumping up skin


However, vacuum massage won't cure cellulite. While massage may improve appearance, results typically don't last for long and in many cases lasting treatments are necessary.

Our vacuum therapy machine device lifts, stretches, and rolls skin to help move fat and reduce cellulite. Meanwhile, we have received ongoing feedback from customers who have shared their experiences of cellulite reduction and disappearance after using our vacuum therapy machine.


Packing List


Main machine

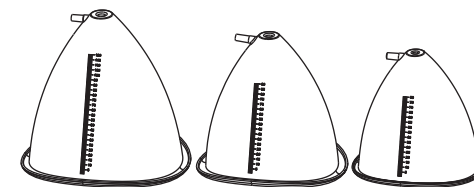
Breast Cups(in 3 sizes) X 6

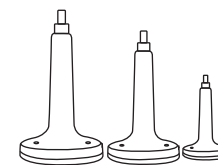
Metal Vacuum Cups x 3

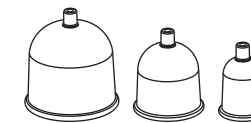Cupping Cups (in 3 sizes) x 6

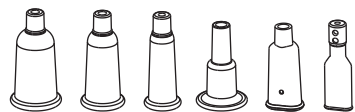
Gas Cups (in 6 specifications) x 12

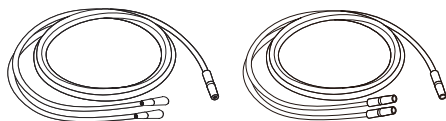
Y-shaped Hose (in 2 sizes) x 2

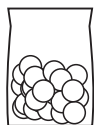spare rubber ring x 2

Big Cotton Filter
X 1 Pack(30PCS)

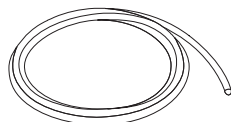
White Hose x 1

Fuse x 4

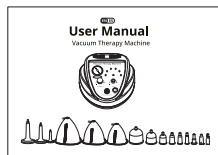
Power Cord x 1

User Manual x 1

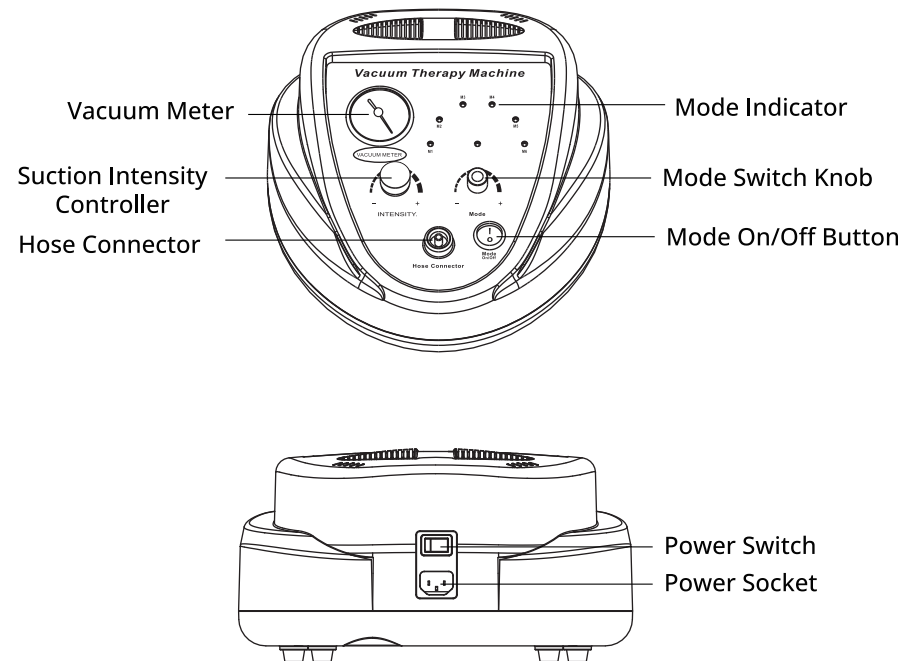
Product Overview

Accessories Overview

1. Cups

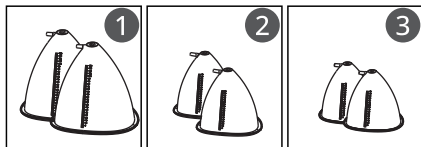1.1. Breast/Butt Cups(in 3 sizes)

Cup No.1 – 3: Choose based on your size for breast augmentation and butt lift.

1.2. Cupping Cups(in 3 sizes)

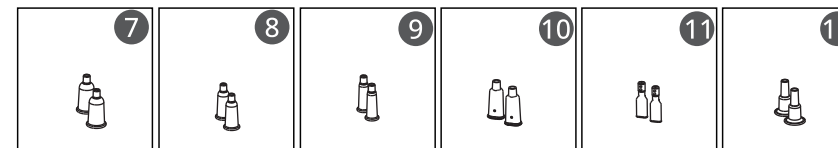
Cup No.4: Used for body cupping therapy on large areas like the buttocks, hips, and abdomen.

Cup No.5: Designed for body cupping therapy and moving cupping therapy on small areas like legs and back.

Cup No.6: Crafted for body cupping therapy and moving cupping therapy on small areas like arms and lower legs.

1.3. Gas Cups(in 6 specifications)

For No.10 cups, locate the mini hole and cover it with your thumb to create suction during use. Remove your thumb to pause.

Cup No.7: Designated for moving cupping therapy on the neck and arms.

Cup No.8: Used for moving cupping therapy and massage on the face and neck.

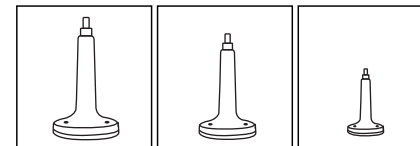
Cup No.9: Intended for facial cleansing and massage.

Cup No.10: Suitable for facial cleansing and lifting.

Cup No.11: Designed for facial cleansing, comedone removal, and eye bag reduction.

Cup No.12: Used for private area massage.

2. Massage Rollers


Large Roller: Applicable for back, legs, and abdomen.

Medium Roller: Designed for the waist, inner thighs, arms, and lower legs.

Small Roller: Crafted for the face, neck, breast, and small body areas.

3. Hoses

<p>Wide Y-Shaped Hose</p> <p>The wide Y-shaped hose is designed for vacuum cup use.</p>		
<p>Narrow Y-Shaped Hose</p> <p>The narrow Y-shaped hose is designed for breast cup use.</p>		
<p>Straight Hose</p> <p>The straight hose is designed for the massage roller.</p>		

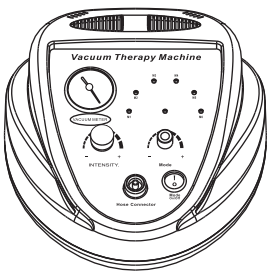	<p>Input Voltage: 110v-220v 50HZ/60HZ</p>	<p>Air flow: 8L/min</p>
	<p>Power: 7W</p>	<p>Noise level: < 65db (30cm away)</p>
	<p>Dimensions: 12" x 11.4" x 6"</p>	<p>Net Weight: 3.5kg</p>
	<p>Vacuum Suction: -50cmHg</p>	<p>Gross Weight: 4.5kg</p>

Suction and Release Pattern of Each Mode

(The default mode is Mode 0, select your desired mode from 1 to 6 using the mode switch knob.)

Mode	Suction Duration	Target area
Mode 0	Only suction	Face, Body
Mode 1	2-sec suction, 2-sec release	Face
Mode 2	3-sec suction, 2-sec release	Face
Mode 3	4-sec suction, 2-sec release	Body
Mode 4	5-sec suction, 2-sec release	Body
Mode 5	6-sec suction, 2-sec release	Breast, Butt
Mode 6	11-sec suction, 2-sec release	Breast, Butt

Working Mode and Pairing Cups

Cup/Mode	Only suction	Mode 1	Mode 2	Mode 3	Mode 4	Mode 5	Mode 6
	✓					✓	✓
	✓						
	✓			✓	✓		
	✓			✓	✓		
			✓	✓			
	✓	✓	✓				
	✓	✓	✓				
	✓						
	✓						
		✓	✓				
Massage Roller	✓			✓	✓		

Instructions

Breast Enlargement(20 minutes):

Step 1: Select the appropriate cups based on your breast size.

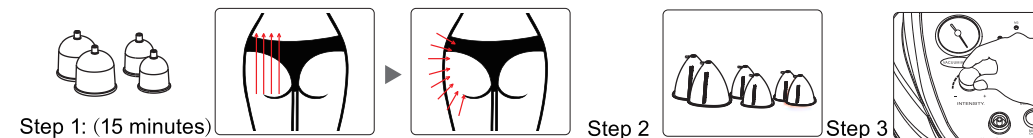
Step 2: Secure the cups firmly on your breasts, ensuring there is no gap to avoid air leakage.

Step 3: Set the appropriate suction intensity.

Step 4: Initiate the treatment.

Butt augmentation:

Follow the illustrations.

Step 1: (15 minutes)

Step 2

Step 3

Step 4: (20 minutes)

Step 1: Massage the buttocks using cups No.4 and No.5, moving from bottom to top and from outside to inside. This step relaxes the muscles and prepares the buttocks for the next vacuum therapy.

Step 2: Select the appropriate cups based on your butt size.

Step 3: Set the appropriate suction intensity.

Step 4: Initiate the treatment.

Face:

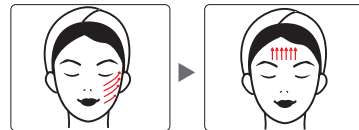
Follow the illustrations.

Step 1: (3 minutes)

Step 2: (10 minutes)

Step 3: (10 minutes)

Step 1: Cleanse your face using cups No.10 and No.11. Firmly apply the cup on your forehead for 1 second, then glide upward to the hairline. (Locate the mini hole on the cup and cover it with your thumb to create suction. Remove your thumb to pause.) Repeat the above operation on your nose, cheeks, and chin.

(3 minutes)

Step 2: Apply cups No.8 or No.9 to your face for firming purposes. (10 minutes)

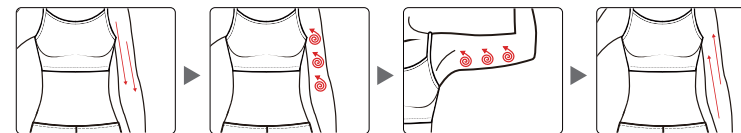
Step 3: End with the small massage roller to enhance blood circulation and detoxification. (10 minutes)

Arm:

Follow the illustrations.

Step 1: (15 minutes)

Step 2: (15 minutes)

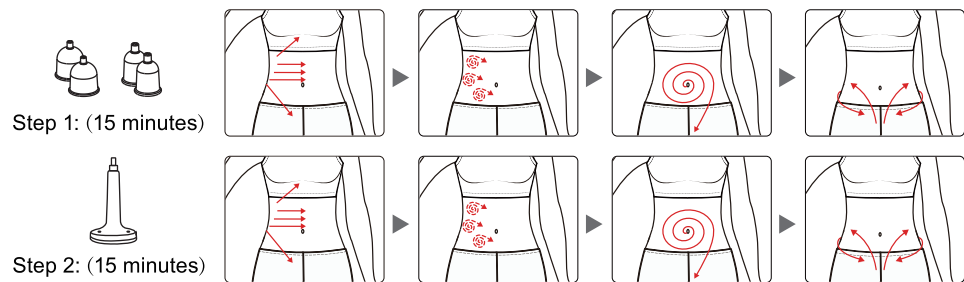
Step 1: Apply cups No.5 or No.6 for clearing dampness and detoxification. (15 minutes)

Step 2: Use the large massage roller to clear channels and enhance blood circulation. (15 minutes)

Belly:

Follow the illustrations.


Step 1: Apply cups No.5 or No.6 for clearing dampness and detoxification. (15 minutes)

Step 2: Use the large massage roller to clear channels and enhance blood circulation. (15 minutes)

Back:

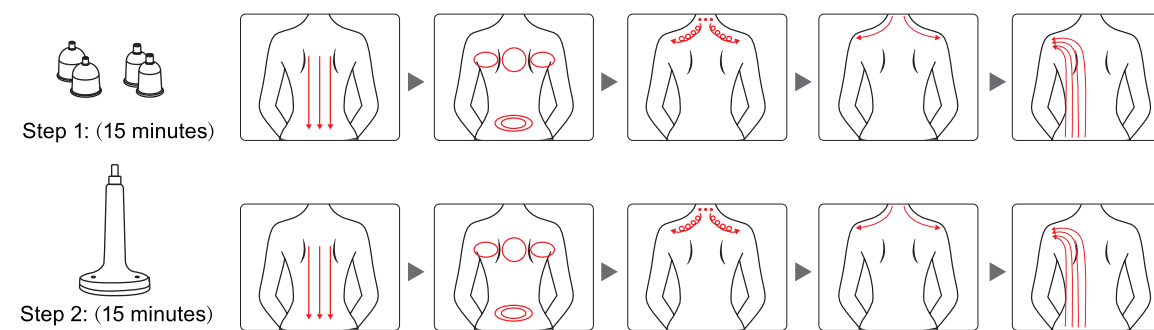Follow the illustrations.

Step 1: Apply cups No.5 or No.6 for clearing dampness and detoxification. (15 minutes)

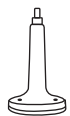
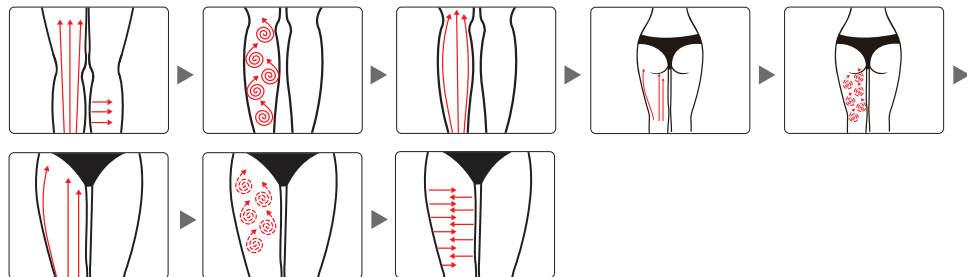
Step 2: Use the large massage roller to clear channels and enhance blood circulation. (15 minutes)

Legs:


Follow the illustrations.

Step 1: (15 minutes)

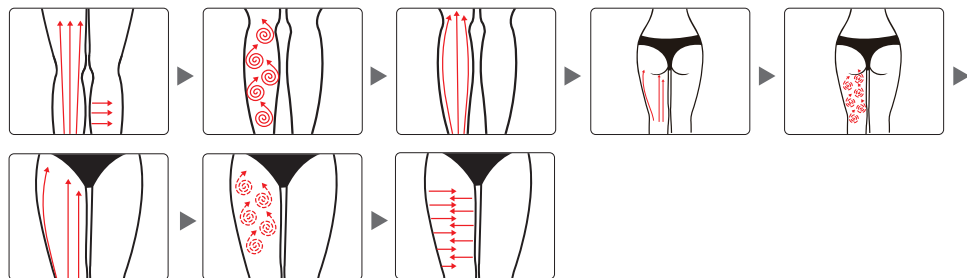Step 2: (15 minutes)

Step 1: Apply cups No.5 or No.6 for clearing dampness and detoxification. (15 minutes)

Step 2: Use the large massage roller to clear channels and enhance blood circulation. (15 minutes)

Troubleshooting

Ask	Cause	Answer
Fail to turn on the instrument.	The power cord is disconnected from the outlet.	Check if the power cord is securely connected to the outlet and ensure that the outlet is functional.
	The fuse is loose or burned.	Replace the fuse on the back of the instrument and restart it.
Weak or no suction from the vacuum cup/massage roller.	The hose is not securely attached.	Ensure the hose is tightly secured to the main unit.
	The cap on the hose connector is loose, or the rubber ring inside is broken or absent.	Firmly secure the cap or replace the rubber ring.
	Check if you have set the suction controller.	Set the suction intensity.

If you have any questions or concerns, please don't hesitate to contact us. We will do our best to ensure your satisfaction.

Warranty and Support

UNOISETION backs up its products with 100% satisfaction guarantee. Deficient products without being improperly used, abused, altered or repaired will be replaced in a timely manner. To validate your warranty, contact us at service.unoisetion@outlook.com WITHIN 12 MONTHS since your purchase date.

1. Provide 12 months warranty service from the date of purchase, but the damage caused by time and use is not within the scope of the service warranty.
2. The warranty period of the product is calculated from the date of purchase by the customer.
3. If the warranty period is exceeded, the consumer needs to bear the cost of parts replacement and maintenance.
4. The warranty of this product does not include any fault caused by:
 - A) Natural disasters, other unpredictable factors, or any damage caused by force majeure.
 - B) Improper use leads to product damage, improper storage leads to product damage, loss and appearance damage.
 - C) Accessories such as diamond microdermabrasion tips and other consumables belong to normal loss and are generally not covered by warranty.
 - D) Damage caused by disassembly, modification and repair of the product by unauthorized personnel of the company.
5. For after-sales service or any questions, please contact us by Email: service.unoisetion@outlook.com

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