

User Manual

SR-NF005

UNOISETION Microdermabrasion Machine

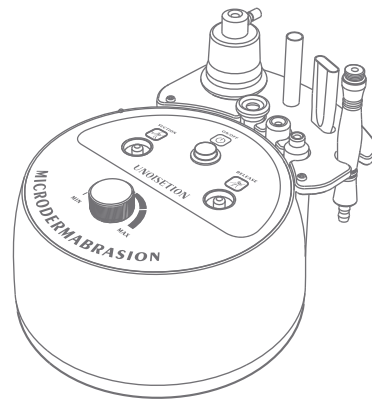


Table of contents

EN

Overview	1
Benefits	1
Packing List	2
Control Panel Instruction	4
User Guide	5
Recommended Application	13
Troubleshooting	17
FAQ	18
Precautions	20
Contraindications	21
Technical Parameters	21
Warranty and Support	22

Overview

Thank you for choosing UNOISETION Microdermabrasion Machine!
Please read this manual carefully before use!

Featuring fine particles on the microdermabrasion tip, this device gently exfoliates the thick, uneven skin on the surface and removes it through suction, promoting skin revitalization. Follow with the glass tube to extract impurities from the pores using negative pressure for a deep cleansing effect. Complete the process with the spray bottle to hydrate and nourish the skin.

Benifits

1. Removes age spots and unclogs pores
2. Improves hyperpigmentation
3. Exfoliates and rejuvenates the skin
4. Reduces fine lines and wrinkles
5. Minimizes or eliminates enlarged pores
6. Reduces comedones

Packing List



Main Machine x 1



Hose x 1



Microdermabrasion Handle x 1



Microdermabrasion Tips x 3



Glass Tubes x 2



Spray Bottle x 1



Rack x 1



Cotton Filter x 1 Pack(50PCS)



Spare Rubber Rings x 3



Screws x 2



L-shaped Hexagon Wrench x 1




Adaptor x 1

Note: ⚠

Please inspect the device and its accessories for completeness before installation.


Control Panel Instruction

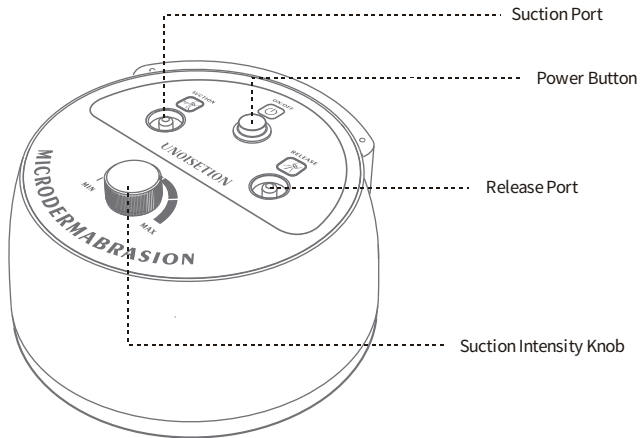


User Guide

1. Installation


1.1 Install the rack.

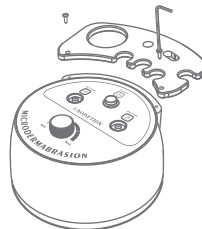

2.1 Connect the power.

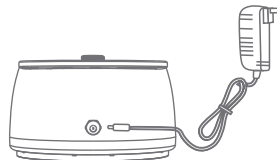


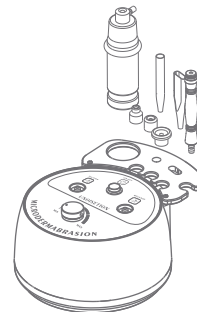


3.1 Place all the components as shown.



2. Microdermabrasion

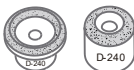2.1 How to Select the Microdermabrasion Tip?


<p>22MM D-240 15MM D-240</p>  <p>For large facial areas.</p> 		
<p>6MM D-240</p>  <p>For small areas around the eyes or nose.</p> 		

Note : ⚠


Do not use the microdermabrasion tip on the same spot for more than 3 seconds.

2.2 How to Install the Microdermabrasion Handle and Tips?


1. Check the microdermabrasion handle for the pre-installed black rubber ring; if it's missing or damaged, replace it.

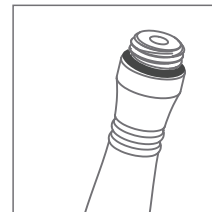


2. Put a cotton filter on the microdermabrasion handle port.

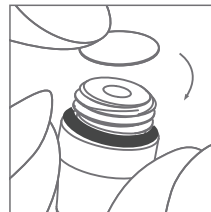

3. Attach the desired tip to the handle.

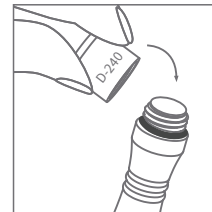

4. Connect one end of the hose to the handle.

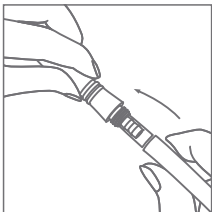


2.3 How to Clean Microdermabrasion Tip?

Clean the tip with water or soapy water after each use. Ensure to remove the microdermabrasion tip from the handle.

Note : ⚠

- ① Replace the tip every 3 months.
- ② Change the cotton filter in the tip every 3 sessions.

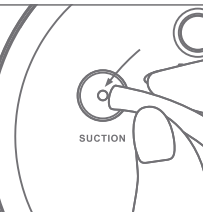


5.1 Connect the other end of the hose to the suction port on the main machine.


3. Glass Tube

3.1 How to Select the Glass Tube?

Flat Glass Tube: For facial massage.

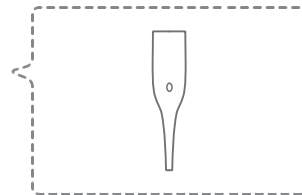

Straight Glass Tube: For clearing severely clogged pores.

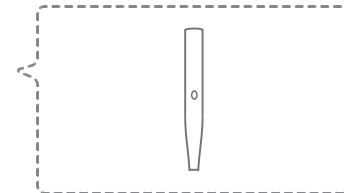


3.2 How to Install the Glass Tube?

1. Connect one end of the hose to the glass tube.




2. Attach the other end to the suction port on the main machine.

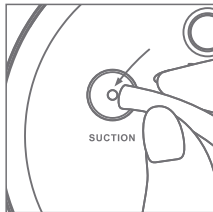

3.3 Caution ⚠

Cover the hole on the glass tube with your finger to prevent air leakage during use.

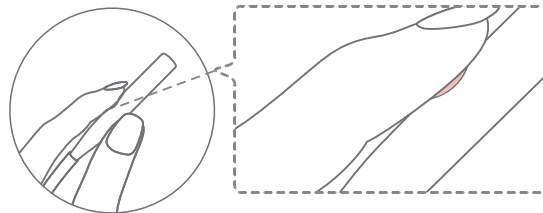


4. Spray Bottle

4.1 How to Install the Spray Bottle?


1. Prepare the spray bottle with pure water, toner, saline, or a mixture of water and serum in a 10:1 ratio.




2. Connect one end of the hose to the release port on the main machine.

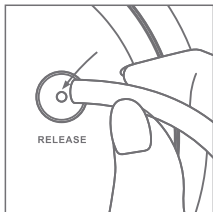

4. Cover the hole at the top of the bottle during use.

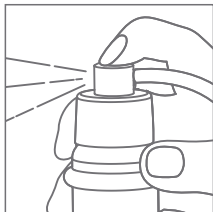


3. Attach the other end to the spray bottle.

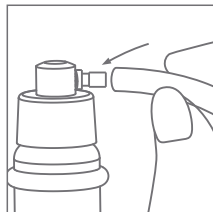


4.2 Precautions: ⚠

Avoid using high-concentration solutions.

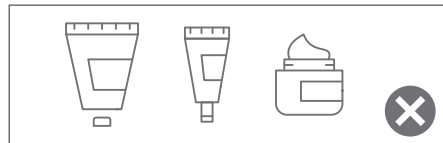


Recommended Application


1. Open pores with a facial steamer or hot towel.

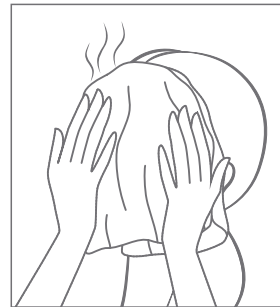


2. Power on the device.
The button lights up red when the microdermabrasion machine is in operation.

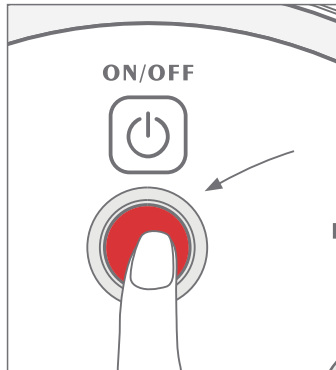


3. Adjust the suction intensity to a comfortable level.

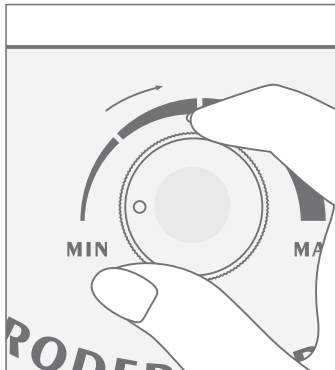

4. Exfoliate and reduce wrinkles with the microdermabrasion tip.

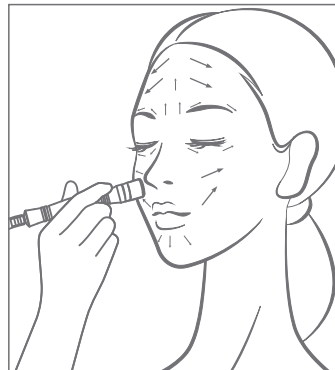


5. Remove comedones using the glass tube.

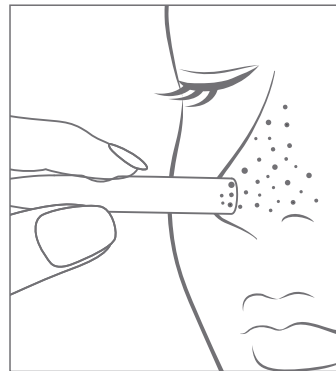

6. Hydrate the skin with the spray bottle.



7. Complete the procedure with a facial mask.

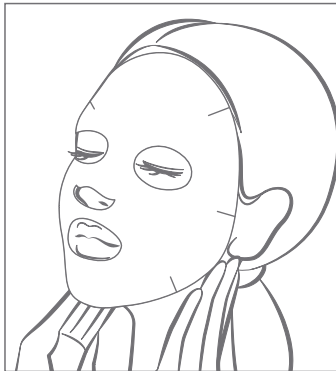

Troubleshooting

1. 1 No suction

1. Check if the hose is properly installed in the suction port.
2. Ensure the hose is tightly connected to the handle.
3. Inspect the rubber ring; if missing or damaged, replace it.
4. Check the hose for any damage; replace if necessary.
5. Verify if the hole on the glass tube is blocked; clean it if needed.

1. 2 No mist from the spray bottle

1. Ensure the liquid is not too concentrated.
2. Press the hole on top with your finger during use.

Note: ⚠

If the problem persists, please contact us for further assistance.

FAQ

1

Q: Why does my face get red after Microdermabrasion?

A: Microdermabrasion works on the principle of negative vacuum pressure. When the microdermabrasion handle or glass tube is applied to your skin, it can cause temporary capillary congestion, resulting in redness. This is normal and should fade within a few hours. It is recommended to moisturize your skin with facial cream after using the microdermabrasion handle or glass tube.

2

Q: How long does it take to see results?

A: After the first session, your skin will feel softer and smoother. With consistent treatment over 8-12 weeks, you may notice visible improvements in fine lines, wrinkles, spots, and hyperpigmentation.

3

Q: How often should I use the Microdermabrasion Machine, and how long should each session be?

A: It is recommended to use the Microdermabrasion Machine 1–2 times per week, with each session lasting 15–30 minutes. Adjust the frequency and duration based on your skin type and personal expectations.

4

Q: Which skin types are suitable for the Microdermabrasion Machine?

A: The Microdermabrasion Machine is suitable for all skin types, including normal, dry, oily, and combination skin. However, individuals with sensitive skin or other skin concerns should consult a professional doctor or therapist for advice.

Precautions

1. Please read this manual carefully before use so you are familiar with the use of this product and related technology.
2. This product is guaranteed by the designated after-sales service unit.
3. Always begin microdermabrasion treatments at the lowest level until your skin adjusts to the treatment.
4. The Microdermabrasion Machine is a high-precision electronic instrument, please keep and store it properly. Do not store it in a place that is damp, exposed to direct sunlight, or around with inflammable and explosive materials to prevent damage.
5. Power off the instrument when not in use, and ensure the general power is turned off after daily use for safety reasons.
6. Avoid using the Microdermabrasion handle and spray bottle at the same time. Disconnect the hose after using the microdermabrasion handle or glass tube, and then connect the spray bottle to the release port with the hose.

Contraindications

Please do not use the machine or consult a doctor or professional before use if you:

1. Have sensitive skin or experience skin irritation.
2. Have facial implant.
3. With extensive facial trauma.
4. Suffer from skin disease or infectious diseases.
5. Recently had injections, such as hyaluronic acid, mesotherapy, botulinum, or plastic surgery.

Technical Parameters

Model:	SR-NF005
Voltage:	110V-240V
Power:	6W
Dimensions:	8.8 x 7 x 4.1 inch (224*177*105mm)
Gross Weight:	1.4kg
Net weight:	1.1kg

Warranty and Support

A one-year warranty is provided from the date of purchase.

If you have any questions or issues with the product, please contact us via email. We are glad to provide you with additional skincare guidance.

Email: mychwayuk@outlook.com

Product Name: Microdermabrasion Machine
Product Model: SR-NF005
Production Date: _____
Quality Inspection: _____

PASS

How to Use SR-NF005:

Scan the QR code to watch the video or search "SR-NF005" on YouTube



How to Install SR-NF005:

Scan the QR code to watch the video or search "SR-NF005" on YouTube

