


## Face Cupping



## Smart Face Cupping Device Product Manual


## Package Includes

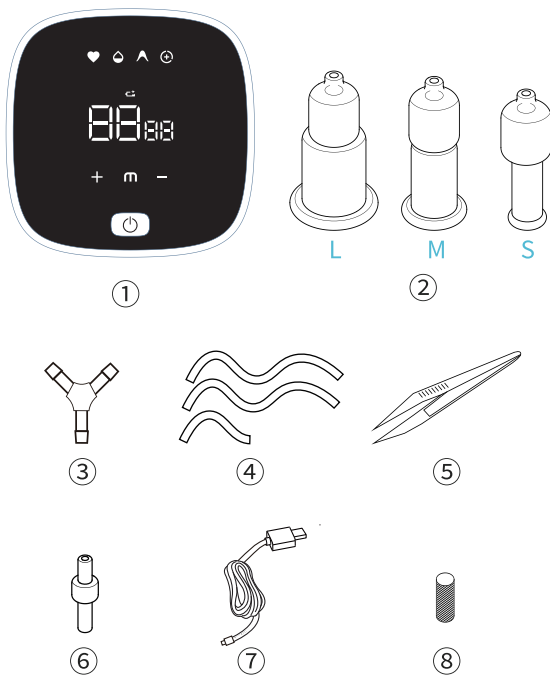

- ① Host \*1
- ② Glass cup\*6
- ③ Trident tube\*1
- ④ Hose\*3
- ⑤ Tweezers\*1
- ⑥ pp adapter\*2
- ⑦ Charging cable\*1
- ⑧ Spare anti-reverse flow cotton\*6

01

## Display and Use

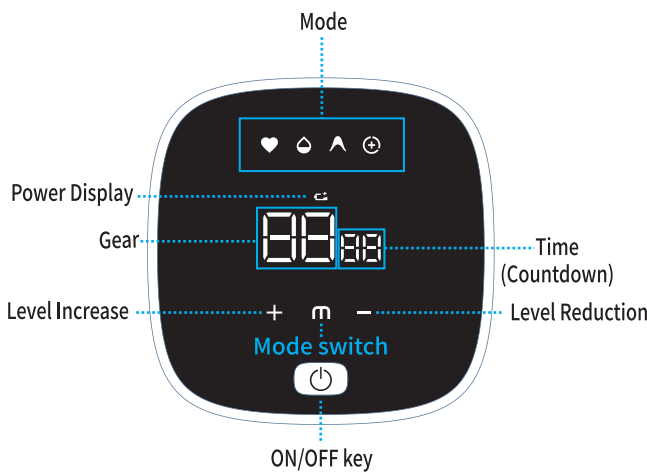

## Mode Introduction

### Cupping mode

Automatically stop after 20 minutes, recommended use for 10 minutes each session.



### Acquaintance mode

Automatically stop after 20 minutes, recommended use for 10 minutes each session.



### Massage mode

Automatically stop after 20 minutes, recommended use for 10 minutes each session.



### strengthen

Automatically stop after 20 minutes, recommended use for 10 minutes each session.



- 1. Cupping mode:** (Default setting is 20 minutes). Use of physical suction long suction skin, stimulate acupuncture points, play a therapeutic role.
- 2. Acquaintance mode:** (Default setting is 20 minutes). Massage the skin with slow-paced suction and alternating releases.
- 3. Massage Mode:** (Default setting is 20 minutes). Massage the skin with fast-paced suction and alternating releases.
- 4. Strengthen Mode:** (Default setting is 20 minutes). Used to strengthen training and prevent rebound.

## Basic Usage

1. Apply cream or nourishing oil to your face.
2. Massage the face and neck, start at the chin and use strokes back along the jaw to the ear.
3. Follow with more strokes from the ear along the sinuses to the nose.
4. Now the forehead, upward strokes start between the eyebrows, followed by outward strokes moving towards and along the hairline, before dropping down to the neck to the collarbone. Tone the skin afterward.
5. Use 2-3 times a week for best results.

03

## Tips

- Please clean the glass cups and keep it dry before use.
- It is recommended to start from the first level. Apply essential oils before use.
- After cupping/scraping, it is normal to appear red marks similar to "bruises" due to skin congestion, and it usually disappears in about 5-7 days.
- Patients with skin allergies, varicose veins, pregnant women, heart disease, abnormal blood pressure, malignant tumors, and cerebrovascular cannot use.

## Disinfection Guidelines

- Straw, trident tube, anti-reverse flow cotton, host, charging cable, and pp adapter, can not be sterilized.
- Glass tips should not be boiled when sterilizing. It is recommended to use UV, alcohol and other non-high temperature sterilization.
- Wash with warm water. The temperature difference between hot and cold does not exceed 20 °C.


04

## Installation Steps

### Single-sided use




### Bilateral use

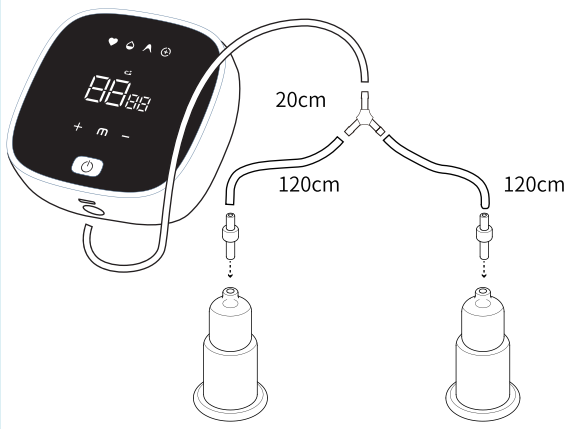

05

## Installation Steps for Spare Anti-reverse Flow Cotton

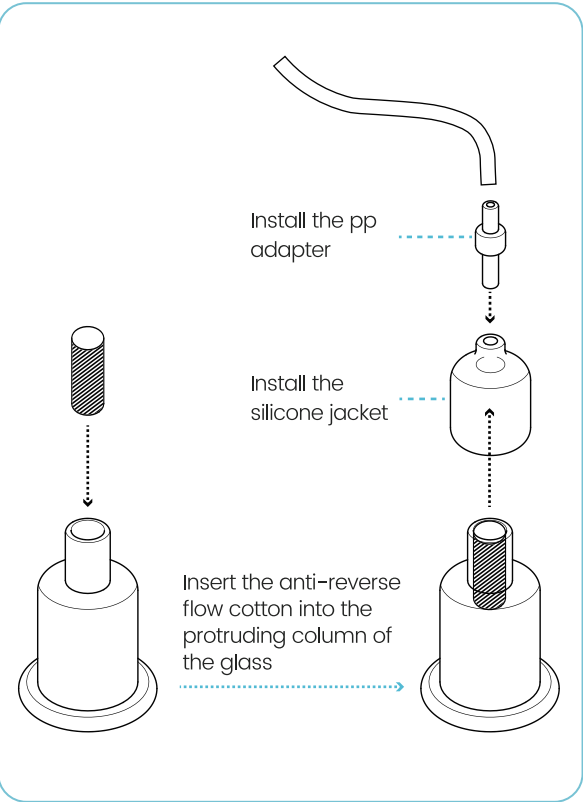

### NOTE:

When you need to apply essential oils for use, you must make sure that anti-reverse flow cotton has been installed inside.

06

## After-sales Service

Our products have a 1-year warranty period. If any quality issues are found within the validity period, please contact us promptly.

07