
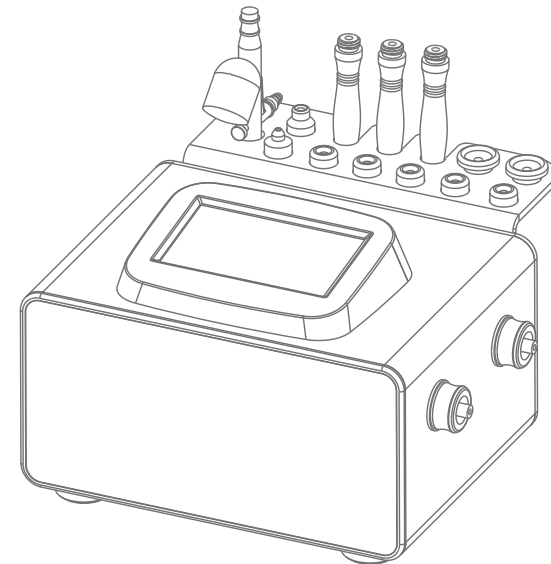


VT-3271

# User Manual

## MultiLift UltraSkin

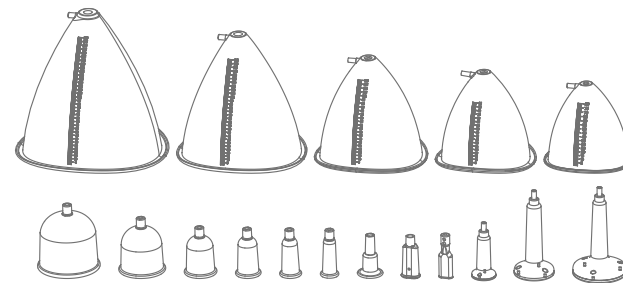

Thank you for purchasing our product.

- Please read this User's Manual carefully for information on how to use the device properly and safely.
- After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
- Make sure to include this User's Manual if you are transferring ownership of the device to another person.

## Contents

Safety Precautions.....	1
Contraindications .....	2
Overview .....	2
Interface Introduction & How to Use...	3
Technical Specification .....	10
Use Precautions.....	11
Maintenance & Cleaning .....	11
FAQ .....	12
Ask & Answer .....	14
Operation Guide.....	15
Warranty & Support .....	21

## Safety Precautions

- 1.Please read this manual carefully before use and familiarize yourself with the operation of this instrument and related technology.
- 2.This instrument is a high-precision electronic device that requires proper storage. Please avoid storing it in humid conditions, direct sunlight, or near flammables and explosives to prevent damage.
- 3.Always use the designated accessory provided or recommended by the manufacturer to ensure performance throughout its service life.
- 4.Do not place the instrument in high-temperature conditions, as this may negatively impact its service life and performance.
- 5.Do not apply the instrument on eyes, thyroid, parathyroid, testicles, pregnant women's abdomen, pacemakers, etc.
- 6.Patients should refrain from using this instrument unless they have obtained permission from their doctor.
- 7.Power off the instrument when not in use to ensure its safety.
- 8.If you are currently taking anti-adiposity drugs, please discontinue their use for 1 - 2 months before using this instrument. Failure to do so may result in an extended course of treatment.
- 9.To maintain hygiene and extend the instrument's service life, sterilize it with saline.

## Contraindications

Before using this instrument, it is essential to consult your doctor or seek advice from professionals. Details below:


1. Pregnant or lactating woman.
2. Individuals with hypertension, heart disease, or a pacemaker installed.
3. Individuals with unhealed wounds or those in the recovery period from surgery.
4. Individuals suffering from epilepsy, severe diabetes, and hyperthyroidism.
5. Individuals suffering from malignant tumors, hemophilia, and hemorrhage.
6. Individuals with skin disease or infectious disease.
7. Individuals who have had a cesarean section within the last six months, a miscarriage within the last three months, or a natural birth within the last two months should refrain from using this instrument.
8. For your safety, refrain from using if you've had plastic surgery or have implants/metal objects near the treatment area within 3 months.



## Overview

Utilizing the most sought-after technique for facial cleansing, hydration, and body sculpting, this instrument guarantees a radiant transformation for both the face and body. It provides the perfect solution for people seeking improvements in skin appearance and body figures without injections, drugs, or side effects. For external use with immediate effects.

## Interface Introduction & How to Use

### Interface Instruction

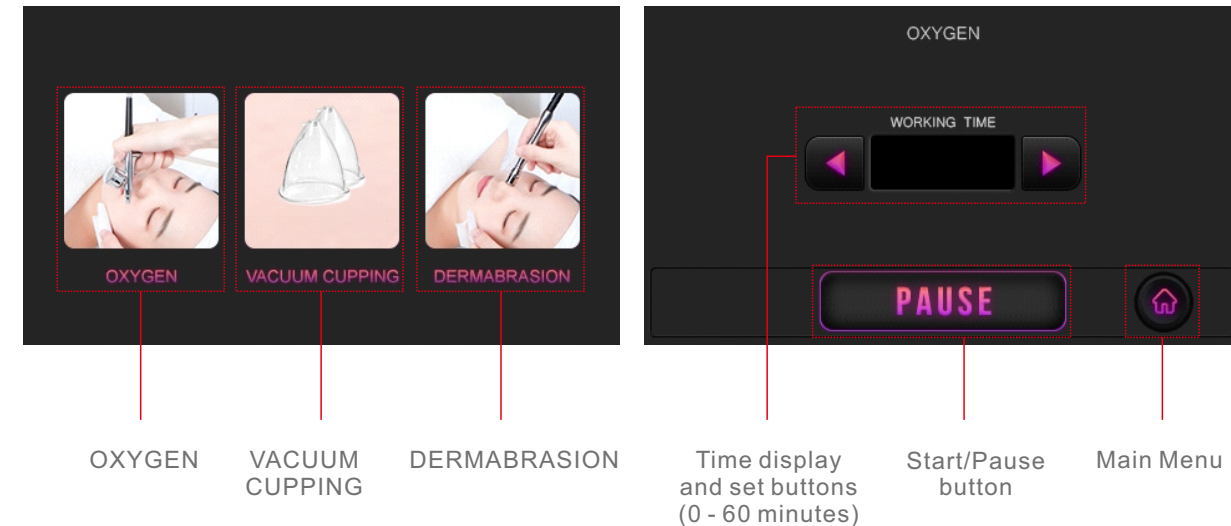

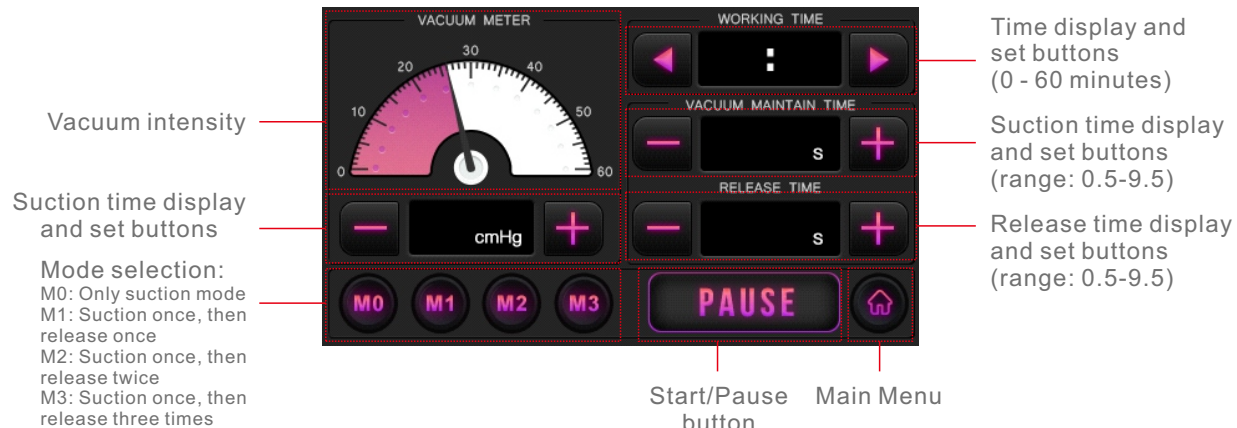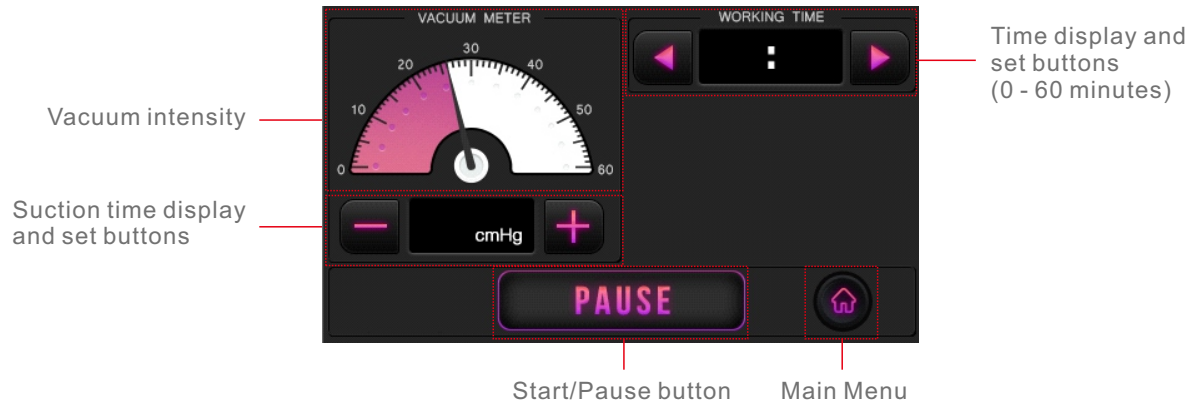

## How to use

### Connect the transparent hose



The long Y-shaped hose is designed for vacuum cup use.



The short Y-shaped hose is designed for breast cup use.




The straight hose is designed for grease cup use.

How to connect the Y-shaped hose: connect one end of the hose to the inlet on the main unit, and connect another two ends with two vacuum cups/breast cups.

(In case only one cup is needed for the treatment, tie a knot on the spare end. And untie the knot if 2 cups are needed for the treatment.)

How to connect the straight hose: connect one end of the hose to the inlet on the main unit, and connect another end with the grease cup.

## How to choose the breast cup/vacuum cup/grease cup

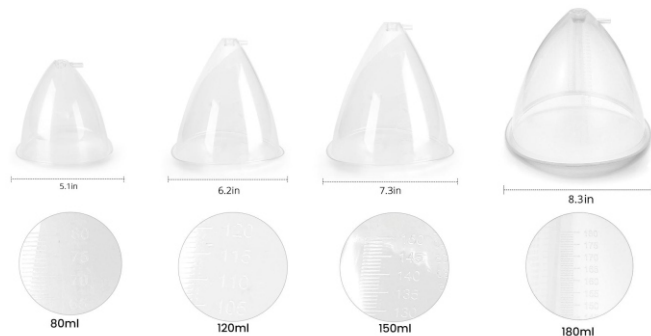

Breast cup, choose according to the size of the breast or buttocks.



- ①Cup No. 1 is used for cupping therapy, moving cupping therapy (for large body areas, e.g., thighs, buttocks, abdomen)
- ②Cup No. 2 is used for cupping therapy, moving cupping therapy (for body areas like legs, and back)

- ③Cup No. 3 is used for cupping therapy, moving cupping therapy (for small body areas, e.g., arms, lower legs)
- ④Cup No. 4 is used for moving cupping therapy (for body areas like the neck, and arms)
- ⑤Cup No. 5 is used for moving cupping therapy, massage (for body areas like the face and neck)
- ⑥Cup No. 6 is used for moving cupping therapy, cleansing (for face area)
- ⑦Cup No. 7 is used for the secret area.
- ⑧Cup No. 8 is used for facial cleansing.



Grease cup

- ①The large cup is used for large body areas, such as arms, abdomen, abdomen, waist, buttocks, legs, and back.
- ②The medium cup is used for the face or parts of the body area, such as arms, neck, face, and lower legs.
- ③The small cup is used for the facial area, such as cheeks, around the eyes, forehead, chin, and neck.

## About Dermabrasion



The black rubber ring must be correctly installed to prevent issues like no suction or weak suction.



Installing the cotton filter is essential as it helps prevent foreign particles from damaging the device and also indicates the cleansing effectiveness.



Select the appropriate working tip based on the size of the treatment area.



Large tips for the large treatment area like the area around the eyes, nose wings, cheeks, and body area.



Round tips are used for the cleansing of the forehead, jaw, nose wings, cheek, and facial wrinkle removal.



Pointed tips are used for blackhead and mite removal in the nasal folds, unclogging pores, and wrinkle removal.


## About OXYGEN Sprayer



(1) Fill the mini-flask up to two-thirds capacity with normal saline and serum, maintaining a ratio of 10:1.

(2) Adjust the spray with the knob at the end of the sprayer.

## Technical Specification

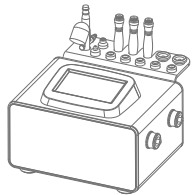	Power supply Input	AC110v-240v 50HZ-60Hz	Net Weight	5.3KG
	Power	14W	Gross Weight	6KG
	Dimensions:13×12.5×7.5inches			
	Input Power:14W			
OXYGEN	Power:14W			
VACUUM CUPPING	Power:14W	Air flow:7LPM		
	Vacuum:-60cmHg	Noise level:<65dB(30cm away)		
DERMABRASION	Power:14W	Air flow:7LPM		
	Vacuum:-60cmHg	Noise level:<65db(30cm away)		

## Use Precautions

- 1.Strictly follow the instructions in the manual to use the instrument and train operator.
- 2.Before using the instrument, use a plug with a grounding pin and ensure the power socket is grounded.
- 3.Make sure the voltage of the instrument is adapted. If the local power supply voltage is unstable, we recommend adding a power-matched regulated power supply between the mains and the instrument.
- 4.The microdermabrasion tip must be connected to the suction hole (not the release hole) to prevent the machine from functioning abnormally.
- 5.The rubber ring must be correctly installed to prevent issues like no suction or weak suction.
- 6.The cotton filter is essential to prevent impurities from entering the tube and damaging the machine.
- 7.Keeping the skin moist contributes to a more comfortable experience.
- 8.Do not use the instrument on an empty abdomen, and take the treatment 1 hour later after the meal.
- 9.Always start with the lowest vacuum intensity.

## Maintenance & Cleaning

- 1.Disconnect the power cord when not in use and store the instrument in a dry, cool, and well-ventilated place.
- 2.After use, rinse the microdermabrasion tip in saline, and then clean it with a brush, allow them to air dry, and store them in a dry condition.
- 3.To clean the vacuum cup and grease cup, wipe them with a wet wipe or saline to clean and sterilize, then dry them with a dry towel.

4. After cleaning the instrument, cover it with a dustcloth or plastic wrap. Before use, ensure to sanitize the handle again for hygienic reasons.

5. Disassembling the instrument and attempting to clean the inside is strictly prohibited. This task should only be performed by professionals.

## FAQ

1. Will the skin get thinner and more sensitive after the microdermabrasion treatment?

A: No. The microdermabrasion treatment is a non-invasive procedure that targets only the aged stratum corneum. Since our skin naturally renews itself every 28 days, the microdermabrasion contributes to smoother, brighter, and finer skin without causing increased sensitivity.

2. Why do I need exfoliation?

A: The stratum corneum is the outermost layer of the epidermis, consisting of 15 – 20 layers of dead cells that can vary from different parts of the human body. The skin's metabolism cycle also varies among individuals of different ages. For example, the cycle is 28 days for teenagers, whereas for adults, it is the sum of their age plus 10. However, for a variety of reasons, the stratum corneum might not shed in a timely manner, resulting in an increased thickness of the stratum corneum, a dull complexion, and an impact on skin metabolism. In such cases, additional exfoliation is necessary.

3. Can I apply makeup immediately after a microdermabrasion treatment?

A: Yes, the treatment contributes to achieving moisturized and plump skin, ensuring a seamless and well-blended makeup application. Microdermabrasion is a non-invasive procedure known for providing quick, safe, and healthy results.

4. Are there any side effects associated with vacuum breast enlargement treatment?

A: No side effects at all. Vacuum breast enlargement therapy is entirely based on physical massage, suction, and release under negative pressure. With the integration of multiple advanced technologies, it provides users with immediate effects without any side effects or rebound effects.

5. When do I expect to see the result of facial treatment with the vacuum cup?

A: After treatment, you'll notice purified pores and smoother skin instantly. We recommend a weekly treatment for best results, with increasing improvements day by day.

6. How does negative pressure help in detoxification?

A: Negative pressure can stimulate capillary congestion and boost cellular vitality. The suction-release process, along with alternating negative pressure, continuously opens and closes local pores, improving skin breathing, enhancing oxygen uptake, and accelerating waste removal during treatment. In summary, this treatment promotes Qi and blood circulation, dispels dampness and cold, clears meridians, detoxifies, and purges heat.



## Ask & Answer

Ask	Cause	Solution
1.Fail to turn on the instrument.	→ The power cord is disconnected from the wall socket.	→ Check if the power cord is well connected to the wall socket and ensure that the socket is functional.
	→ The fuse is loose or burned.	→ Replace the fuse on the back of the instrument and restart it.
2.The instrument is not working and can't be restarted.	→ The fuse may have burned out due to unstable voltage or prolonged use.	→ Please replace the fuse and restart the instrument 5 minutes later.
3.The screen shows an error message.	→ Data error.	→ 1. Turn off the instrument, disconnect the power cord, and wait for 1 minute to restart the instrument. 2. If the above method doesn't resolve the problem, please contact your dealer.
4.The dermabrasion tip has a loss of suction.	→ The rubber ring is missing.	→ Put in a rubber ring.
	→ The rubber ring on the tip is excessively twisted, causing air leakage.	→ Loosen it slightly.

## Operation guide

### Face :

#### DERMABRASION tip



Step 1: (3 minutes)



#### Parameters setting

Recommended working time: 2 minutes


Recommended vacuum intensity: 20cmHg

Product needed: Wet cotton pad

#### Vacuum cup for facial cleansing



Step 2: (3 minutes)

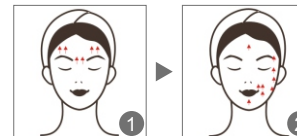

#### Parameters setting

Recommended working time: 3 minutes

Recommended vacuum intensity: 20cmHg

Recommended suction time: 4.5


Recommended release time: 1.5

Product needed: Wet cotton pad

## OXYGEN Spray



Step 3: (10 minutes)

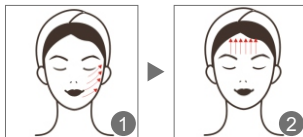

Parameters setting

Recommended working time: 10 minutes

Product needed: Normal saline and serum in a ratio of 10:1.

Spray adjustment: Adjust with the knob on the end of the sprayer.


## Vacuum cup/Grease cup



Step 4: (10 minutes)



Step 5: (10 minutes)

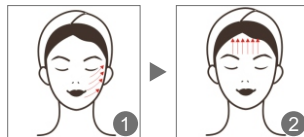

Parameters setting

Recommended working time: 10 minutes

Recommended vacuum intensity: 20cmHg

Recommended Mode: M0


Product needed: Facial essential oil

## Breast enlargement treatment

### Breast cup



Step 1: (20 minutes)

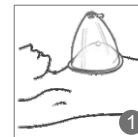

Parameters setting

Recommended working time: 20 minutes.

Recommended vacuum intensity: 20cmHg

Recommended Mode: M0


Product needed: essential oil

## Butt augmentation

### Vacuum cup



Step 1: (15 minutes)

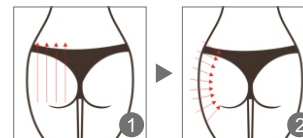

Parameters setting

Recommended working time: 15 minutes

Recommended vacuum intensity: 20cmHg

Recommended suction time: 4.5


Recommended release time: 1.5

Product needed: essential oil

## Breast cup



Step 2: (20 minutes)

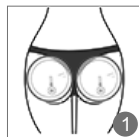

### Parameters setting

Recommended working time: 20 minutes.

Recommended vacuum intensity: 20cmHg

Recommended Mode: M0

Product needed: essential oil

## Arm


### Vacuum cup/Grease cup



Step 1: (15 minutes)



Step 2: (15 minutes)

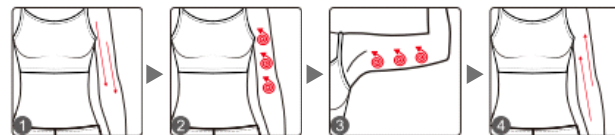

### Parameters setting

Recommended working time: 15 minutes

Recommended vacuum intensity: 20cmHg

Recommended suction time: 4.5

Recommended release time: 1.5

Product needed: essential oil

## BACK


### Vacuum cup/Grease cup



Step 1: (15 minutes)



Step 2: (15 minutes)

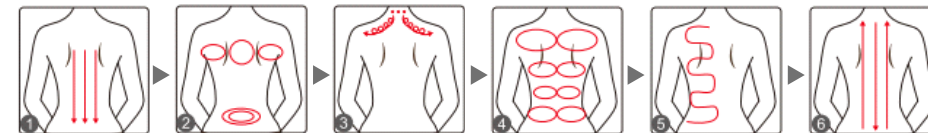

### Parameters setting

Recommended working time: 15 minutes

Recommended vacuum intensity: 20cmHg

Recommended Mode: M0

Product needed: essential oil

## Leg


### Vacuum cup/Grease cup



Step 1: (15 minutes)



Step 2: (15 minutes)

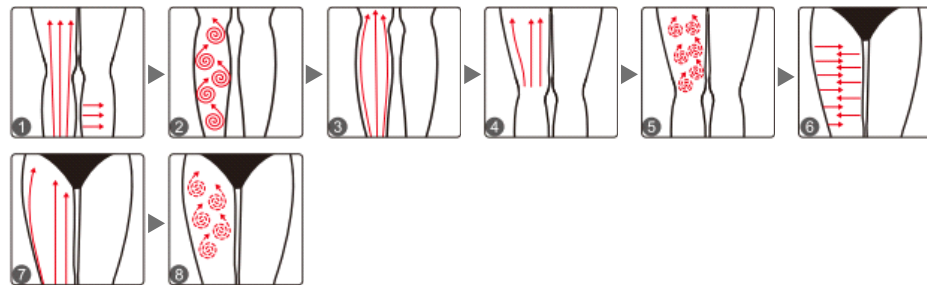

#### Parameters setting

Recommended working time: 15 minutes

Recommended vacuum intensity: 20cmHg

Recommended suction time: 4.5

Recommended release time: 1.5

Product needed: essential oil

## Warranty & Support

We offer a one-year warranty on the main unit and a three-month warranty on components.

Product Name: Versatile Beauty Device

Product Model: VT-3271

Production Date: \_\_\_\_\_

Quality Inspection: \_\_\_\_\_

PASS