


User Manual

MultiLift UltraSkin

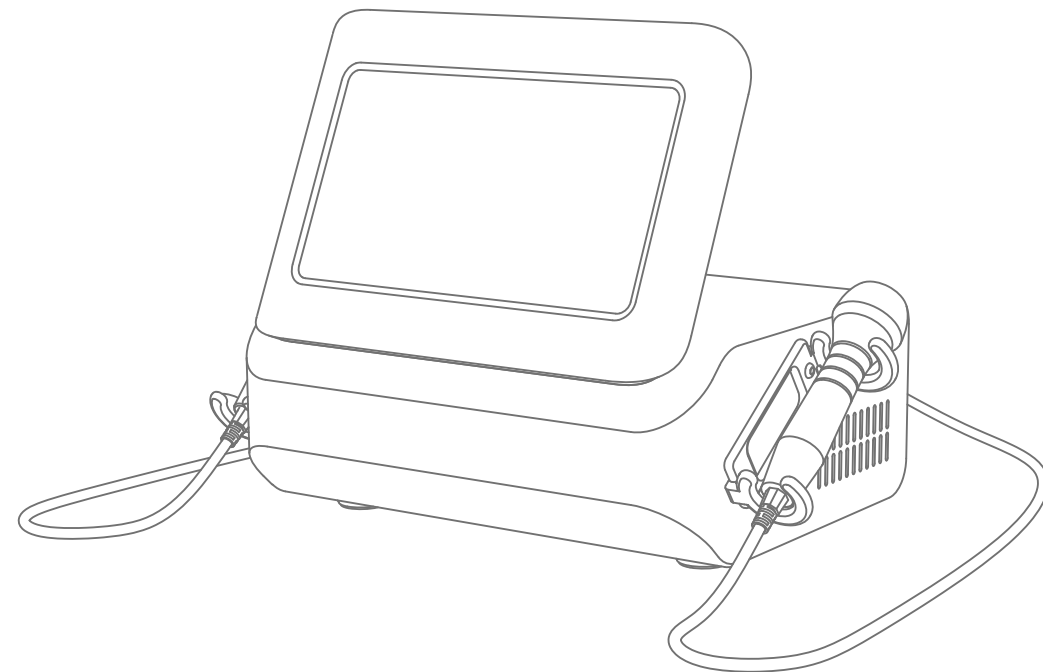

Contents

Part I

1. Product Introduction	1
1.1 Product Overview	1
1.2 Product Components	2
2. Precautions and Safety Instructions	3
2.1 For Installation	3
2.2 For Use	5
2.3 Warning	5
2.4 Contraindications	6
3. Product Features	7
3.1 Product Specifications	7
3.2 Product Configuration	8
3.3 Interface Operation	10
4. How to use	17
4.1 Preparation	17
4.2 Installation	18
4.3 Treatment Application	18
5. Care after treatment	43
6. Skin care error	43
7. FAQ	45
8. Warranty & Support	46


1.Product Introduction

1.1 Product Overview

With the cutting-edge and secure dual-frequency ultrasound technique, this beauty device not only introduces essential nutrients to the skin but also targets connective tissue. This process promotes skin absorption and skin regeneration, resulting in a range of skin enhancements such as improved firmness, a brighter complexion, and effective skin repair. And associated massage relaxes the body and promotes circulation, ensuring a more vibrant skin and figure. It provides the perfect solution for people seeking improvements in skin appearance and body figures without injections, drugs, or side effects. For external use with visible effects.

When compared to traditional single-frequency ultrasound, dual-frequency ultrasound can effectively target multiple skin layers and deliver energy to deeper layers of the skin. The 3MHz/10MHz ultrasound targets skin issues, it covers the epidermal layer and corium layer, efficiently activating and repairing the skin's regenerative capacity. This process significantly alleviates issues such as aging, sagging, wrinkles, sensitivity, and more, resulting in a radiant and youthful skin appearance. The 1MHz/3MHz ultrasound, tailored for body treatments, effectively penetrates from the skin surface to the subcutaneous tissue. Utilizing massage, it triggers collagen production and restructuring, promoting circulation, and resulting in notable body toning and sculpting.

1.2 Product Components

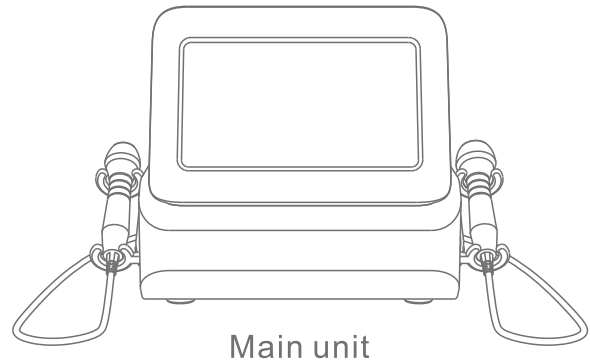

Applicator

Skin care (massage)		
Applicator		
Working area	FACE	BODY

2.Precautions and Safety Instructions

2.1 For Installation

(1) For installation safety

- a.Do not place down the device or strike it, as this may cause internal damage.
- b.Ensure to install and operate the device on a flat and stable surface.
- c.Please install the device in a dry and cool environment. Exposure to moisture or water can lead to insulation degradation.
- d.Strictly follow the instructions in the manual to install the device.
- e.Ensure that the socket is grounded, and add an additional ground connection if the socket is not grounded.
- f.This device only operates at the rated voltage, so please check the voltage before installation.
- g.Always keep the power cord away from water.
- h.Do not use a damaged power cord.
- i.Never share the same power socket with other appliances. As this may cause abnormal heat generation, leading to a fire.

(2) For user safety

- a.This device should only be operated by trained individuals.
- b.This device is not intended for use by pregnant women or children.
- c.Individuals with the following conditions should refrain from using this device:
pacemaker, malignant tumor, hyperpyrexia, a history of heart attack, severe heart disease, pregnancy,

acute disease, skin disease at the treatment area, and artificial internal organs or bones.

d.This device should only be used on the human body.

e.Before operation, do not use cream, foundation, powder, etc. Cleanse your face with a mild facial cleanser after the operation.

f.Before operation, check if there are any dirt or foreign bodies on the treatment head, and clean it with gauze and saline.

g.Always start with the lowest energy level on areas with thinner skin, and gradually increase the energy level as you become accustomed to the treatment.

2.2For Use

a.Please read this manual carefully and follow the instructions in the manual to operate the device.

b.Make sure to familiarize yourself with the application method before use.

c.Before operating, ensure that the treatment head is in good condition. Stop using it and replace it immediately if it is broken.

d.During the treatment, it's important to apply gel or serum to the skin.

e.If the device malfunctions, stop using it immediately and power it off.

f.If there is a power failure, immediately power off the device and reset all the buttons.

g.Do not use the device for any other purpose.

2.2For Use

a.Please read this manual carefully and follow the instructions in the manual to operate the device.

b.Make sure to familiarize yourself with the application method before use.

c.Before operating, ensure that the treatment head is in good condition. Stop using it and replace it immediately if it is broken.

d.During the treatment, it's important to apply gel or serum to the skin.

e.If the device malfunctions, stop using it immediately and power it off.

f.If there is a power failure, immediately power off the device and reset all the buttons.

g.Do not use the device for any other purpose.

2.3Warning

(1)This device should not be used in combination with the following electronic medical devices.

a.Implantable electronic medical devices like a pacemaker should be avoided, as interference may cause malfunctions. If individuals with a pacemaker wish to use this device, they should consult the manufacturer of the pacemaker or a cardiologist and obtain permission before using this device.

b.Life-sustaining electronic systems such as artificial hearts or lungs.

c.Portable or mounted electronic medical devices, such as electrocardiographs.

d.Do not use this device with other electrical stimulators.

(2)General Conditions

- a.It is strictly prohibited to disassemble or modify the device in any way.
- b.Pregnant women should consult their doctor before using this device.
- c.Always keep the manual readily available near the device for reference.
- d.Do not place down the device or strike it.
- e.Individuals with infectious diseases or injuries should avoid contact with this device.
- f.Stop using the device if there is a crack on the applicator, as this may cause skin injuries.

2.4Contraindications

Before using this instrument, it is essential to consult your doctor or seek advice from professionals. Details below:

- 1.Pregnant or lactating woman.
- 2.Individuals with hypertension, heart disease, or a pacemaker installed.
- 3.Individuals with unhealed wounds or those in the recovery period from surgery.
- 4.Individuals suffering from epilepsy, severe diabetes, and hyperthyroidism.
- 5.Individuals suffering from malignant tumors, hemophilia, and hemorrhage.
- 6.Individuals with infectious disease.
- 7.Individuals have had a cesarean section within six months, a miscarriage within three months, or a natural birth within two months.
- 8.People who had plastic surgery, with implants or metal objects at the treatment area in 3 months.

3.Product Features


3.1Product Specifications

Power supply Input: AC110V-260v 50HZ-60Hz	Net Weight: 6KG	
Power: 19W	Gross Weight: 7KG	
Dimensions: 420 × 350 × 280mm		
Input Power: 19W		
1M&3M body	1M Frequency: 1MHz	Power: 15W
	3M Frequency: 3MHz	
3M&10M face	3M Frequency: 3MHz	Power: 19W
	10M Frequency: 10MHz	
Noise level: < 60db (30cm away)		

3.2Product Configuration

Appearance Description


a.Main unit

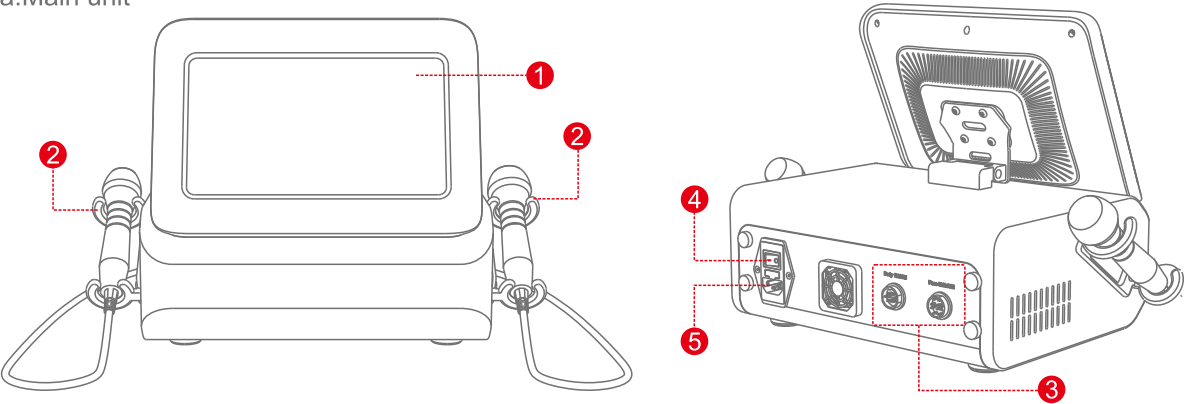

No.	Name	Working area
①	Monitor & Operator Panel	The touchscreen displays the operating interface with various control and display functions, including options for starting and stopping, setting the treatment time, adjusting the energy level, etc.
②	Applicator holder	The place where the applicator is positioned.
③	Socket for the applicator	The left socket is designated for the body applicator, while the right one is intended for the face applicator.


No.	Name	Working area
④	Power button.	To power on/off the device.
⑤	Socket for the power cord	The socket for the power cord.

b. Applicator

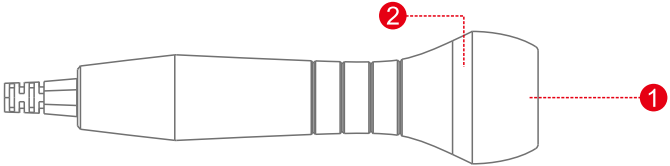


No.	Name	Working area
①	Treatment head	Apply this part of the applicator to the skin during operation.
②	LED	The LED lights blue during operation.

3.3Interface Operation

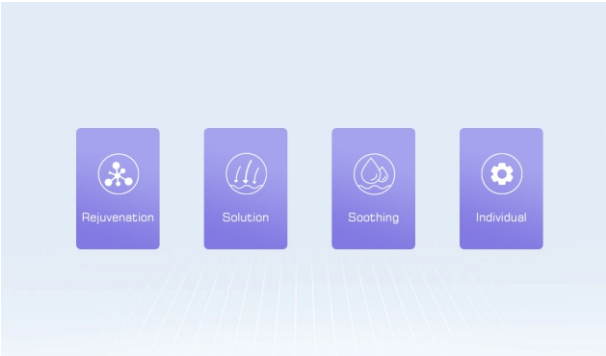

Function Selection Interface (Rejuvenation, Solution, Soothing, Individual)


a.Select Rejuvenation

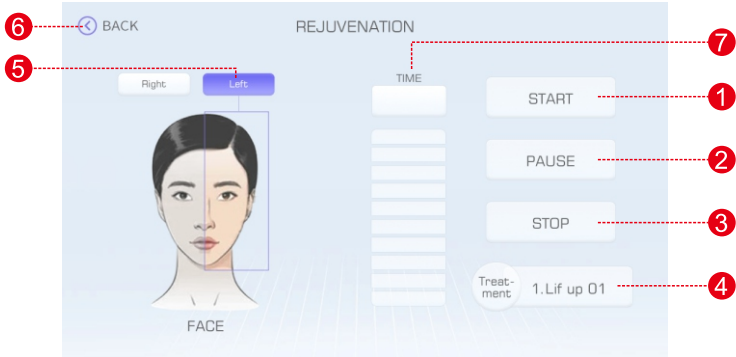

- 1 START button**
After selecting "Rejuvenation," the interface shown above will appear. You can then tap the "START" button to begin the treatment.
- 2 PAUSE button**
Tap this button to pause the treatment in progress, and the timer will also pause. To resume the treatment, tap the "START" button.
- 3 STOP**
Tap this button to terminate the treatment in progress, and the working time will reset to the initial value.
- 4 Treatment**
The default treatment is "1. Lif up 01", tap the "treatment" button and you may switch to other treatments as listed below.



For example, select "2.Lif up 02", and then the treatment is switched from "1. Lif up 01" to "2. Lif up 02". For more details on these treatments, please refer to section 4.3 "Treatment Application."

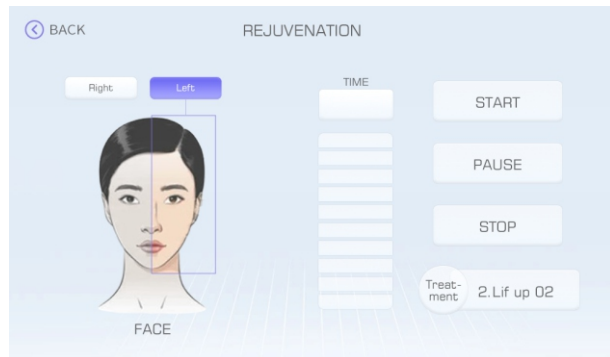

5 Treatment area

The current treatment area is displayed on the screen.

6 Back to the previous page

7 The remaining working time

The timer shows the remaining working time after tapping the "START" button.


b. Select Solution



Tap the "treatment" button and you may select from the treatments listed below. For more details on these treatments, please refer to section 4.3 "Treatment Application."



c. Select Soothing

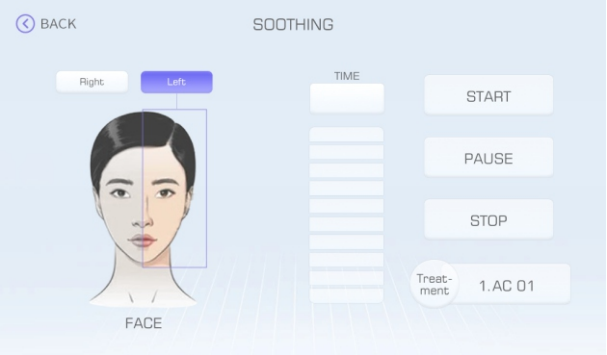


(1) The default treatment is "1. Lif up 01", tap the "treatment" button and you may switch to other treatments as listed below. For more details on these treatments, please refer to section 4.3 "Treatment Application."




d. Select Individual (customized treatment)

Dual frequency, page (1)

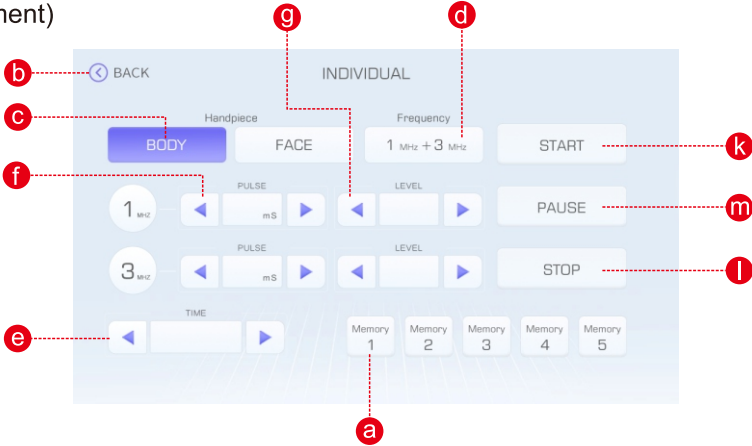

Single frequency, page (2)

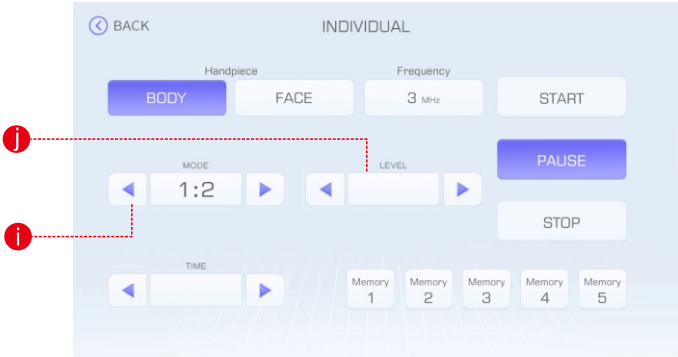

a Parameter Save Button

Long press the button to save the parameters, while a short press allows you to read them.

b Back to the previous page

c Select the applicator

Only one applicator is available, BODY or FACE.

d Frequency set button

Select from the available frequencies listed below.

BODY: 1MHz, 3MHz(Single), 1/3MHz(Dual) FACE: 3MHz, 10MHz(Single), 3/10MHz(Dual)

The screen displays page (1) when you select dual-frequency, and page (2) is displayed when you select single-frequency.

e Working time set button

With each press, the working time will increase/decrease by 1 minute, and the time is displayed in the e area.

f Set the output time for dual-frequency

In dual-frequency mode, you can adjust the output time within the range of 1ms to 20ms, increasing or decreasing it by 1ms with each press.

g Set the output intensity for dual-frequency

In dual-frequency mode, you can set the output intensity within the range of 1 to 10 levels, increasing or decreasing it by 1 level with each press.

i Single-frequency mode set button

Tap the button to switch the mode from 1:2, 1:5, 1:10, and Continuous.

j Set the output intensity for single-frequency

For modes 1:2, 1:5, and 1:10, there are 1 to 10 levels of output intensity available. For mode Continuous, there are 1 to 5 levels available.

k Start button

l Stop button

m Pause button

4. How to use

4.1. Preparation


1) Power supply check

Check if the power socket is grounded.

Check if all the cords are properly connected.

Before turning on the device, make sure the applicators are securely connected to the main unit.

4.2 Installation

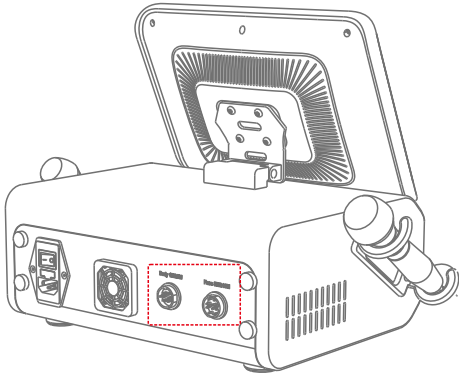

How to connect the applicators to the main unit.

Find the sockets on the back of the main unit, and connect the FACE applicator to the left socket and the BODY applicator to the right socket.

4.3.Treatment Application
A. Treatments of the "Rejuvenation"

1.Lif up 01(Wrinkles)

Duration: 12 minutes

Treatment area: Face

Frequency: 3/10MHz

Applicator: FACE

Recommended treatment sessions: 10 - 12 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator					
Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M					
Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHZ)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	5	3	Left side	Applicator 2
3/10MHz	3ms/7ms	3MHz—5 10MHz-5	3	Left side	Applicator 2
10MHz	Cont	5	3	Right side	Applicator 2
3/10MHz	3ms/7ms	3MHz—5 10MHz-5	3	Right side	Applicator 2

Application: Begin the treatment in the area indicated on the screen, and switch to the other side after 6 minutes when you hear a beep.

Note: This treatment is designed for crow's feet. It not only improves complexion but also alleviates deep wrinkles. The operating area should be larger than the treatment area, and you can use retinoid products simultaneously.

2.Lif up 02 (Folds)

Duration: 22 minutes

Treatment area: Face

Frequency: 1/3/10MHz

Applicator: FACE/BODY

Recommended treatment sessions: 10 - 12 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator					
Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M					
Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	5	3	Left side	Applicator 2
3/10MHz	3ms/7ms	3MHz—5 10MHz-5	3	Left side	Applicator 2
10MHz	Cont	5	3	Right side	Applicator 2
3/10MHz	3ms/7ms	3MHz—5 10MHz-5	3	Right side	Applicator 2
1/3MHz	3ms/7ms	1MHz—5 3MHz-5	5	Left side	Applicator 1
1/3MHz	3ms/7ms	1MHz—5 3MHz-5	5	Right side	Applicator 1

Application: Start the treatment on the left side of your face with the applicator 2(3/10MHz). Switch to the right side after 6 minutes when you hear a beep. Use the applicator 1(1/3MHz) on the left side of

your face after another 6 minutes when you hear a beep. Switch to the right side of your face after 5 minutes when you hear a beep. (Left-Right-Left-Right)

Note:

The operating area should be larger than the treatment area. You can prolong the treatment time on areas with more pronounced wrinkles, such as folds. Avoid application on the bony area, and you can use retinoid products simultaneously.

3.Lif up 03 (Connective tissue lift)

Duration: 24 minutes

Treatment area: Face

Frequency: 1/3/10MHz

Applicator: FACE/BODY

Recommended treatment sessions: 10 - 12 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHZ)	Mode	Levels	Time (min)	Application area	Applicator
3/10MHz	1ms/1ms	3MHz—5 10MHz-5	3	Left side	Applicator 2
3/10MHz	1ms/1ms	3MHz—5 10MHz-5	3	Left side	Applicator 2
3/10MHz	1ms/1ms	3MHz—5 10MHz-5	3	Right side	Applicator 2
3/10MHz	1ms/1ms	3MHz—5 10MHz-5	3	Right side	Applicator 2
1/3MHz	1ms/1ms	1MHz—4 3MHz-4	3	Left side	Applicator 1
1/3MHz	1ms/1ms	1MHz—4 3MHz-4	3	Left side	Applicator 1
1/3MHz	1ms/1ms	1MHz—4 3MHz-4	3	Right side	Applicator 1
1/3MHz	1ms/1ms	1MHz—4 3MHz-4	3	Right side	Applicator 1

Application:
Start the treatment on the left side of your face with the applicator 2 (3/10MHz), and switch to the right side after 6 minutes when you hear a beep. After another 6 minutes, switch to the applicator 1 (1/3MHz) and apply it on the left side of your face. And switch to the right side in the 19th minute when you hear a beep (Left-Right-Left-Right).

- Note:
- a.Evenly apply the applicator to the entire face.
 - b.Your face may experience slight edema and redness after the treatment.

c.When treating bony areas like the forehead, it's important to reduce the energy level and increase the number of applications.

d.For individuals who smoke or are frequently exposed to sunlight, the treatment may be less effective compared to others. Therefore, it's recommended that they undergo more frequent treatments and additional sessions. A customized treatment plan is advisable for better results.

4.Tight 01 (Skin tightening)

Duration: 12 minutes
Treatment area: Face
Frequency: 3/10MHz
Applicator: FACE
Recommended treatment sessions: 10 - 12 sessions
Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	4	3	Left side	Applicator 2
3/10MHz	3ms/7ms	3MHz—4 10MHz-5	2	Left side	Applicator 2
3/10MHz	1ms/1ms	3MHz—4 10MHz-5	1	Left side	Applicator 2

Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	1.5	3	Right side	Applicator 2
3/10MHz	3ms/7ms	3MHz—4 10MHz-5	2	Right side	Applicator 2
3/10MHz	1ms/1ms	3MHz—5 10MHz-5	1	Right side	Applicator 2

Application: Switch to the right side of your face after 6 minutes when you hear a beep. (Left-Right)

Note: The operating area should be larger than the treatment area. During one session, continue the application until the face appears soft and thermal puffiness. A maintenance treatment is recommended every 2 to 3 months for optimal results.

5.Tight 02 (Skin tightening body)

Duration: 22 minutes

Treatment area: Body

Frequency: 1/3/10MHz

Applicator: FACE/BODY

Recommended treatment sessions: 10 - 12 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator					
Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M					
Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHZ)	Mode	Levels	Time (min)	Application area	Applicator
3/10MHz	3ms/7ms	3MHz—4 10MHz-4	5	Left side	Applicator 2
3/10MHz	1ms/1ms	3MHz—4 10MHz-4	1	Left side	Applicator 2
3/10MHz	3ms/7ms	3MHz—4 10MHz-4	5	Right side	Applicator 2
3/10MHz	1ms/1ms	3MHz—4 10MHz-4	1	Right side	Applicator 2
1/3MHz	1ms/9ms	1MHz—4 3MHz-4	3	Left side	Applicator 1
1/3MHz	1ms/1ms	1MHz—4 3MHz-4	2	Left side	Applicator 1
1/3MHz	1ms/9ms	1MHz—4 3MHz-4	3	Right side	Applicator 1
1/3MHz	1ms/1ms	1MHz—4 3MHz-4	2	Right side	Applicator 1

Application: Start the treatment on the left side of your face with the applicator 2(3/10MHz). Switch to the right side after 6 minutes when you hear a beep. Use the applicator 1(1/3MHz) on the left side of your face after another 6 minutes when you hear a beep. Switch to the right side after 5 minutes when you hear a beep. (Left-Right-Left-Right)

Note: During one session, continue the application to ensure the repeated appearance of soft and thermal puffiness on the face. Take a treatment every 2-3 months as needed.

6.Tight 03(Skin rejuvenation)

Duration: 12 minutes

Treatment area: Face

Frequency: 3/10MHz

Applicator: FACE

Recommended treatment sessions: 10 - 12 sessions

Recommended treatment frequency: 2 sessions per week

Applicator					
Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M					
Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	2	3	Left side	Applicator 2
3/10MHz	3ms/7ms	3MHz—2 10MHz-2	3	Left side	Applicator 2
10MHz	Cont	2	3	Right side	Applicator 2
3/10MHz	3ms/7ms	3MHz—2 10MHz-2	3	Right side	Applicator 2

Application: Switch to the right side after 6 minutes when you hear a beep. (Left-Right)

a.This treatment is designed to address skin aging issues, such as skin laxity.

b.Evenly apply the applicator to the entire face.

c.Your face may experience slight edema and redness after the treatment.

d.When treating bony areas like the forehead, it's important to reduce the energy level and increase the number of applications.

e.For individuals who smoke or are frequently exposed to sunlight, the treatment may be less effective compared to others. Therefore, it's recommended that they undergo more frequent treatments and additional sessions. A customized treatment plan is advisable for better results.)

f.It is recommended to take vitamins A, C, and E for optimal results.

g.This treatment can be combined with peeling, tightening (RF), or skin rejuvenation procedures.

7.Rejuve 01(Anti-aging)

Duration: 12 minutes

Treatment area: face, neck

Frequency: 3/10MHz

Applicator: FACE

Recommended treatment sessions: 10 - 12 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	2	3	Left side	Applicator 2
3/10MHz	1ms/1ms	3MHz—2 10MHz-2	3	Left side	Applicator 2
10MHz	Cont	2	3	Right side	Applicator 2
3/10MHz	1ms/1ms	3MHz—2 10MHz-2	3	Right side	Applicator 2

Application: Switch to the right side of your face after 6 minutes when you hear a beep. (Left-Right)

Note:

- a.The operating area should be larger than the treatment area and the thermal puffiness should be avoided.
- b.Take a treatment every 2-4 months as needed, and avoid sun exposure during the 2nd and 4th sessions.

8.Rejuve 02(Stretch mark)

Duration: 22 minutes

Treatment area: Body

Frequency: 1/3/10MHz

Applicator: FACE/BODY

Recommended treatment sessions: 10 - 12 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	2	3	Left side	Applicator 2
3/10MHz	1ms/9ms	3MHz—2 10MHz-2	3	Left side	Applicator 2
10MHz	Cont	2	3	Right side	Applicator 2
3/10MHz	1ms/9ms	3MHz—2 10MHz-2	3	Right side	Applicator 2
1/3MHz	1ms/1ms	1MHz—2 3MHz-4	2	Left side	Applicator 1
1/3MHz	3ms/7ms	1MHz—2 3MHz-4	3	Left side	Applicator 1
1/3MHz	1ms/1ms	1MHz—2 3MHz-4	2	Right side	Applicator 1
1/3MHz	3ms/7ms	1MHz—2 3MHz-5	3	Right side	Applicator 1

Application: Switch to the right side of your face after 6 minutes when you hear a beep. When you hear the 2nd beep after another 6 minutes, switch to applicator 1 (1/3MHz) and apply it on your left side. After 5 minutes, switch to the right side when you hear the 3rd beep. (Left-Right-Left-Right)

Note: The operating area should only cover the treatment area. During one session keep repeating the application. When combined with other treatments such as MTS and laser, it can help achieve improved results.

9.Rejuve 03(Cellulite)

Duration: 30 minutes

Treatment area: Body

Frequency: 1/3MHz

Applicator: BODY

Recommended treatment sessions: 10 - 12 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
1MHz	Cont	5	5	Left side	Applicator 1
1/3MHz	1ms/3ms	1MHz—5 3MHz-7	5	Left side	Applicator 1
1/3MHz	1ms/1ms	1MHz—7 3MHz-7	5	Left side	Applicator 1
1MHz	Cont	5	5	Right side	Applicator 1
1/3MHz	1ms/3ms	1MHz—5 3MHz-7	5	Right side	Applicator 1
1/3MHz	1ms/1ms	1MHz—7 3MHz-7	5	Right side	Applicator 1

Application: Switch to the right side of the body after 15 minutes when you hear a beep (Left-Right).

Note:

a.The operating area should be larger than the treatment area. During one session keep repeating the application.

b.Take a treatment every 2-3 months as needed.

B.Treatments of the "SOLUTION"

Sono-P F(Sonophoresis face)

Duration: 12 minutes

Treatment area: Face

Frequency: 3/10MHz

Applicator: FACE

Recommended treatment sessions: As required

Recommended treatment frequency: As required

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHZ)	Mode	Levels	Time (min)	Application area	Applicator
3/10MHz	3ms/7ms	3MHz—4 10MHz-4	3	Left side	Applicator 2
3/10MHz	1ms/1ms	3MHz—4 10MHz-4	3	Left side	Applicator 2
3/10MHz	3ms/7ms	3MHz—4 10MHz-4	3	Right side	Applicator 2
3/10MHz	1ms/1ms	3MHz—4 10MHz-4	3	Right side	Applicator 2

Application: Switch to the right side of your face after 6 minutes when you hear a beep (Left-Right).

Note:

- a.The operating area should be larger than the treatment area.
- b.This beauty treatment can be applied on both the face and neck.
- c.Medications should be selected by a qualified medical professional to avoid unpredictable results.

2.Sono-P B(Sonophoresis BODY)

Duration: 20 minutes.

Treatment area: Body

Frequency: 1/3MHz

Applicator: BODY

Recommended treatment sessions: As required

Recommended treatment frequency: As required

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHZ)	Mode	Levels	Time (min)	Application area	Applicator
1/3MHz	7ms/3ms	1MHz—3 3MHz-5	10	Left side	Applicator 1
1/3MHz	1ms/1ms	1MHz—4 3MHz-4	10	Left side	Applicator 1

Application: The customized treatment set in the "Individual" mode is recommended for addressing larger treatment areas.

Note: Evenly apply the applicator to the treatment area (Medications should be selected by a qualified medical professional to avoid unpredictable results.)

C.Treatments of the "SOOTHING"

1.AC 01(Acne comedonica)

Duration: 12 minutes

Treatment area: Face

Frequency: 3/10MHz

Applicator: FACE

Recommended treatment sessions: 10 - 12 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	2	3	Left side	Applicator 2
3/10MHz	3ms/7ms	3MHz—2 10MHz-4	3	Left side	Applicator 2
10MHz	Cont	2	3	Right side	Applicator 2
3/10MHz	3ms/7ms	3MHz—2 10MHz-4	3	Right side	Applicator 2

Application: Switch to the right side of your face after 6 minutes when you hear a beep (Left-Right).

Note: A maintenance treatment is recommended every 1 to 3 months for optimal results.

2.AC 02(Acne Papulopustulosa)

Duration: 22 minutes

Treatment area: face, back

Frequency: 1/3/10MHz

Applicator: FACE/BODY

Recommended treatment sessions: 10 - 12 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	2	3	Left side	Applicator 2
3/10MHz	3ms/7ms	3MHz—2 10MHz-4	3	Left side	Applicator 2
10MHz	Cont	2	3	Right side	Applicator 2
3/10MHz	3ms/7ms	3MHz—2 10MHz-4	3	Right side	Applicator 2
1/3MHz	4ms/6ms	1MHz—2 3MHz-4	5	Left side	Applicator 1
1/3MHz	4ms/6ms	1MHz—2 3MHz-4	5	Right side	Applicator 1

Application: Start the treatment with the FACE applicator. Switch to the right side of your face after 6 minutes when you hear a beep. Use the BODY applicator on the left side after another 6 minutes when you hear a beep. Switch to the right side of your face after 5 minutes when you hear a beep. (Left-Right-Left-Right)

Note: A maintenance treatment is recommended every 1 to 3 months for optimal results. This treatment should only be taken after squeezing out the acne.

3.Sensitive S 01 (Atopic dermatitis)
Duration: 7 minutes
Treatment area: Face
Frequency: 3/10MHz
Applicator: FACE
Recommended treatment sessions: 4-8 sessions
Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	2	2	Treatment area	Applicator 2
10MHz	Cont	2	2	Treatment area	Applicator 2
3/10MHz	3ms/7ms	3MHz—2 10MHz-2	3	Treatment area	Applicator 2

Application: For treatment on larger areas of the body, repeat the treatment several times or customize the treatment in the "Individual" setting.
Note: Take a maintenance treatment as required.
Consistent treatments lead to improved results.

4.Sensitive S 02(Perioral Dermatitis)
Duration: 6 minutes
Treatment area: Face
Frequency: 10MHz
Applicator: FACE
Recommended treatment sessions: 4-8 sessions
Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	2	3	Treatment area	Applicator 2
10MHz	Cont	2	3	Treatment area	Applicator 2

Application: For treatment on larger areas of the body, repeat the treatment several times or customize the treatment in the "Individual" setting.
Note: Take a maintenance treatment as required.

5.Sensitive S 03(Eczemas)

Duration: 6 minutes

Treatment area: face, body

Frequency: 3/10MHz

Applicator: FACE

Recommended treatment sessions: 4-8 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	2	3	Treatment area	Applicator 2
3/10MHz	3ms/7ms	3MHz—2 10MHz-2	3	Treatment area	Applicator 2

Application: For treatment on larger areas of the body, repeat the treatment several times or customize the treatment in the "Individual" setting.

Note: Take a maintenance treatment as required.

6.Scar 01 (Hypertrophic scars/keloids)

Duration: 7 minutes

Treatment area: Scar area

Frequency: 3/10MHz

Applicator: FACE

Recommended treatment sessions: 6-8 sessions

Recommended treatment frequency: 1 session per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	5	2	Treatment area	Applicator 2
10MHz	Cont	4	2	Treatment area	Applicator 2
3/10MHz	3ms/7ms	3MHz—4 10MHz-7	3	Treatment area	Applicator 2

Application: For treatment on larger areas of the body, repeat the treatment several times or customize the treatment in the "Individual" setting.

Note: Take a maintenance treatment as required (The results of the treatment will become noticeable after 6 months.).

7.Scar 02 (Acne Scars)

Duration: 14 minutes

Treatment area: face, back

Frequency: 3/10MHz

Applicator: FACE

Recommended treatment sessions: 6-8 sessions

Recommended treatment frequency: Once a week or once every two weeks.

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	4	2	Left side	Applicator 2
10MHz	Cont	4	2	Left side	Applicator 2
3/10MHz	3ms/7ms	3MHz—4 10MHz-4	3	Left side	Applicator 2
10MHz	Cont	4	2	Right side	Applicator 2
10MHz	Cont	4	2	Right side	Applicator 2
3/10MHz	3ms/7ms	3MHz—4 10MHz-4	3	Right side	Applicator 2

Application: Switch to the right side of your face after 7 minutes when you hear a beep. (Left-Right)

Note: Take a maintenance treatment as required. When combined with other treatments such as MTS and RF Microneedling, this can help achieve improved results.

8.Scar 03(Chronic leg ulcer/ulcuscruis)

Duration: 6 minutes

Treatment area: Applicable area

Frequency: 10MHz

Applicator: FACE

Recommended treatment sessions: 4-10 sessions

Recommended treatment frequency: 1 - 2 sessions per week

Applicator Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHZ)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	2	6	Treatment area	Applicator 2

Application: This treatment area may cover the entire body, so for larger areas with scars or multiple lesion sites, customize the treatment using the "Individual" setting.

Note: Take a maintenance treatment anytime as required (It is recommended to use the protection film for severe scars, please contact the dealer for more details).

9.Redness 01(Lupus erythematosus)

Duration: 12 minutes

Treatment area: Face

Frequency: 10MHz

Applicator: FACE

Recommended treatment sessions: 8 - 10 sessions

Recommended treatment frequency: Once a week or once every two weeks.

Applicator					
Applicator 1 (BODY applicator) 1M/3M Ultrasonic applicator - 1M& 3M					
Applicator 2 (FACE applicator) 3M/10M Ultrasonic applicator - 3M& 10M					
Frequency (MHz)	Mode	Levels	Time (min)	Application area	Applicator
10MHz	Cont	2	6	Left side	Applicator 2
10MHz	Cont	2	6	Right side	Applicator 2

Application: Switch to the right side of your face after 6 minutes when you hear a beep (Left-Right).

Note: Take a maintenance treatment anytime as required.

5.Care after treatment

Check if the treatment areas show signs of swelling, redness, blisters, etc.


Cleanse the treatment area with cold water or a mild facial cleanser.


Avoid irritating the skin with makeup or other products, as the skin is sensitive after treatment.

Avoid direct sunlight on the treated area. It is recommended to use sunscreen with an SPF of 30 or higher.

6.Skin care error

Error	
Details	Abnormally high temperature in the FACE/BODY applicator's treatment head.
Measure	Pause the treatment for 5 minutes. If the treatment head cools down during this time, you can then resume the treatment.

Error	
Details	The FACE/BODY applicator is not connected to the main unit properly.
Measure	Turn off the device and then reconnect the FACE/BODY applicator to the main unit (as the main unit may fail to detect the applicator when it's in operation)

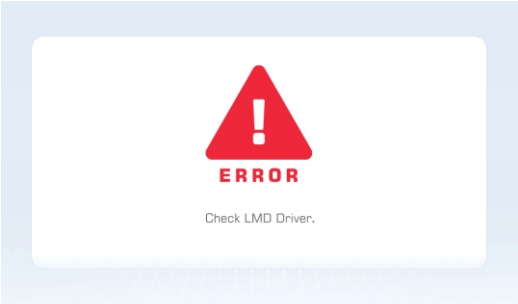Error	
-------	--

Details	Invalid data
Measure	Restart the device, or consult the dealer.

7.FAQ

1.Any skincare precautions after the treatment?

A: It is advisable to employ a facial mask for a span of three uninterrupted days while incorporating a night cream to achieve the most favorable outcomes. Additionally, pay attention to moisturization and sun protection.

2.Can I apply makeup after the treatment?

A: You may apply makeup after treatment if your skin is free from inflammation. However, if inflammation or skin irritation is present, it is advisable to refrain from applying makeup until your skin has stabilized.

3.Why does my skin become dry after the treatment?

A: This treatment elevates the hyaluronic acid content, which functions to retain moisture. Skin dryness is a positive sign of the treatment's effectiveness. It is advisable to apply a mask for three consecutive days to alleviate this effect.

8.Warranty & Support

We provide a one-year warranty and a three-month return and replacement policy from the date of purchase.

Product
Name:

MultiLift UltraSkin

Product
Model:

FC-23K1

Production
Date:

Quality
Inspection:

PASS