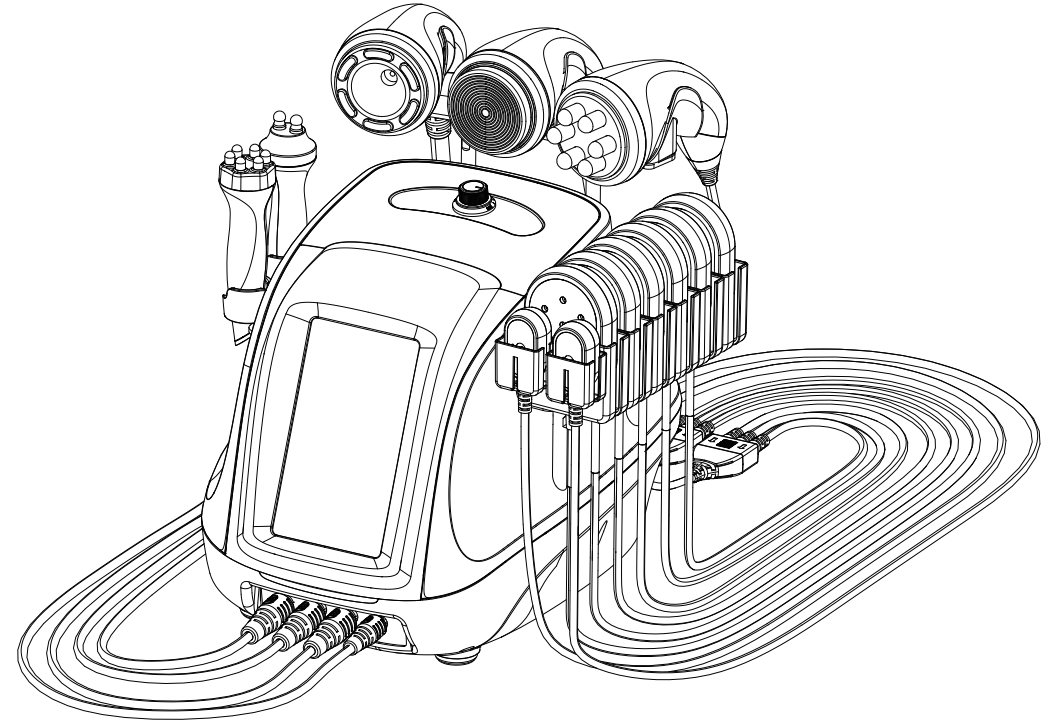


User Manual

LS-78D2SB



Content

1.Safety precautions	1
2.Contraindications	1
3.Overview	2
4.Interface function	3
5.Technical specification	8
6.Precautions	9
7.Maintenance & Cleaning	10
8.FAQ	10
9.Ask & Answer	14
10.Operation guide	16
11.Warranty and support	19

1.Warranty and support

- 1.Please read this manual carefully before use and familiarize yourself with the operation of this instrument and related technology.
- 2.This instrument is a high-precision electronic device that requires proper storage. Please avoid storing it in humid conditions, direct sunlight, or near flammables and explosives to prevent damage.
- 3.Always use the designated accessory provided or recommended by the manufacturer to ensure performance throughout its service life.
- 4.Do not place the instrument in high-temperature conditions, as this may negatively impact its service life and performance.
- 5.Do not apply the instrument on eyes, thyroid, parathyroid, testicles, pregnant women's abdomen, pacemakers, etc.
- 6.Patients should refrain from using this instrument unless they have obtained permission from their doctor.
- 7.Ensure to apply the instrument on moist skin and avoid the operation on dry skin.
- 8.Power off the instrument when not in use to ensure its safety.
- 9.If you are currently taking weight-loss drugs, please discontinue their use 1 - 2 months before using this instrument. Failure to do so may result in an extended course of treatment.
- 10.To maintain hygiene and extend the instrument's service life, sterilize it with saline.

2.Contraindications

Before using this instrument, it is essential to consult your doctor or seek advice from professionals.

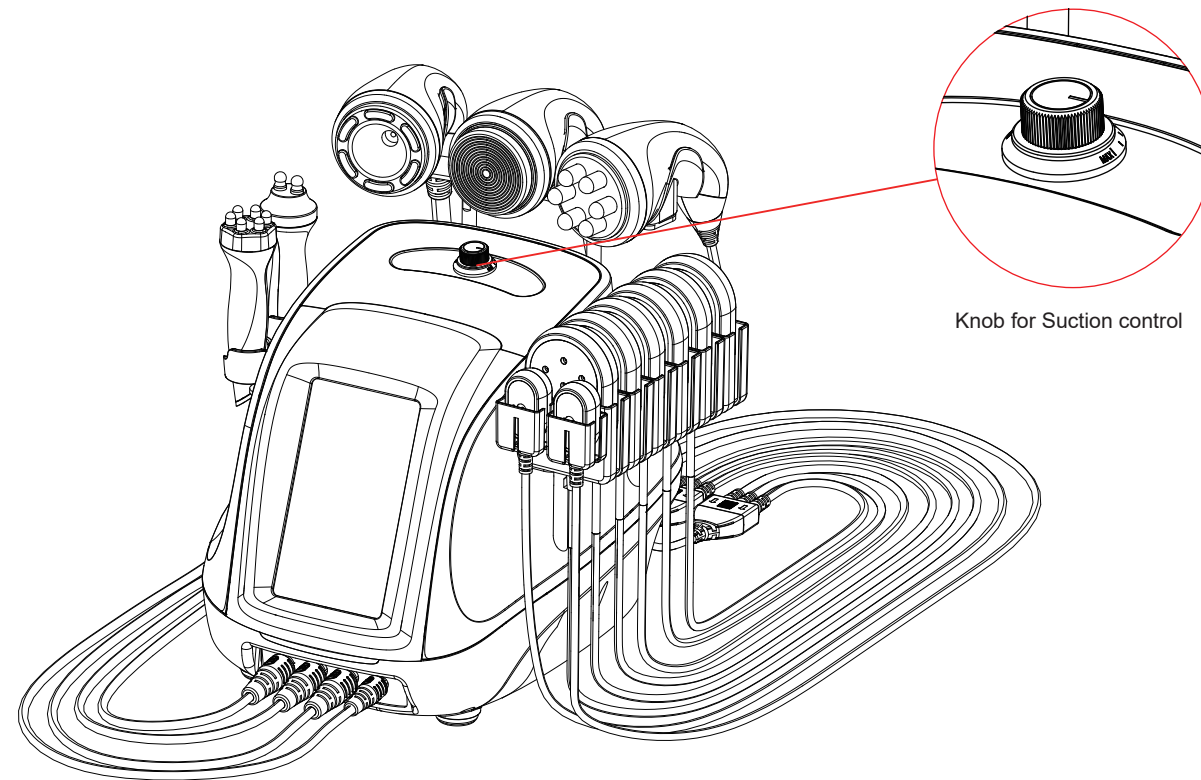
Details as below:

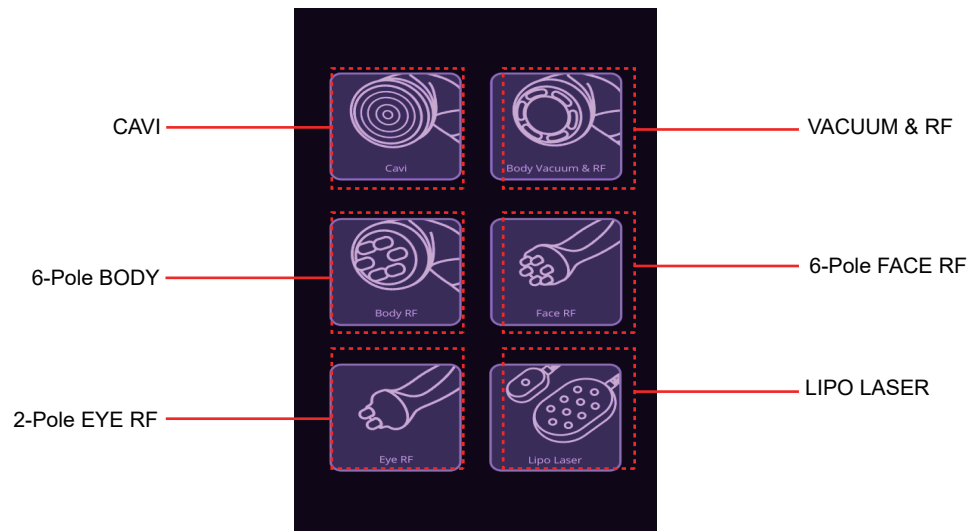
1. Pregnant or lactating woman.
2. Individuals with hypertension, heart disease, or a pacemaker.
3. Individuals with unhealed wounds or those in the recovery period from surgery.
4. Individuals suffering from epilepsy, severe diabetes, and hyperthyroidism.
5. Individuals suffering from malignant tumors, hemophilia, and hemorrhage.
6. Individuals with skin disease or infectious disease.
7. Individuals have had a cesarean section within six months, a miscarriage within three months, or a natural birth within two months.
8. For your safety, refrain from using if you've had plastic surgery or have implants/metal objects near the treatment area within 3 months.

3. Overview

With the most popular and safest technology for anti-aging, firming skin, and body shaping, this 6-in-1 beauty instrument ensures you achieve your desired body figure. It provides the perfect solution for people seeking improvements in skin appearance and body figure without injections, drugs, or side effects. For external use with immediate effects.

4. Interface function






Time Setting (0-60 mins)

Energy Setting (0-10 levels)

Mode Setting
 ∞: Continuous working
 ∩: Intermittent working

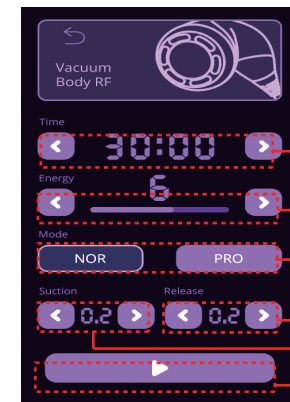Start/Stop



Time Setting (0-60 mins)

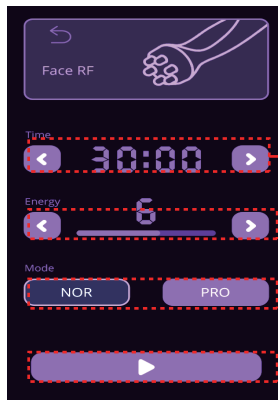
Energy Setting (0-10 levels)

Mode Setting
 NOR: The handle heats up directly upon contact with the skin.
 PRO: The handle gradually heats up upon contact with the skin.

Release Time Setting

Suction Time Setting

Start/Stop




Time Setting (0-60 mins)

Energy Setting (0-10 levels)

Mode Setting

NOR: The handle heats up directly upon contact with the skin.
PRO: The handle gradually heats up upon contact with the skin.

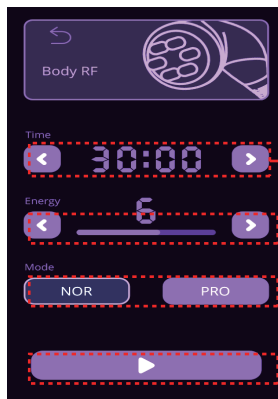Start/Stop



Time Setting (0-60 mins)

Energy Setting (0-10 levels)

Mode Setting

NOR: The handle heats up directly upon contact with the skin.
PRO: The handle gradually heats up upon contact with the skin.

Start/Stop



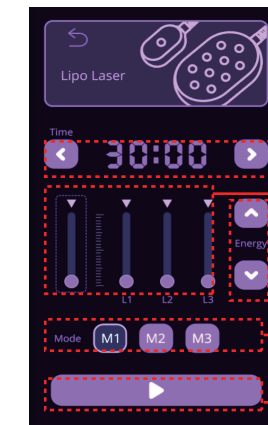
Time Setting (0-60 mins)

Energy Setting (0-10 levels)

Mode Setting

NOR: The handle heats up directly upon contact with the skin.
PRO: The handle gradually heats up upon contact with the skin.

Start/Stop



Time Setting (0-60 mins)

Display of laser energy level

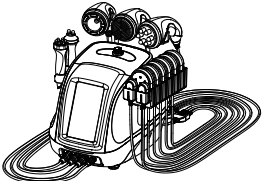
Energy Setting (0-10 levels)




Mode Setting


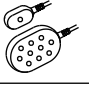
M1: Continuous working
M2: Flash at a medium frequency
M3: Flash at a fast frequency


Start/Stop

5. Technical specification

	Power supply Input: AC110V-220v 50-60Hz		Net Weight: 6.70KG	
	Power: 75W		Gross Weight: 7.74KG	
	Dimensions: 43.5 × 32.5 × 31.5 cm (14.96 × 15.75 × 11.81 cm)			
	Input Power: 75W			

	Frequency: 40KHz		Frequency: 1MHz		Frequency: 1MHz
	Power: 30W		Power: 65W		Power: 55W

	Frequency: 1 MHz		Wavelength: 650nm5mw	
	Power: 32W		Power: 15W	

	Frequency: 1MHz	Air flow: 10LPM	
	Power: 75 W	Noise level: <65db (30cm away)	
	Vacuum: -73KPa		

6. Precautions

1. Strictly follow the instructions in the manual to use the instrument and train operator.
2. Before using the instrument, you must use a plug with a grounding pin and ensure that the power socket is grounded.
3. Make sure the voltage of the instrument is adapted. If the local power supply voltage is unstable, we advise you to add a power-matched regulated power supply between the mains and the instrument.
4. Remove all the metal objects inside the body to avoid causing accidents or affecting treatment effects.
5. Stop the intake of weight drugs 1 - 2 months prior to treatment, or extend the treatment course if you want to take the treatment immediately.
6. Do not use the instrument on an empty abdomen, and take the treatment 1 hour later after the meal.
7. Ensure full skin contact with the treatment head throughout the session to avoid uneven heat distribution.
8. Refrain from using the CAVI handle on areas such as the head, chest, breast, heart, or back.
9. Always start with the lowest vacuum intensity.
10. Ensure the filter cotton is installed inside the vacuum handle, as the foreign body may block the handle.
11. Fix the laser pads onto the client's targeted area, then power on the device. Gradually adjust the energy level from low to high to ensure the client's comfort and warmth.
12. Maintain complete contact between the laser pad surface and skin, avoiding partial contact.
13. Combining the laser pad with the massage or other handles enhances the treatment's effectiveness for optimal results. You can also use the laser pads alone. Simply fix them in the treatment area and wait for 20 minutes.

7. Maintenance & Cleaning

1. Disconnect the power cord when not in use and store the instrument in a dry, cool, and well-ventilated place.
2. Clean and sanitize the handle's surface with a damp cloth or saline, and then dry it thoroughly with a dry towel to ensure it remains dry.
3. Before cleaning the Vacuum & RF handle, ensure to remove the filter cotton to prevent the growth of bacteria.
4. After cleaning the instrument, cover it with a dustcloth or plastic wrap. Before use, ensure to sanitize the handle again for hygienic reasons.
5. Disassembling the instrument and attempting to clean the inside is strictly prohibited. This task should only be performed by professionals.

8. FAQ

1. How long does it take to see results?

A: You might experience skin tightness right after the treatment, as the heat stimulation tightens the collagen tissue, and the RF stimulates the subcutaneous tissue the brings the continuous regeneration of collagen, leading to more noticeable results within one week.

2. Is RF harmful to the skin?

A: RF firming and anti-wrinkle treatment is a non-surgical method that stimulates collagen regeneration

and metabolism in the deeper layers of the skin without causing any harm. It's normal to experience temporary heat and redness in the treated area, which indicates increased blood circulation and should fade within a few hours. There's no need to be concerned about this reaction.

3. What is the duration of a body sculpting treatment?

A: Approximately 60 minutes, our skilled techniques and advanced instrument guarantee impressive results.

4. Liposuction or our device?

A: Liposuction surgery utilizes negative pressure to remove excess fat from specific areas, resulting in a slimming effect. It provides you with immediate results, but it also carries risks, requires anesthesia, and involves downtime. While our device offers a risk-free solution with no side effects. It sculpts your body, stimulates collagen regeneration for lifted skin, and creates the ideal figure for an overall improved appearance.

5. Is it safe to use ultrasound on the entire body?

A: NO. Ultrasound can break down fat cells, causing the cell walls to rupture and releasing fat, which is then absorbed and metabolized by the lymphatic system. The power of ultrasound is quite strong, and the heart is particularly sensitive to sound waves due to being a hollow organ. Sound waves can cause vibrations in the heart, and there can be variations in how sound waves are conducted through muscle tissue and blood within the heart, resulting in back-and-forth reflections. As a consequence, the valve may become dislodged from the heart muscle. Certainly, ultrasound treatment should be strictly avoided around the eye area to prevent the risk of retinal detachment, as well as areas such as the waist, back, and chest, as well as areas such as the waist, back, and chest.

6. Are there any side effects of ultrasound on the human body?

A: High-frequency sound wave is a non-surgical, non-invasive procedure that does not require incision or anesthesia. How does the ultrasound work: ultrasound (mechanical wave) - cavitation blasting - ultrasound focussing. The effect: fat breaking - fat dissolving - skin tightening and body sculpting, therefore, it only targets the low-density fat tissue and protects the high-density tissue such as blood vessels and nerves, ensuring no adverse effects on the human body. You may notice a slight noise in your ears during the session, which is a normal occurrence and requires no cause for concern.

7. Why do I get tinnitus during the treatment?

A: When the ultrasound with concentrated energy vibrates strongly and has a frequency higher than 20Khz, it can penetrate the skin to reach the fat layer that is 20mm deep. The energy causes high-speed friction and heat among the fat cells, which breaks them down and emulsifies them. However, this process can create noise that may be audible in your ears.

8. What is collagen?

A: Collagen is a biological polymer substance, a white opaque fibrous protein without branched chains, which can replenish essential nutrients needed by all layers of the skin and enhance the collagen activity in the skin, with the effects of locking moisture, moisturizing the skin, delaying aging, improving facial laxity and hair care. Collagen is an essential nutrient for promoting anti-aging effects. As you age, the gradual loss of collagen in your skin becomes evident. Typically, a woman's skin starts showing signs of aging around the age of 20 due to collagen depletion, and the rate of loss peaks at 25. Research suggests that by the age of 40, collagen levels are only about half of what they were at 18. The reduction in collagen and moisture content can contribute to the development of facial folds and

wrinkles. This is primarily due to the breakdown of collagen fibers and elastic networks that support the skin. This loss also leads to oxidation, atrophy, and collapse of skin tissue, resulting in various signs of aging such as dehydration, wrinkles, looseness, and lack of elasticity.

9. Will the reduced fat come back after the lipo laser treatment?

A: The number of adipose cells in an adult remains constant, and weight gain is typically due to an increase in the size of these cells. This treatment employs a specific laser to stimulate adipose cells, leading them to release stored fatty acids and triglycerides, thereby reducing cell volume. Consequently, the occurrence of fat returning is highly unlikely. Additionally, maintaining a healthy diet would greatly contribute to preventing this.

9. Ask & Answer

Ask	Cause	Answer
Instrument fails to turn on?	The power cord is disconnected to the wall socket.	Check if the power cord is well connected to the wall socket and ensure that the socket is functional.
	The fuse is loose or burn?	Replace the fuse on the back of the instrument and restart it.
The instrument has stopped functioning and cannot be restarted?	The fuse may have burned out due to unstable voltage or prolonged operation.	Please replace the fuse and restart the instrument 5 minutes later.
The screen shows an error message.	Data error.	<ol style="list-style-type: none"> 1. Turn off the instrument and unplug the power cord, and restart the instrument 1 minutes later. 2. If the above method doesn't work, please contact the dealer for a replacement.

Ask	Cause	Answer
Why is the handle not generating RF energy?	The handle is not properly connected.	Please confirm if the handle plug is well connected to the main unit.
	Oil on the treatment area could result in insufficient contact between the handle and the body.	Apply some mediums on the treatment area.
	The handle is not contact with the body properly.	Ensure complete contact between the treatment head and the skin.
CAVI handle is not functioning.	You might mistakenly think that the handle is not functioning properly due to its low-noise design.	To check if the handle is functioning, please refer to the CAVI handle pretest.
	You may have forgotten to apply the gel or have used a product that is not designed for use with the instrument.	Please use the professional gel for CAVI treatment.
	Insufficient massage strength.	Apply more pressure and massage with greater force.

10.Operation guide

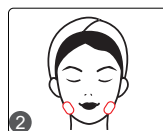
Face:



Step 1: (5 minutes)




Step 4: (20 minutes)



Eye:



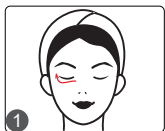Step 1: (10 minutes)



Arm:



Step 1: (15minutes)



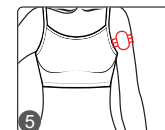
Step 2: (15 minutes)



Step 3: (15 minutes)



Step 4: (20 minutes)



Belly:



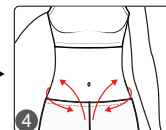
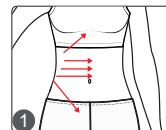
Step 1: (15 minutes)




Step 2: (15 minutes)



Step 3: (15 minutes)



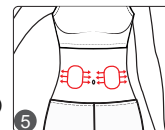Step 4: (20 minutes)




Back:



Step 1: (15 minutes)



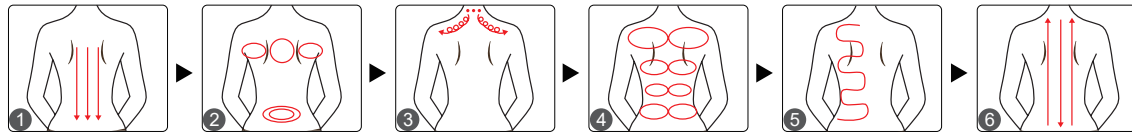Step 2: (15 minutes)



Leg:



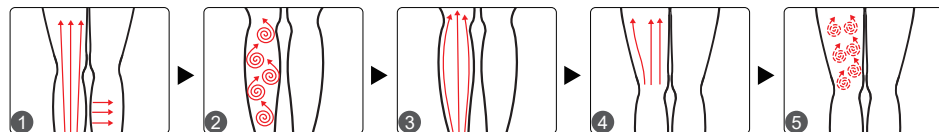
Step 1: (15 minutes)



Step 2: (15 minutes)



Step 3: (15 minutes)



Step 4: (20 minutes)



11. Warranty and support

We offer a one-year warranty and a three-month period for free returns and replacements from the date of purchase.

Product Name: Versatile Beauty Device

Product Model: LS-78D2SB

Production

Date: _____

Quality

Inspection: _____

PASS