

# Usage Demonstration



1 . Product connect to the controller.



2 . Plug in the power cord ( plug not included )



3 . Long press the switch to turn it on or off



4 . Put the hat on the head



I II III IV V VI VII  
✓ APPLICABLE ✗ NOT APPLICABLE

Levels	3 levels adjustable pulse mode
L0 : No Pulse	Suitable for short-term continuous treatment
L1 : 10HZ	Helps induce brainwave entrainment
L2 : 40HZ	Helps restore the brainwaves to a natural order

Product name	Red light therapy hat
Wavelength	660nm&830nm
Number of led light	60pcs*660nm,60pcs*830nm Total 120 pcs led lights
Product perimeter	25.6in/65cm
Product diamete	7.5in/19cm
Product weight	4.8ounce/140g
Due to manual measurement error between 0.4-0.8 inch	

## Functional Features

### Different lighting effects :

Various data show that different wavelengths of light can achieve different penetration requirements . The product has both 660nm and 830nm wavelengths , which can meet the needs of different parts . 60pcs 660nm red light beads . 60pcs 830nm infrared light beads , a total of 120 light beads.

### Premium-quality and light:

Low light wane and long usage lifespan . Foldable and portable, light weight , only 140g ( 4.8 ounce ),no burden to use.

### Optical Design of Convex:

Convex lamp beads can make the light more concentrated so that the light travels in a straight line . It can reach 1-3mm hair follicle root activate the hair follicle and let it re-enter the healthy growth phase.Pure physical therapy with no side effects.

### Automatic timing function:

Once turned on , it will automatically turn off after 20 minutes , no worry about over-treatment . If you need to pause during use, you can long press the switch to turn it off . It is recommended to use 20 minutes each time . 3-4 times a week.

### Function:

Red light therapy hat can promote the proliferation of brain cells, speed up cell energy metabolism , and relax the brain . It also can reduce oil production , decrease shedding , promote hair growth . restore hair growth healthy cycle , and repair damaged hair.

## Caution


- 1.Do not washed and avoid water contact with plugs and connectors. It is recommended to use professional cleaners for wiping.
- 2.Please use the matching power cord , it is recommended to use a 5V1A plug for charging.
- 3.When used with a power bank , the product will be limited by the power of the power bank device , and the power and brightness will be significantly reduced
- 4.Avoid squeezing and over-bending the product when using it.

## Quality and Technical Certification

- 1.The product is produced by Shenzhen Yihong Lighting Co , Ltd. entrusted by Scienlodic Ltd.
- 2.FDA owner number is 10082798.
- 3.The product has passed the European patent , and the US patent is pending.
- 4.Do not misappropriate and imitate without authorization.

## After Sales

If you have any problems in use , please seek help from the seller or professional guidance in time.  
All Scienlodic products with quality problems can be reissued or replaced free of charge within 2 years . Please check and use the product in time , so as not to miss the post-sale period.

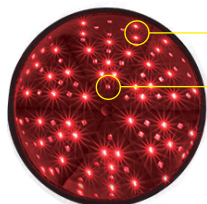

60pcs visible red light

60pcs invisible infrared light

Tips : Visible light for human eyes is about 400-800nm .  
Half of the light on the hat is 830nm invisible infrared light.  
You can use a power meter to test whether the light is bright.

# Instructional Manual

Thank you for purchasing our light,we hope you can reap benefits of red light therapy at home.

