

Thank you for choosing our red light therapy panel. We appreciate your support and hope you'll enjoy this device.

F© (€

Red Light Therapy Pad Instructional Manual



Thank you for purchasing our light, we hope you can reap the benefits of red light therapy at home.

Say hello to your red light!

Before starting your first session, please make sure to read this manual as it may explain important operating procedures.

Notes: Each LED inside with three chips . Near Infrared 850nm is invisible to the human naked eves.

The red light device uses a combo

of red 660nm and Infrared 850nm

LEDs which is chosen to give the

ideal wavelengths found in the

decades.

tissue.

660nm of red light

eyes and targets skin cells.

850nm of near-infrared light

This light cannot be seen by the

human eyes and targets deep

scientific research of the last three

This light can be seen by the human

This package includes

1*Red light Therapy Sleep Bag

2*Power Cord

2*Adaptor

2*Connecting wire

1*User Manual


Specifications

| Model: | Red Light therapy Sleep Bag | Wave lenath : | 660nm / 850nm |
|--------------|---|---------------|---------------------------|
| Size : | 180 * 80*20cm | Function : | Body slim / Pain Relief |
| Color : | Black | Custom : | Logo / Shape / Wavelength |
| Power: | 512W | Power Bank: | Not |
| Led Number : | 2560pcs (2048pcs 660nm leds, 512pcs 850nm leds) | | |


Red Light

Therapy Sleep bag

2048Red Lights + 512 Infrared Lights , Unique pulse function , the effect is more obvious

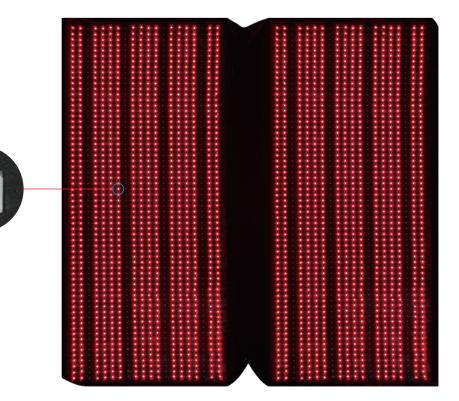


CONTROL PANEL

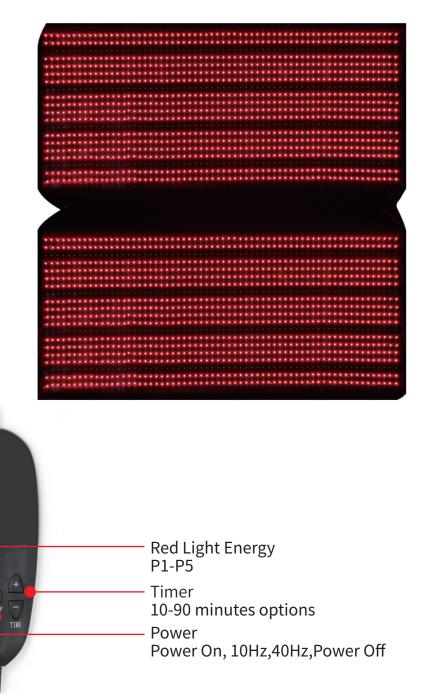


Benefits backed by science:

Research shows that the wavelengths in the spectrum of 660-900nm show optimal biotical response.

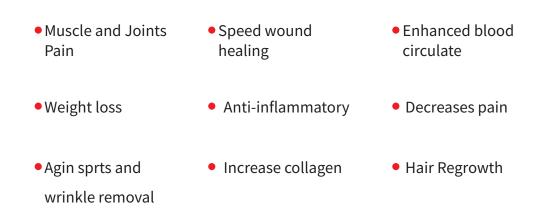

Even though this technology doesn't have side effects, we suggest taking the following precautions:

- ▶ This device is not waterproof!
- ▶ Not intended for use by children and pregnant women.
- ▶ Do not exceed 30 Min. on one same treatment.

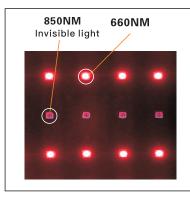



Features:


1 : Perfect performance for body treatments like weight loss and pain relief . 2 Protable design so you can use it in different situations.

3. 2560 pcs leds with wavelength 660nm red light and 850nm infrared light (this light can not see by our own eyes , but they can be lighten up , you can check with you phone camera)
4. We provide lyear warranty for our pad.

NEAR INFRARED – INFRARED A (IR A) Usually use 660nm–850nm

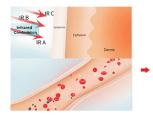

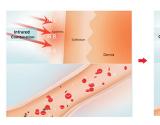


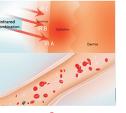

Invisible near infrared energy are absorbed by muscle, subcutaneous tissue which results in THERMAL EFFECT your body responds to this radiation heat with increased blood flow, allowing the parts of the body to receive the oxygen and nutrients they need in order to function more effectively, regeneration is stimulated.

Measured value



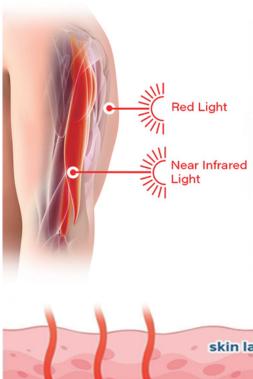


Red Light at 660nm

This specific wavelength is absorbed at a higher proportion by theskin tissue, making it especially effective for boosting skin healthand increasing collagen production.

Near Infrared Light at 880nm

This wavelength is actually outside of the visual spectrum, so it mayappear as if the NIR lights are not working, butthey just can't beenseen by the naked eye. Near infraredlight has a particularly strongability to penetrate deeperinto the tissue, organs, and joints, makingit ideal forenhancing muscle recovery and reducing joint pain

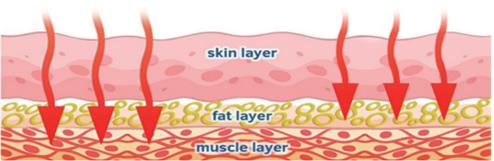

Application

Here show you some examples of how to use the red light device. You and your pets can enjoy the red light in different ways like shoulder, back, feet knee, etc.

LED THERAPY LIGHT PAD FOR JOINTS PAIN RELIEF