


User Manual

Home SPA Instrument

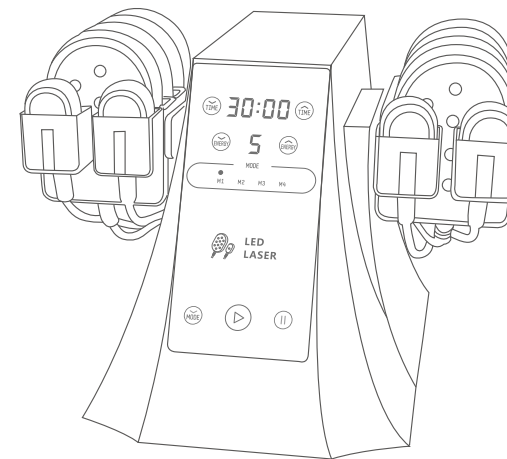

Thank you for purchasing our product.

- Please read this User's Manual carefully for information on how to use the device properly and safely.
- After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
- Make sure to include this User's Manual if you are transferring ownership of the device to another person.

Contents

Part I

Safety Precautions	1
Contraindications	2
Overview	2
Interface Function	3
Technical Specification	4
Precautions	5
Maintenance & Cleaning	6
FAQ	6
Ask & Answer	7
Operation Guide	8
Warranty & Support	9

Safety Precautions

- 1.Please read this manual carefully before use and familiarize yourself with the operation of this instrument and related technology.
- 2.This instrument is a high-precision electronic device that requires proper storage. Please avoid storing it in humid conditions, direct sunlight, or near flammables and explosives to prevent damage.
- 3.Always use the designated accessory provided or recommended by the manufacturer to ensure performance throughout its service life.
- 4.Do not place the instrument in high-temperature conditions, as this may negatively impact its service life and performance.
- 5.Patients should refrain from using this instrument unless they have obtained permission from their doctor.
- 6.Power off the instrument when not in use to ensure its safety.
- 7.If you are currently taking weight-loss drugs, please discontinue their use 1 - 2 months before using this instrument. Failure to do so may result in an extended course of treatment.
- 8.Do not apply the instrument on eyes, thyroid, parathyroid, testicles, pregnant women's abdomen, pacemakers, etc.

Contraindications


Before using this instrument, it is essential to consult your doctor or seek advice from professionals.
Details as below:

1. Pregnant or lactating woman.
2. Individuals with hypertension, heart disease, or a pacemaker installed.
3. Individuals with unhealed wounds or those in the recovery period from surgery.
4. Individuals suffering from epilepsy, severe diabetes, and hyperthyroidism.
5. Individuals suffering from malignant tumors, hemophilia, and hemorrhage.
6. Individuals with skin disease or infectious disease.
7. Individuals have had a cesarean section within six months, a miscarriage within three months, or a natural birth within two months.
8. People who had plastic surgery, with implants or metal objects at the treatment area in 3 months.



Overview

With the most popular and safest technology for fat removal, metabolism enhancement, and body shaping, this machine promises you the desired body figure and glowing face. It provides the perfect solution for people seeking improvements in facial skin appearance and body figures without injections, drugs, or side effects. For external use with immediate effects.

Interface Function

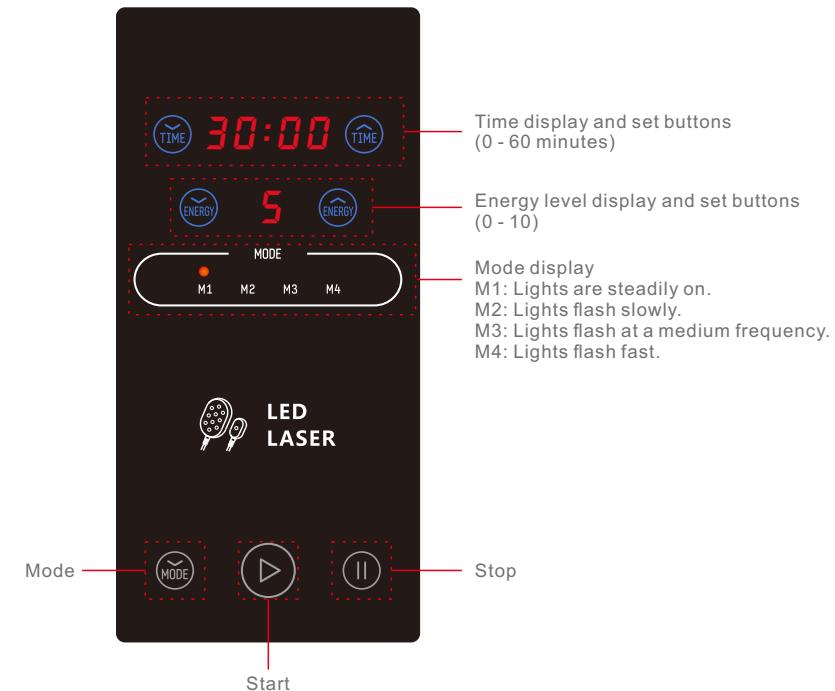

Technical Specification

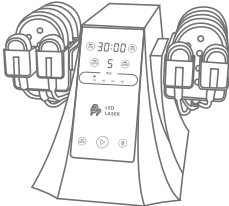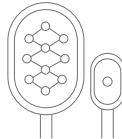	Power supply Input: AC110V-220V 50-60Hz		Net Weight: 3.5kg
	Power:	16W	Gross Weight: 4.8kg
	Dimensions:	11×13.5×8.2inches	
	Input Power:	16W	
	Wavelength: 650nm 5mw		
	Power: 16W		

Precautions

- 1.Strictly follow the instructions in the manual to use the instrument and train operator.
- 2.Before using the instrument, use a plug with a grounding pin and ensure the power socket is grounded.
- 3.Make sure the voltage of the instrument is adapted. If the local power supply voltage is unstable, we recommend adding a power-matched regulated power supply between the mains and the instrument.
- 4.Remove all the metal objects inside the body to avoid causing accidents or affecting treatment effects.
- 5.Do not use the instrument on an empty abdomen, and take the treatment 1 hour later after the meal.
- 6.Fix the laser pads on the treated individual first, then switch on the instrument, and gradually adjust the energy intensity from low to high to ensure a warm and comfortable experience for the user.
- 7.Make sure the laser pads are in full contact with the skin during the procedure, and avoid partial contact.
- 8.You can enhance the effectiveness of this instrument by combining it with other beauty devices. However, if you use this instrument alone, simply fix the laser pad on the treatment area, and wait for 20 minutes.

Maintenance & Cleaning

- 1.Disconnect the power cord when not in use and store the instrument in a dry, cool, and well-ventilated place.
- 2.After cleaning the instrument, cover it with a dustcloth or plastic wrap. Before use, ensure to sanitize the handle again for hygienic reasons.
- 3.Disassembling the instrument and attempting to clean the inside is strictly prohibited. This task should only be performed by professionals.

FAQs

1.What is the function of the laser?

A: The laser is a special form of light that possesses the ability to penetrate the skin's surface. It operates by effectively eliminating fat deposits, stimulating metabolic processes, rejuvenating the soft tissue within the epidermis, and triggering the regeneration of collagen deep within the skin. These actions collectively result in skin that is notably smoother and significantly more elastic in texture.

2.When can I expect to see the results?

A: Typically, you can observe the effects immediately following your treatment. Lipo laser treatments work by dissolving fat and boosting metabolism. The dissolved fat is subsequently eliminated from the body through urine and sweat. This process does take some time, meaning that the more treatments you take, the more pronounced the effects will become.

3.Will I regain weight after the lipo laser treatment?

A: The number of adipose cells in an adult remains constant, and weight gain is typically due to an increase in the size of these cells. This treatment employs a specific laser to stimulate adipose cells, leading them to release stored fatty acids and triglycerides, thereby reducing cell volume. Consequently, the occurrence of fat returning is highly unlikely. Additionally, maintaining a healthy diet would greatly contribute to preventing it.

Ask & Answer

Ask		Cause		Answer
Fail to turn on the instrument.	→	The power cord is disconnected to the wall socket.	→	Check if the power cord is well connected to the wall socket and ensure that the socket is functional.
	→	The fuse is loose or burned.	→	Replace the fuse on the back of the instrument and restart it.
The instrument is not working and can't be restarted?	→	The fuse may have burned out due to unstable voltage or prolonged use.	→	Please replace the fuse and restart the instrument 5 minutes later.
The screen shows an error message.	→	Data error?	→	1.Turn off the instrument and unplug the power cord, and restart the instrument 1 minutes later. 2.If the above method doesn't work, please contact with the dealer for a replacement.

Operational Guidance

Face




20mins



Arm




20mins

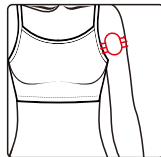

Abdomen




20mins

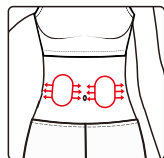

Leg



20mins

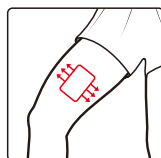

Warranty & Support

We offer a one-year warranty on the main unit and a three-month warranty on components.

Product

Name: Lipo Laser Body Sculptor

Product

Model: LY-12101J

Production

Date: _____

Quality

Inspection: _____

PASS