
User Manual

Home Beauty Instrument

MS-11R9

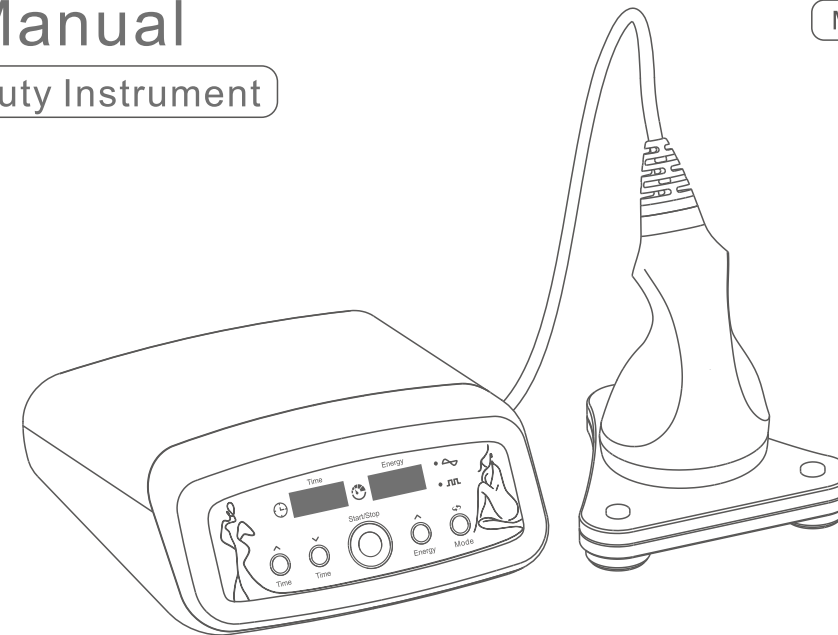Thank you for purchasing our product.

- Please read this User's Manual carefully for information on how to use the device properly and safely.
- After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
- Make sure to include this User's Manual if you are transferring ownership of the device to another person.

Table of contents

Introduction	1
Efficacy	1
Precautions	1
Contraindications	2
Product parameters	3
Packing list	4
User guide	5
After treatment	7
Steps	8
Treatment frequency recommended	11
Troubleshooting	12
FAQ	12
After-sales guarantee and support	14

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Introduction

Thank you for choosing Home Beauty Instrument MS-11R9, please read this manual carefully before you use the instrument!

The instrument adopts the principle of high frequency sound waves, using strong acoustic wave energy to hit human body fat cells, causing friction movement among fat cells to achieve the effect of body shaping. The instrument can quickly consume body calories and cellular water, and shrink fat cells, thus achieving the effect of fat removal. The instrument is currently widely used in body shaping field and is well received by consumers.

Efficacy

1. Improve the waist and abdomen flab, love handles, and belly fat.
2. Improve the bingo wings, thick arms, arm fat accumulation.
3. Tighten excess leg flab, say goodbye to thick legs.
4. Accelerate metabolism, improve constipation and bowel movement.
5. Improve the obesity caused by overeating.

Precautions ⚠

1. Please read this manual carefully before use so you are familiar with the use of this instrument and related technology.
2. The beauty instrument is a high-precision electronic instrument, so it needs to be stored properly. Do not store it in a place that is damp, exposed to direct sunlight, or around with inflammable and explosive materials to prevent damage.

- 3. Always use the designated accessories provided or recommended by the original manufacturer to ensure the performance of the instrument within normal service life.
- 4. Do not place the instrument near strong heat, as this may affect its service life and performance.
- 5. Do not use instrument on the eyes, head, chest, back, thyroid, parathyroid, testicles, abdomen of pregnant women and pacemaker.
- 6. The instrument is not intended for use by patients, unless they have the permission of a doctor.
- 7. During the operation, the operating part should stay moist, avoid operation on dry skin.
- 8. Please turn off the power of the instrument when not in use, and make sure to turn off the general power of the instrument before leaving after daily use to ensure the safety of instrument.
- 9. If you are taking other weight-loss drugs, it is recommended to stop the drug for 1-2 months before using the instrument, otherwise the course of treatment needs to be extended.
- 10. Clean the instrument with wet wipe after operation to ensure the cleanliness and hygiene, and prolong its service life.

Contraindications


Please do not use the instrument or consult your doctor before use if you:

- 1. Are menstruating, pregnant and lactating women.
- 2. Are suffering from severe disease, such as diabetes, heart disease, history of scarring, hyperthyroidism, etc.
- 3. Have recently had surgery or are recovering from a surgery.
- 4. Have malignant tumors, hemophilia or severe bleeding.
- 5. Have skin diseases and infectious diseases.
- 6. Are wearing pacemakers (heart rate units and other internal electronic devices), vascular stents or metal implants.


- 7. Have had surgery or cesarean section within 6 months, miscarriage within 3 months, or natural birth within 2 months.
- 8. Are drunk within 12 hours, full within half an hour, fasting, thirsty, and overworked.
- 9. Are a minor or with reduced sensory or mental capabilities.

Product parameters

Product name	Home Beauty Instrument
Product model	MS-11R9
Frequency	40kHz
Power	30W
Product size	18.5x19.5x8cm (7.3x7.7x3.2inch)
Product weight	0.9kg (1.98lb) (excluding 40k head)

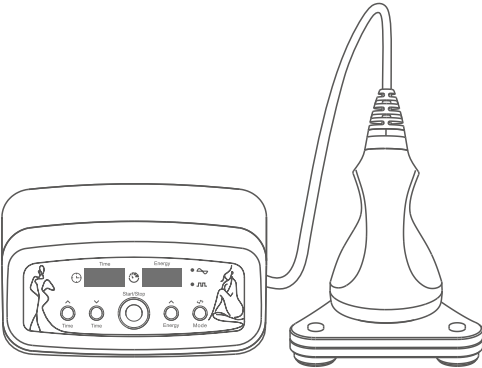
Packing list

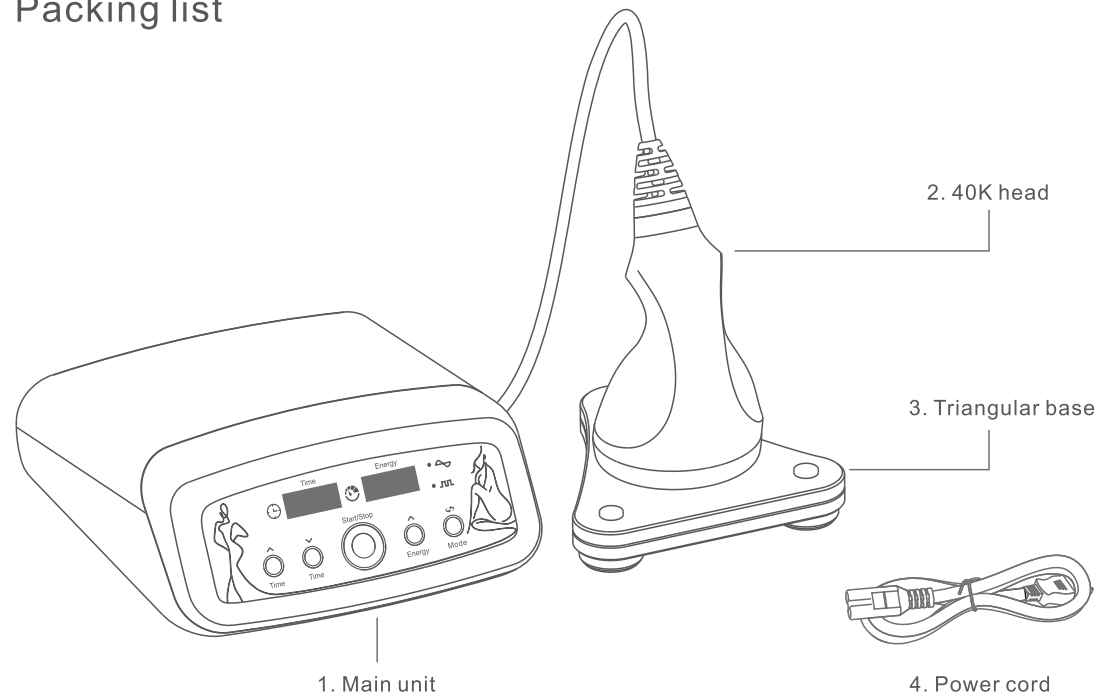
Note: Please check the instrument and its accessories or attachments before use, and remove all packaging materials after confirming that they are complete.

User guide

1. Head installation

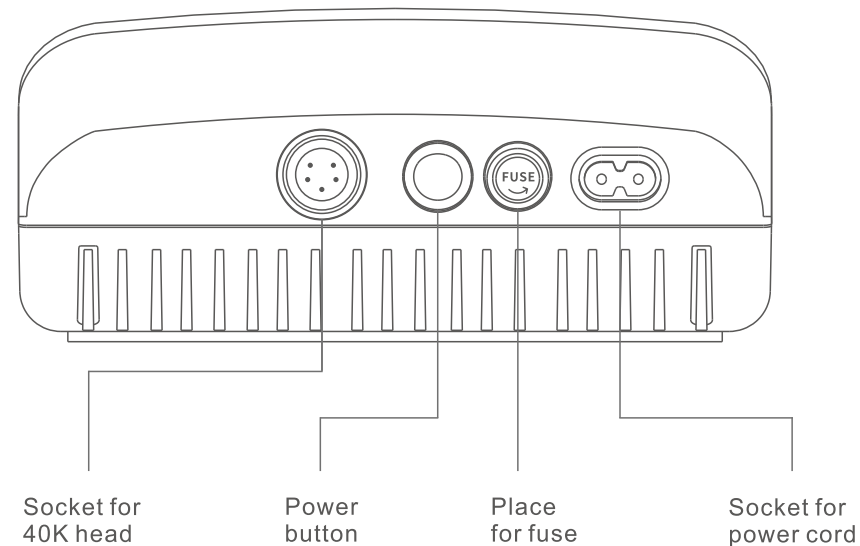2. Interface setting

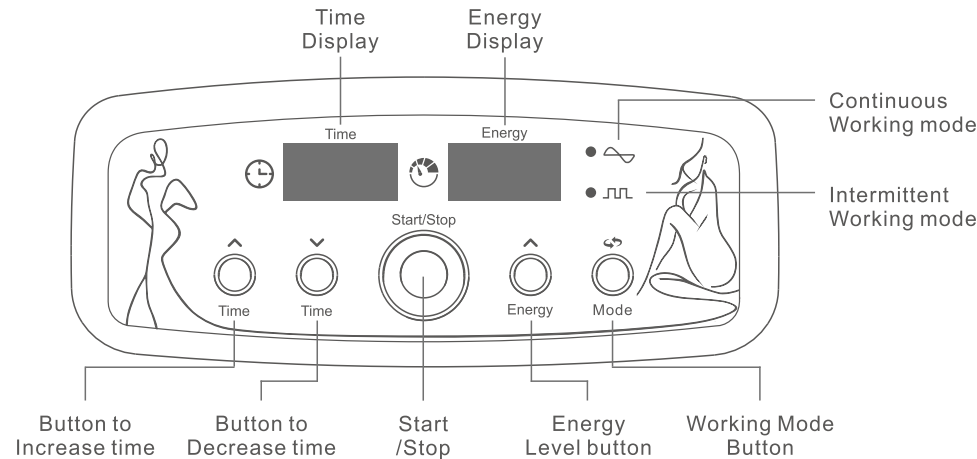Energy level button

(if you need to reduce the energy, press the energy button one more time after setting to the highest level, it will automatically skip to the lowest level.)

Continuous working mode

Suitable for people with more fat accumulation and thicker fat.

Intermittent working mode

Suitable for people who want to improve body shape and local obesity.

3. 40K head test

1. Pour some water into the triangular base.

2. Turn on the instrument and set the energy to between 6 and 9, set to continuous working mode [continuous wave icon], then press the start button.

3. Immerse the 40K head in water for about 1mm, but do not fully immerse the head (metal part) in water. When you hear the sound of "zzz", the head is working properly. If the head does not work properly, please contact the seller by email to get a new head replacement.

After treatment ⚠

1. Keep warm and do not eat cold food after the operation with 24 hours after treatment.
2. Do not take a shower within 4 hours after treatment, avoid catching a chill by the wind.
3. Avoid spicy and greasy food.
4. Drink more warm water to accelerate metabolism.
5. Avoid staying up late, smoking and alcohol.
6. Avoid sauna or hot springs in the short term.

Steps


Please prepare the gel and towel before operation.

1. Abdomen(Parameter setting suggestions)


Duration: 30 minutes

Energy level: between 3-7


Mode: continuous working mode

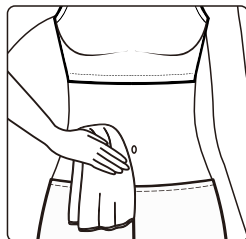
1. Wipe the treatment area with a wet towel.

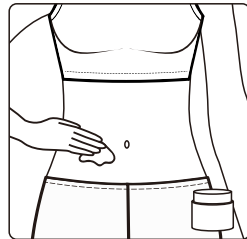
2. Apply the gel evenly to the treatment area.

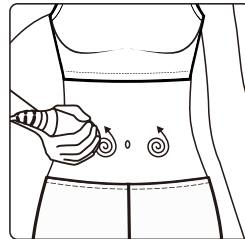
3. Move the head in small circles evenly on your abdomen for 5 minutes.

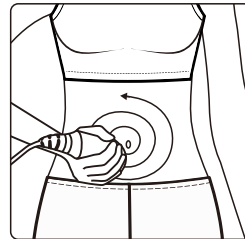
4. Move the head in large circles evenly around your navel for 5 minutes.

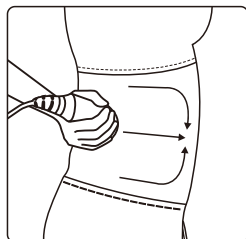
5. Move from the flank (side of the waist) towards the middle of your abdomen (3 minutes each right/left side).

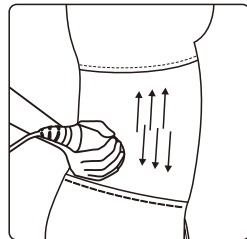6. Move upwards and downwards on the both sides of the waist and abdomen (3 minutes each right/left side).

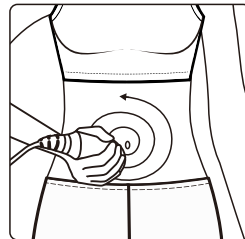7. Move the head in circles on your abdomen for 8 minutes for reinforcement.

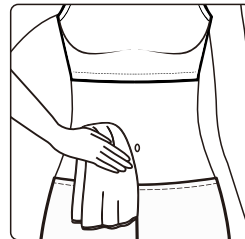
8. Wipe the treatment area.

2. Thighs(Parameter setting suggestions)

Duration: 30 minutes

Energy level: between 3-7


Mode: continuous working mode

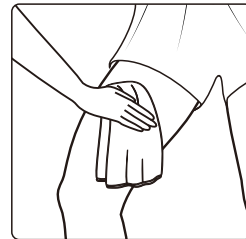
1. Wipe the treatment area with a wet towel.


2. Apply the gel evenly to the treatment area.

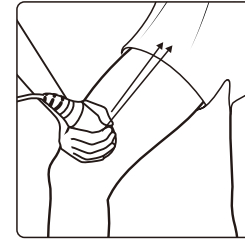
3. Move upwards on the front thigh from your knee towards your groin (5 minutes each right/left side).

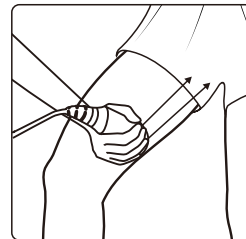4. Move upwards on the inner thigh from your knee towards the top of your thigh (5 minutes each right/left side).

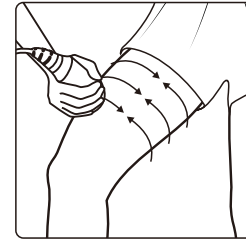5. Move from the sides of your thigh towards the middle (5 minutes each right/left side).

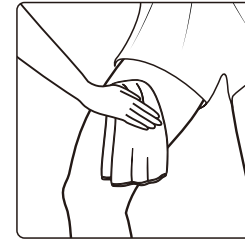6. Wipe the treatment area.

3. Upper Arms(Parameter setting suggestions)


Duration: 20 minutes Energy level: between 3-7 Mode: continuous working mode


1. Wipe the treatment area with a wet towel.

2. Apply the gel evenly to the treatment area.

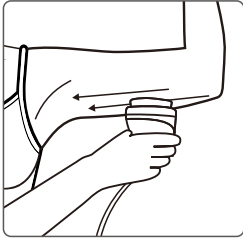3. Move from elbow towards underarm or outside of upper arm (5 minutes each right/left side).

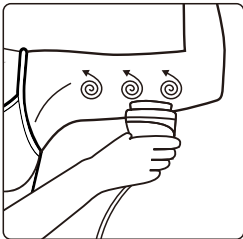4. Move the head in small circles on upper arm (5 minutes each right/left side).

5. Wipe the treatment area.

Treatment frequency recommended

Operate 3-4 times a week, leave an interval of 1 day between treatments (stop the treatment during menstruation). During the treatment period, try to avoid high-calorie food and cooperate with moderate exercise to achieve better results! In the later stage, operate once in 3-5 days for consolidation!

SUN	MON	TUE	WED	THU	FRI	SAT
30	1 ✓	2	3 ✓	4	5 ✓	6
7 ✓	8	9 ✓	10	11 ✓	12	13 ✓
14	15 ✓	16	17 ✓	18	19 ✓	20
21 ✓	22	23 ✓	24	25 ✓	26	27 ✓
28	29 ✓	30	1 ✓	2	3 ✓	4

Troubleshooting

Problem	Possible causes and solutions
1. The instrument cannot be turned on?	→ A. Please check if the power cord is connected to a valid socket.
	→ B. Please replace the fuse with the spare fuse in the fuse box.
2. The main unit can be turned on, but the 40K head does not work?	→ A. Please check whether the 40k head plug into main unit's socket tightly.
	→ B. The 40k head is damaged, then cause it does not work.
3. Why is the treatment head getting warm after a period of operation?	→ A. Check if gel or other media products are applied during operation.
	→ B. Please switch on the instrument after the treatment head touches the skin.
	→ C. The operation time should not exceed 30 minutes.

FAQ

1. Why is there tinnitus during 40K operation?
A: The high frequency sound wave power is very strong; the sound waves with a frequency higher than 20kHz act on the fat layer 20mm deep under human skin through the high frequency sound wave focusing. The high frequency sound wave with concentrated energy effectively causes the high-speed friction and heat among the fat cells in the focal area to break and emulsify them, this process emits sound, and that's where the tinnitus phenomenon comes from.

2. Will I regain weight after treatment?
A: It is not easy for you to regain weight after taking treatment. This treatment reduces fat cells, not water, and the formation of fat cells takes a long time to accumulate, so it is not easy to gain it back.

3. Do I need to control diet during the body shaping period?
A: You need to control diet a bit. After taking the treatment, fat-blasting will accelerate metabolism, avoid eating spicy, greasy, and fried foods, so as not to affect the metabolism. Meanwhile appropriate exercise can make the excess fat smoothly excreted from your body, and the body shaping effect will be more obvious.

4. Can I use 40K head all over the body?
A: Do not use 40K head on heart and area around eyes. The power of high frequency sound wave is very strong, and the heart is pretty sensitive to sound waves since the it is a hollow organ, after being vibrated by sound waves, due to the different conduction of sound waves by muscle tissue and blood in the heart, back and forth reflections will occur, this force can pull the valve away from the heart muscle, if it is aimed at the eye, it will cause retinal detachment (Head, low back, chest is prohibited.)

5. Does 40K head perform any side effects on the human body?
A: High frequency sound wave is a non-surgical, non-invasive procedure that does not require incision or anesthesia. The principle of high frequency sound wave: the focused strong sound waves enter the body and produce a strong impact on human fat cells to make them burst.
The effect: fat breaking – fat dissolution – skin tightening – body shaping.
Therefore, it only targets the low-density fat tissue and protects the high-density tissue such as blood vessels and nerves, and has no side effects on the human body. There will be slight tinnitus during the operation, which is normal and need not be worried.

6. Can I set the frequency of the 40K head to between 30kHz and 40kHz? Or is it constantly 40k?
A: The frequency of 40K head cannot be adjusted, it is a head with a fixed 40kHz frequency.

7. How often do I need to take the treatment? How long does each treatment take?
A: Take the treatment once every other day, and for 30 – 40 minutes on a single treatment area, and at least 3 times a week.

8. How long did it take to get these results?
A: Generally, you can feel the slight changes in the operating parts for about a month, and the effect will become more and more obvious if you stick to 2 - 3 months.

9. Can I use the 40K head under chin? Is that dangerous?

A: This instrument is forbidden to be used on the head and near the heart, the chin is also part of the head, please do not use it on the head!

10. What products can I use with the 40k head for treatment?

A: Gel, cream, lotion, etc. are all available, but relatively speaking, gel is more moisturizing, and it is not easy to get dry during operation. As a medium, gel is easier to receive energy and transmit it to the body, which is more suitable as a companion product for 40K head!

11. How to maintain my beauty instrument and clean my 40k head?

Turn off the power, unplug the power adapter, wipe the 40k head clean with a wet tissue. (Note: Do not use alkaline cleaners, alcohol, gasoline, nail polish remover and other solvents to clean)

After-sales guarantee and support

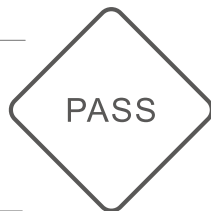
Free return is provided within 90 days starting from the date of purchase.

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