


User Manual

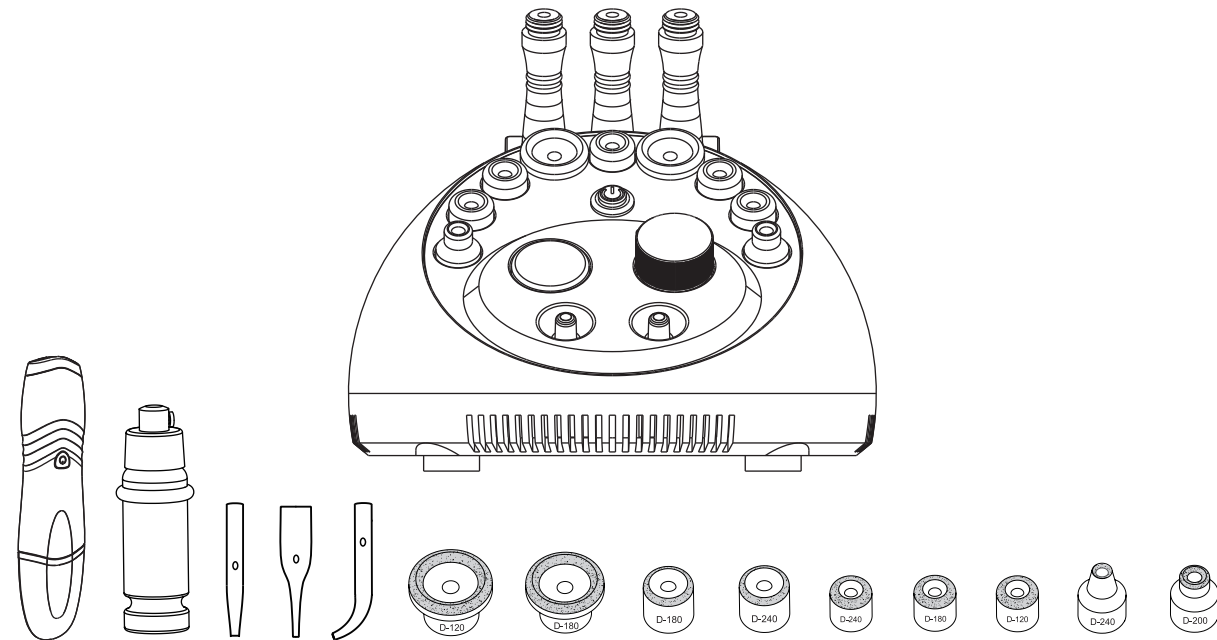

Catalogs

Part I

Security Precautions	1
Product Overview	3
User guide	5
FAQ	16
Specification Information	17

Thank you for purchasing our product and we hope you have a pleasant experience using it. We strongly recommend that you read this user manual carefully before you start using this product to get the best experience.

Security Precautions

contraindication ⚠



Not all users are suitable to use this diamond microdermabrasion machine, do not use this machine or consult your doctor if you have the following conditions:

- 1.Pregnant or lactating women.
- 2.People with heart disease or pacemakers.
- 3.People whose wounds have not yet healed and who are recovering from surgery.
- 4.Patients with epilepsy, severe diabetes and hyperthyroidism.
- 5.Those with malignant tumors, hemophilia and hemorrhage.
- 6.People with skin diseases, infectious diseases.
- 7.The machine cannot be used by those who have had a planter birth within six months, a miscarriage within three months, or a natural birth within two months.
- 8.People who have just had plastic surgery within three months, implantation of prosthesis, metallic substances in the treatment area, etc.


Precautions

1. Follow the instructions in this manual to operate the unit.
2. When using this device, please keep the target area moist and avoid using it on dry skin.
3. If you are taking weight loss medication, it is recommended that you stop taking the medication for 1-2 months before using the device, otherwise the treatment cycle may be prolonged.
4. Do not disassemble or replace any parts or accessories of this equipment by yourself.
5. Do not attempt to repair or maintain this equipment yourself.
6. When in use, place the unit on a flat and stable location.
7. Make sure the power is turned off before cleaning the unit.
8. Avoid excessive use of essential oils or other cosmetics when using this device.
9. After use, please adjust the air pressure to the lowest level.
10. It is not recommended to use the unit for scrubbing operations for more than 30 minutes at a time.
11. Please keep this equipment away from humid environment, direct sunlight or flammable and explosive materials to avoid damage.
12. Do not place this equipment near high temperatures as this may affect its life and performance.
13. Do not use this device in sensitive areas such as the eyes, thyroid, parathyroid or testicles, on the abdomen of a pregnant woman or near a pacemaker.
14. Outlet and suction holes must be clearly distinguishable, do not insert the diamond micro-carving into the outlet hole, otherwise there is no way to work properly.


Product Overview

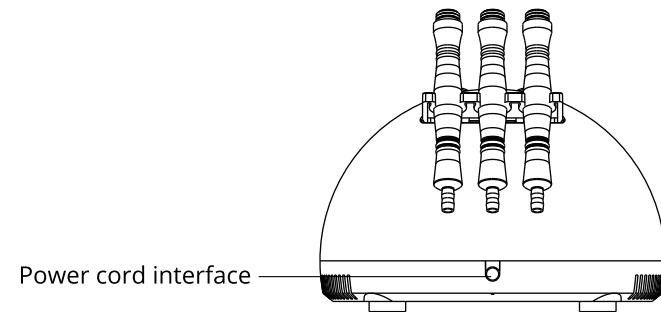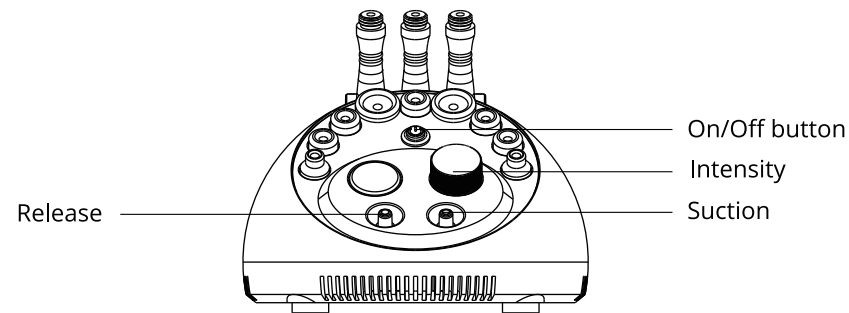


How to Choose the Right Diamond Tips?

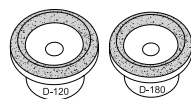


For knee / elbow areas and striae of pregnancy.

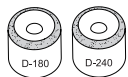

Face and for the larger blemishes.

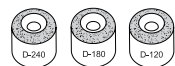

For Scars, fine lines and imperfections.

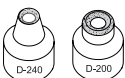

For Eyes and both sides of nose, for wrinkles, oil teblets, blackheads.

Note:

1-The higher the D number,the finer the grit.

2-Clean diamond tips with Alcohol,then sterilize them with UV light or high temperature.


3-Replace the cotton filter inside the diamond tip or gas nipple after 5-6 customers treatment.

User Guide


1. Microdermabrasion

1.1 How to install the Microdermabrasion Wand?


1. Put a cotton filter in to the replacement tip.

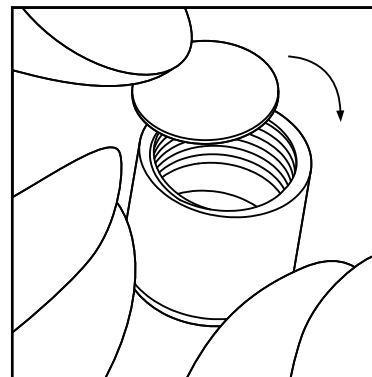


2. Install the replacement tip on the stainless steel wand.

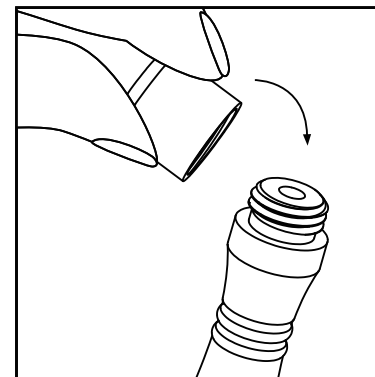


3. Unscrew the cap on the other end of stainless steel wand.

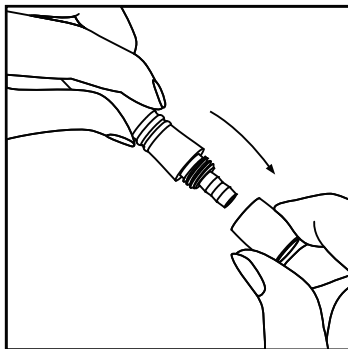

4. Connect the connecting hose and stainless steel wand and screw the cap.

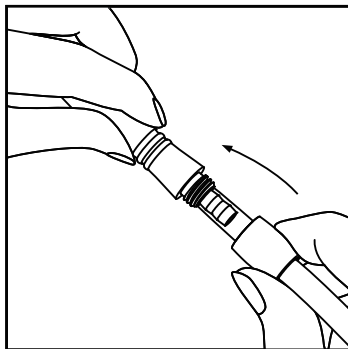


5. Connect the other end of connecting hose and suction port.

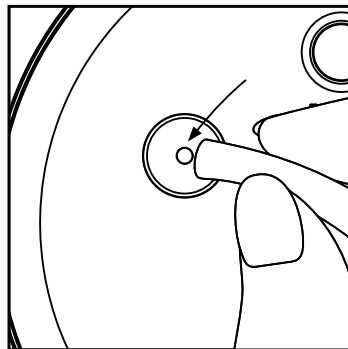


1.2 How to select diamond replacement tip?


D-180 & D-120 for body large areas or body joints.

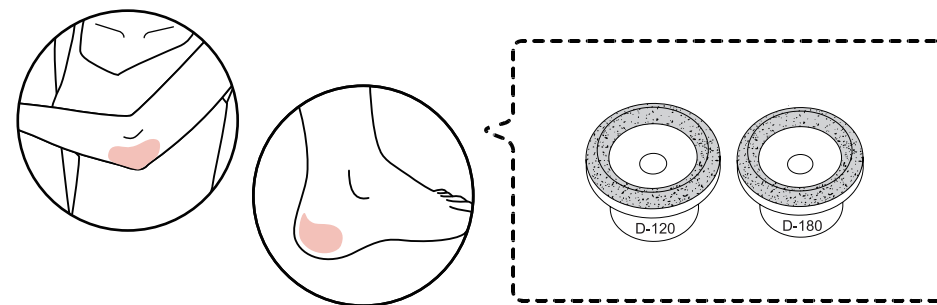

D-240 & D-180 & D-120 for facial large areas.

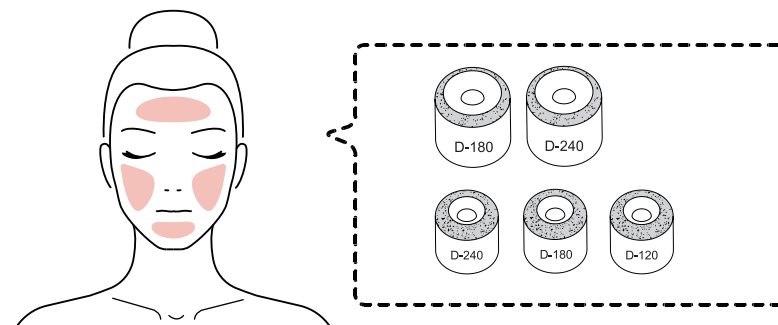

D-200 & D-240 for small areas around eyes or nose on face.

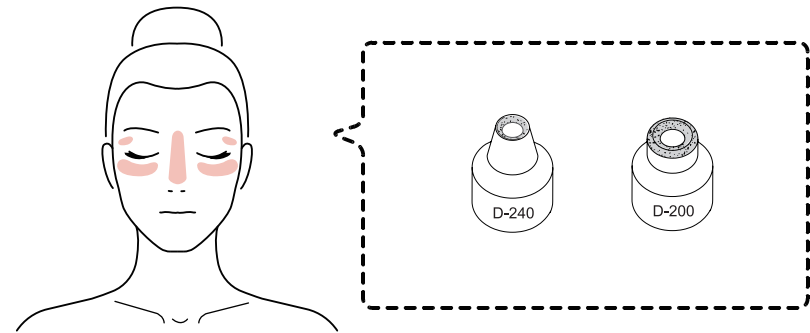


1.3 How to clean diamond replacement tip?



Wipe it with water or soapy water after use.

Caution :

- ① Replace diamond tip every 3 months.
- ② Replace the cotton filter in the diamond tip after each use of the microdermabrasion wand.

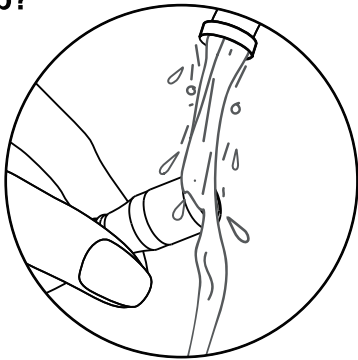

2.Ultrasonic Skin Scrubber

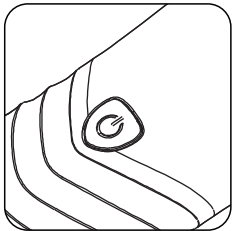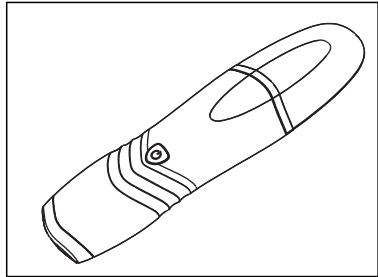

Power and Mode Switch Button


2.1 3 Modes Available

Mode 1	Click once and the skin scrubber will work continuously. This mode can solve the skin problems of enlarged pores, clogged and rough skin.
Mode 2	Click twice and the skin scrubber will work slowly and intermittently. This mode is suitable for people with clogged pores and excessive oil secretion.
Mode 3	Click three times and the skin scrubber will work quickly and intermittently. This mode can exfoliate the skin and tighten your skin.


3. Glass Tube

3.1 How to select Glass Tube


The flat glass tube for facial massage.

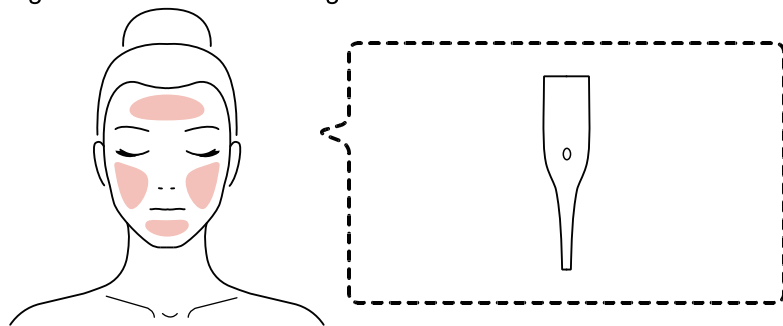

The straight glass tube for cleaning up severe clogged pores.

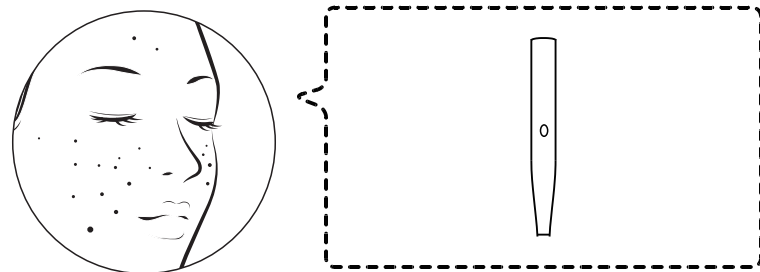


The bend glass tube for blackheads removal.

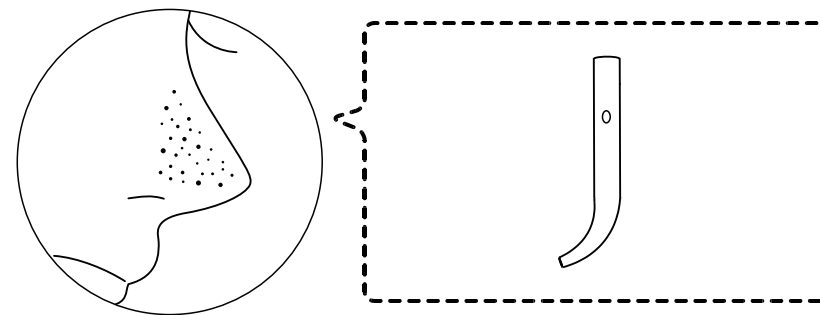

3.2 Caution ⚠


Press the hole on the glass tube with your finger to avoid air leakage when using.

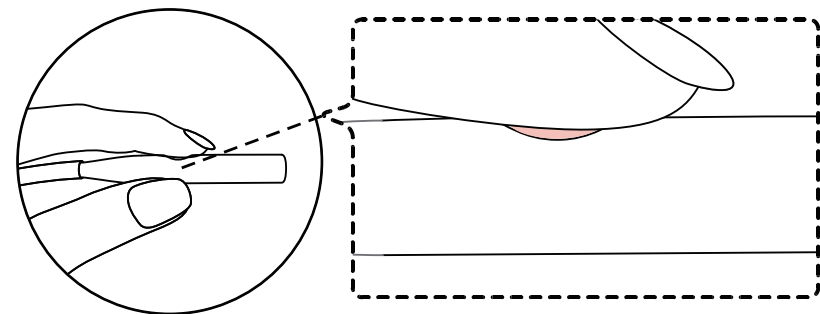


4. Spray Bottle

Caution : ⚠


1. No high concentration solution.

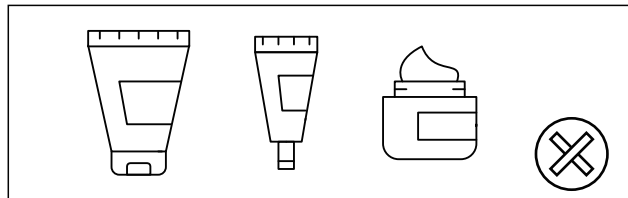

2. Press the hole on top with your finger when using.

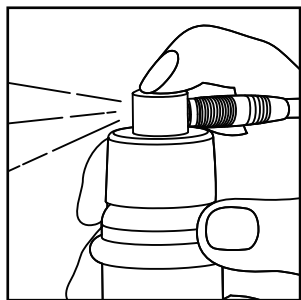


3. Connecting hose and release port should stay connected.

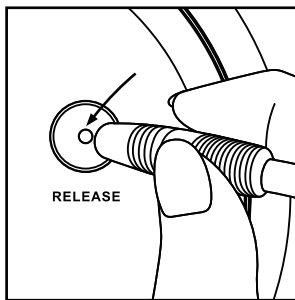

Recommended Application

1. Open pores with a facial steamer or hot towel.

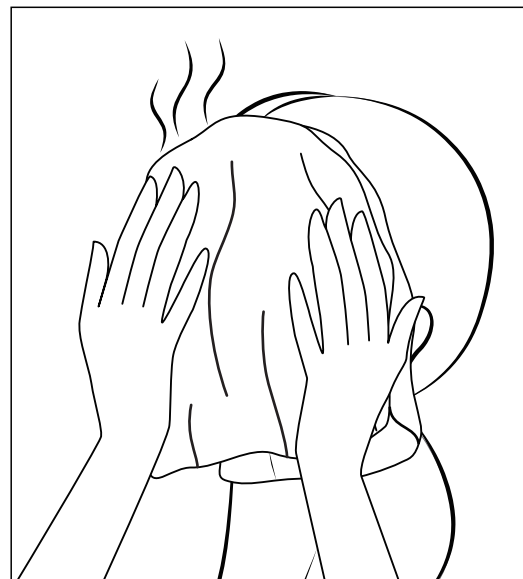


2. Exfoliate and reduce wrinkles with microdermabrasion wand.




3. Remove blackheads and whiteheads with glass tube.



4. Hydrate facial skin with spray bottle.

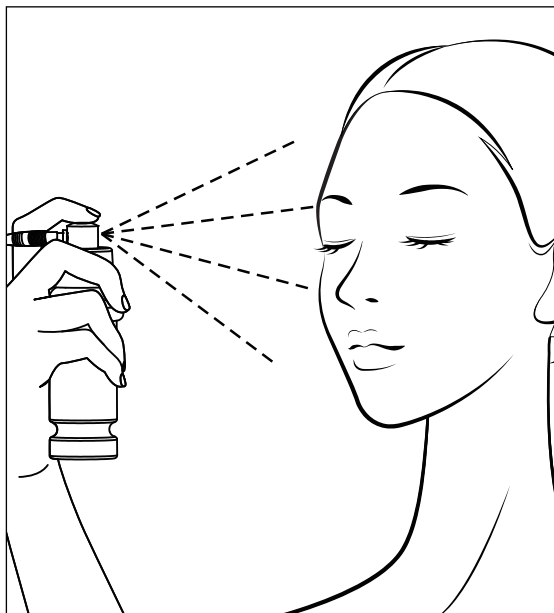


5. Apply a mask.

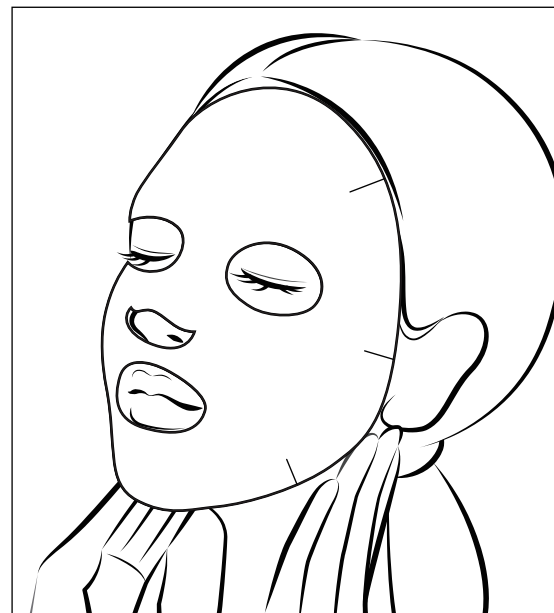

Note: ⚠

After Microdermabrasion Treatment, you can use soothing products, collagen and nourishing facial mask for facial massage, and finally use moisturizing cream with high protection coefficient. In the first 48 hours, you should not expose yourself to the sun, and do not make up after about four hours.

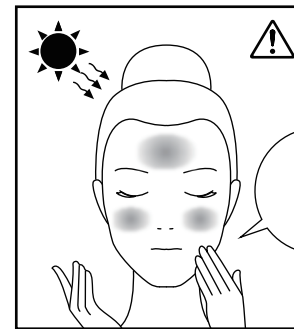

FAQ

1

Q: Why does my face get redness after Microdermabrasion?

A: Microdermabrasion Machine adopts the principle of negative vacuum pressure, when the microdermabrasion wand or glass tube is applied to your face, the capillaries of your face will be congested with blood, so the appearance of skin redness is normal. It should disappear in few hours.

2

Q: How long does it take to see results?


A: After the first treatment, you skin will be softer and smoother. After 8-12 weeks of continuous treatment, you may see a reduction in the appearance of fine lines, wrinkles, dark spots and hyperpigmentation.

3

Q: Why the Spray Bottle fails to spray?

A: ① It is possible that the concentration of the liquid in the bottle is too high, which causing the blockage.
② You are not pressing the top hole with your finger during use.

Specification Information

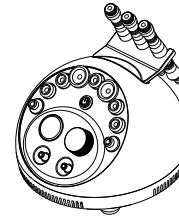

Power supply Input:
AC110V-240v 50HZ-60Hz

Net Weight: 0.9kg

Power: 3W

Gross Weight: 1.3kg

Dimensions: 6.7 x 6.7 x 4.7in

Input Power: 3W