

User Manual

Microdermabrasion Machine

MS-33P3N

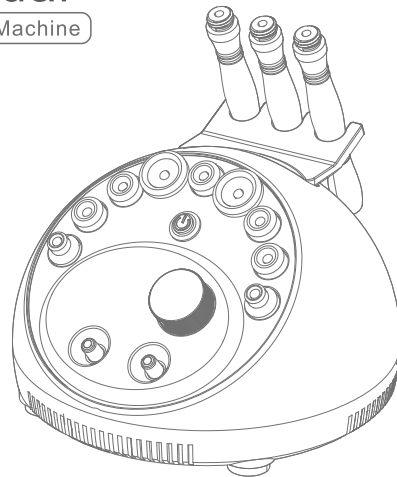


Table of contents **EN**

Safety Precautions	1
Contraindications	2
Overview	2
Interface Introduction	3
Usage Details	4
Technical Specification	7
Use Precautions	8
Maintenance & Cleaning	9
FAQ	9
Ask & Answer	11
Operation Guide	12
Warranty & Support	13

Safety Precautions

1. Please read this manual carefully before use and familiarize yourself with the operation of this instrument and related technology.
2. This instrument is a high-precision electronic device that requires proper storage. Please avoid storing it in humid conditions, direct sunlight, or near flammables and explosives to prevent damage.
3. Always use the designated accessory provided or recommended by the manufacturer to ensure performance throughout its service life.
4. Do not place the instrument in high-temperature conditions, as this may negatively impact its service life and performance.
5. Do not apply the instrument on eyes, thyroid, parathyroid, testicles, pregnant women's abdomen, pacemakers, etc.
6. Patients should refrain from using this instrument unless they have obtained permission from their doctor.
7. Ensure that the instrument is applied on moist skin and avoid the operation on dry skin.
8. Power off the instrument when not in use to ensure its safety.
9. To maintain hygiene and extend the device's service life, sterilize it with saline.

Contraindications

Before using this instrument, it is essential to consult your doctor or seek advice from professionals.


Details below:

1. Pregnant or lactating woman.
2. Individuals with hypertension, heart disease, or a pacemaker installed.
3. Individuals with unhealed wounds or those in the recovery period from surgery.
4. Individuals suffering from epilepsy, severe diabetes, and hyperthyroidism.
5. Individuals suffering from malignant tumors, hemophilia, and hemorrhage.
6. Individuals with skin disease or infectious disease.
7. For your safety, refrain from using if you've had plastic surgery or have implants/metal objects near the treatment area within 3 months.

Overview


This 4-in-1 beauty machine harnesses cutting-edge technology for thorough cleansing and skin brightening, ensuring users achieve a radiant and vibrant complexion. It provides the perfect solution for people seeking improvements in skin appearance without injections, drugs, or side effects. This external-use device delivers noticeable improvements day by day.

Interface Introduction

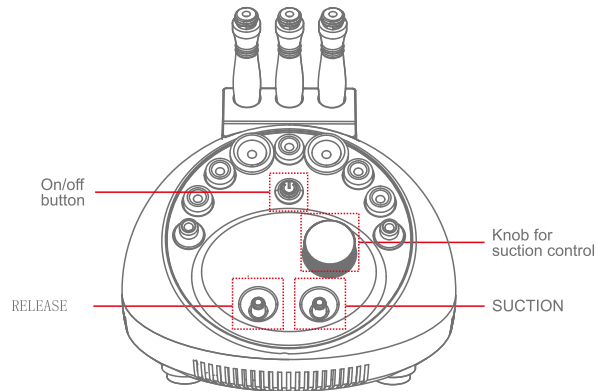

Usage Details

Scrubber:

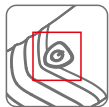

Power and Mode Switch button


3 modes available:

Mode 1: Cleansing mode


Mode 2: Introducing Mode

Mode 3: Lifting mode


Microdermabrasion:

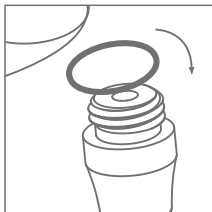

The black rubber ring must be correctly installed to prevent issues like no suction or weak suction.

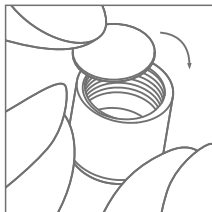

Installing the cotton filter is essential as it helps prevent foreign particles from damaging the device and also indicates the cleansing effectiveness.

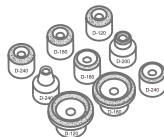

Select the appropriate working tip based on the size of the treatment area.



Recommended range of application: Nose wings, cheeks, and body




Recommended range of application: Forehead, chin, nose wings, and cheeks



Recommended range of application: blackhead and mite removal in the nasal folds, unclogging pores.

Transparent glass tube:

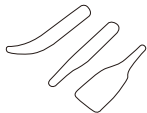

There are 3 transparent glass applicators. Insert one end of the hose into the suction hole on the main machine and the other end into the handle. The applicator starts drawing in while blocking the hole with your thumb and stops drawing in while unblocking the hole.




The applicator indicated by an arrow is used for cleansing blackheads in small treatment areas.



The applicator indicated by an arrow is used for cleansing clogged pores in small treatment areas.



The applicator indicated by an arrow is used for cleansing smaller areas, such as the cheeks.

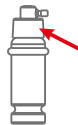


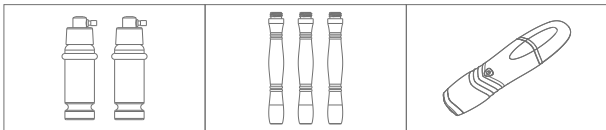
Unscrew the black cap and pour the liquid into the white bottle. The liquid can be pure water, toner, saline, watery serum, or any other suitable option based on your skin condition. It is recommended to fill the bottle halfway.



Insert one end of the hose into the release hole of the main unit and the other end into the location pointed by the arrow. Avoid inserting the hose into the suction hole, as this may cause damage to the machine.

Technical Specification

	Power supply Input: AC110V-240v 50HZ-60Hz	Net Weight : 0.9KG
	Power : 3W	Gross Weight: 1.3KG
	Dimensions: 6.7×6.7×4.7inches	
	Input Power: 3W	



Use Precautions

1. Strictly follow the instructions in the manual to use the instrument and train operator.
2. Before using the instrument, use a plug with a grounding pin and ensure the power socket is grounded.
3. Make sure the voltage of the instrument is adapted. If the local power supply voltage is unstable, we recommend adding a power-matched regulated power supply between the mains and the instrument.
4. The microdermabrasion tip must be connected to the suction hole (not the release hole) to prevent the machine from functioning abnormally.
5. The rubber ring must be correctly installed to prevent issues like no suction or weak suction.
6. The cotton filter is essential to prevent impurities from entering the tube and damaging the machine.
7. Keeping the skin moist contributes to a more comfortable experience.

Maintenance & Cleaning

1. Disconnect the power cord when not in use and store the instrument in a dry, cool, and well-ventilated place.
2. After use, rinse the microdermabrasion tip in saline, clean the transparent tube with a brush, allow them to air dry, and store them in a dry condition.
3. After cleaning the instrument, cover it with a dustcloth or plastic wrap. Before use, ensure to sanitize the handle again for hygienic reasons.
4. Disassembling the instrument and attempting to clean the inside is strictly prohibited. This task should only be performed by professionals.

FAQ

1. Will the skin get thinner and more sensitive after the microdermabrasion treatment?

A: No, the microdermabrasion treatment is a non-invasive procedure that targets only the aged stratum corneum. Since our skin naturally renews itself every 28 days, microdermabrasion contributes to smoother, brighter, and finer skin without causing increased sensitivity.

2. Why do I need exfoliation?

A: The stratum corneum is the outermost layer of the epidermis, consisting of 15 – 20 layers of dead cells that can vary from different parts of the human body. The skin's metabolism cycle also varies among individuals of different ages. For example, the cycle is 28 days for teenagers, whereas for adults, it is the sum of their age plus 10. However, for a variety of reasons, the stratum corneum might not shed promptly, resulting in an increased thickness of the stratum corneum, a dull complexion, and an impact on skin metabolism. In such cases, additional exfoliation is necessary.

3.Can I apply makeup immediately after a microdermabrasion treatment?

A: Yes, the treatment contributes to achieving moisturized and plump skin, ensuring a seamless and well-blended makeup application. Microdermabrasion is a non-invasive procedure that delivers fast, safe, and healthy results.


4.When can I expect to see results?

A: Immediately after the initial session, you will observe a noticeable improvement in your skin's clarity and transparency. Multi-sessions contribute to more visible effects including deeply hydrated skin, brightened complexion, reduced fine lines, and prevention of the aging process.

5.Is the microdermabrasion addicting?

A: No. Microdermabrasion is a skincare procedure that enhances skin condition and provides deep cleansing. It goes beyond mere treatment; it helps the skin absorb nutrients and addresses concerns like dullness, dehydration, rough texture, and wrinkles. Discontinuing it will simply allow the skin to return to its original state and aging process without any addictive effects.

Ask & Answer

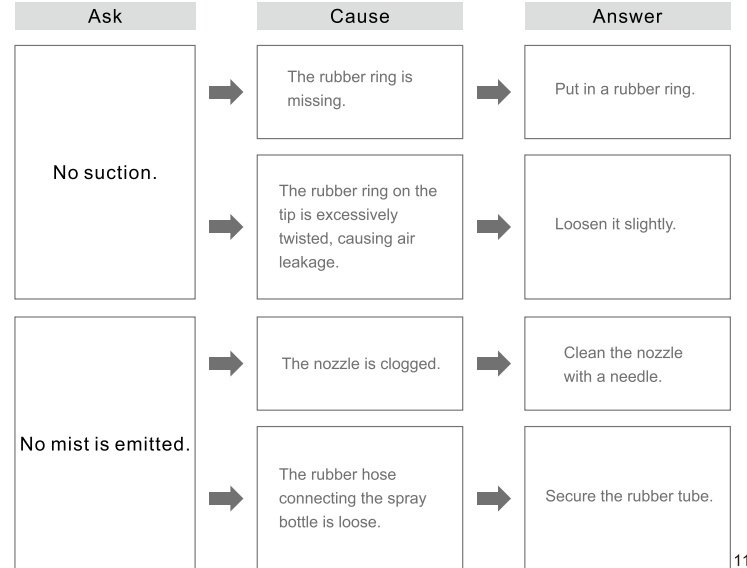

Operation Guide

Face

Microdermabrasion (Step1:  3mins > Step2:  2mins)

Adjust the suction intensity according to individual endurance.

Product required: Wet cotton pad



Sprayer (> Step4:  10mins)
Product required: Serum + Saline (1: 10)



Scrubber (> Step3:  2mins)

Recommended Mode: Cleansing mode

Product required: Wet cotton pad



Warranty & Support

We provide a one-year warranty and a three-month return and replacement policy from the date of purchase.

Product Name: Cleansing & Hydrating Machine
Product Model: MS-33P3N
Production Date: _____
Quality Inspection: _____

