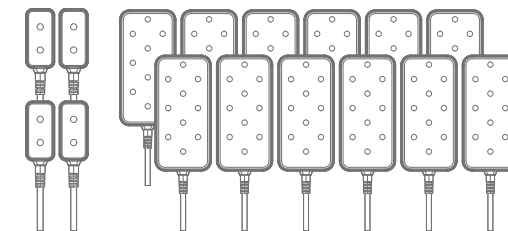
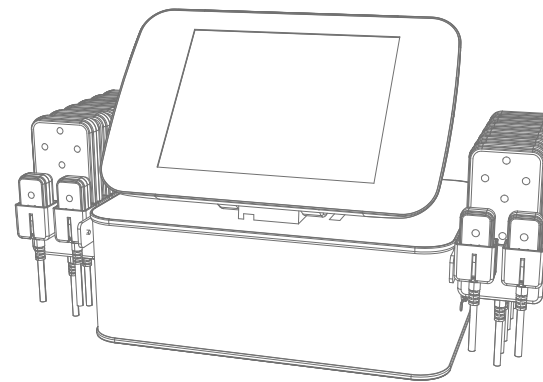


User Manual

Home SPA Instrument

MY-16151



Thank you for purchasing our product.

- Please read this User's Manual carefully for information on how to use the device properly and safely.
- After reading the contents hereof, be sure to keep this User's Manual in an accessible location.
- Make sure to include this User's Manual if you are transferring ownership of the device to another person.

Contents

Part I

Safety precautions	1
Contraindications	2
Overview	2
Interface function	3
Technical specification	4
Use Precautions	5
Maintenance & Cleaning	5
FAQ	6
Ask & Answer	7
Operation guide	8
Warranty & Support	9

Safety precautions

- 1.Please read this manual carefully before use and familiarize yourself with the operation of this instrument and related technology.
- 2.This instrument is a high-precision electronic device that requires proper storage. Please avoid storing it in humid conditions, direct sunlight, or near flammables and explosives to prevent damage.
- 3.Always use the designated accessory provided or recommended by the manufacturer to ensure performance throughout its service life.
- 4.Do not place the instrument in high-temperature conditions, as this may negatively impact its service life and performance.
- 5.Do not apply the instrument on eyes, thyroid, parathyroid, testicles, pregnant women's abdomen, pacemakers, etc.
- 6.Patients should refrain from using this instrument unless they have obtained permission from their doctor.
- 7.Power off the instrument when not in use to ensure its safety.
- 8.If you are currently taking weight-loss drugs, please discontinue their use 1 - 2 months before using this instrument. Failure to do so may result in an extended course of treatment.
- 9.To maintain hygiene and extend the instrument's service life, sterilize it with saline.

Contraindications

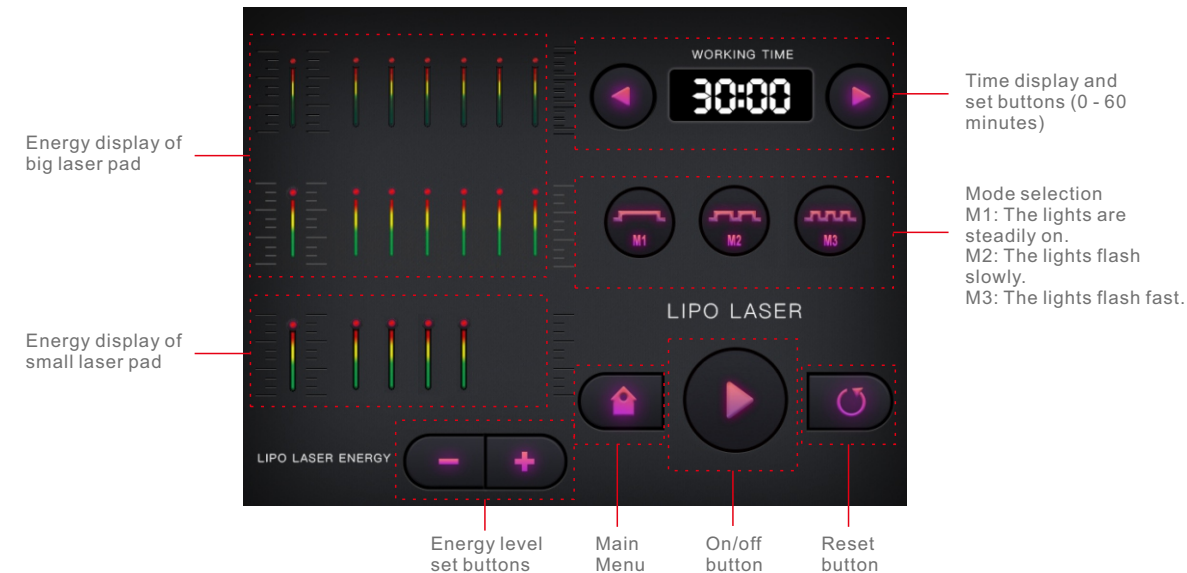
Before using this instrument, it is essential to consult your doctor or seek advice from professionals.
Details as below:

1. Pregnant or lactating woman.
2. Individuals with hypertension, heart disease, or a pacemaker installed.
3. Individuals with unhealed wounds or those in the recovery period from surgery.
4. Individuals suffering from epilepsy, severe diabetes, and hyperthyroidism.
5. Individuals suffering from malignant tumors, hemophilia, and hemorrhage.
6. Individuals with skin disease or infectious disease.
7. Individuals have had a cesarean section within six months, a miscarriage within three months, or a natural birth within two months.
8. For your safety, refrain from using if you've had plastic surgery or have implants/metal objects near the treatment area within 3 months.


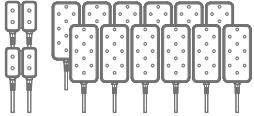
Overview

With the most popular and safest technology for anti-aging, firming skin, and body shaping, this LED beauty instrument ensures you achieve your desired body figure and radiant complexion. It provides the perfect solution for people seeking improvements in skin appearance and body figures without injections, drugs, or side effects. For external use with immediate effects.

Interface function



Technical specification

	Power supply : 110V-220V 50HZ/60HZ	Net Weight: 9.5KG
	Power: 40W	Gross Weight:10.6KG
	Dimensions: 14.6 × 9.84 × 13.4 inches	
	Input Power: 40W	
	Wavelength: 160MW	
	Power: 40W	

Use Precautions

- 1.Strictly follow the instructions in the manual to use the instrument and train operator.
- 2.Before using the instrument, use a plug with a grounding pin and ensure the power socket is grounded.
- 3.Make sure the voltage of the instrument is adapted. If the local power supply voltage is unstable, we recommend adding a power-matched regulated power supply between the mains and the instrument.
- 4.Remove all the metal objects inside the body to avoid causing accidents or affecting treatment effects.
- 5.Stop the intake of weight drugs 1 - 2 months before the treatment, or extend the treatment course if you want to take the treatment immediately.
- 6.Do not use the instrument on an empty abdomen, and take the treatment 1 hour later after the meal.

Maintenance & Cleaning

- 1.Disconnect the power cord when not in use and store the instrument in a dry, cool, and well-ventilated place.
- 2.Clean and sanitize the handle's surface with a damp cloth or saline, and then dry it thoroughly with a dry towel to ensure it remains dry.
- 3.After cleaning the instrument, cover it with a dustcloth or plastic wrap. Before use, ensure to sanitize the handle again for hygienic reasons.
- 4.Disassembling the instrument and attempting to clean the inside is strictly prohibited. This task should only be performed by professionals.

FAQ

1.LED fat removal treatment or liposuction surgery?

A: Compared to traditional liposuction, LED treatment scans the fat area one by one, without causing unevenness in the skin, and conversely, leaving the skin smooth and tight without sagging. In addition, the LED energy used in LED fat removal treatment only works on the set fat layer, not hurting the surface skin while accelerating the self-repair of skin elastic fibers, leaving the skin of the treatment area tight and smooth.

2.Will I regain weight after LED fat removal treatment?

A: The number of fat cells in the body is fixed and remains constant as the age grows to adulthood. The increase in body weight is due to the increase in the volume of fat cells. LED treatment removes the local fat structure and reduces the volume of fat, however, the fact that weight regaining is rare after the treatment doesn't mean it will never occur. Additionally, maintaining a healthy diet would greatly contribute to preventing it.

3.Do I need to diet and exercise during LED fat removal treatments?

A: The effect of LED treatment is significant, but we cannot only rely on it. To ensure lasting results, maintaining a healthy diet after reaching your desired weight post-treatment is crucial, as failure to do so may lead to weight regain. It is important to know that no weight loss product in the world guarantees permanent effect. So, the LED treatment is an auxiliary weight loss method. When the weight reaches the ideal range, you need to control your diet and take regular exercise after the treatment.

Ask & Answer

Ask		Cause		Answer
Fail to turn on the instrument.	→	The power cord is disconnected to the wall socket.	→	Check if the power cord is well connected to the wall socket and ensure that the socket is functional.
	→	The fuse is loose or burned.	→	Replace the fuse on the back of the instrument and restart it.
The instrument is not working and can't be restarted?	→	The fuse may have burned out due to unstable voltage or prolonged use.	→	Please replace the fuse and restart the instrument 5 minutes later.
The screen shows an error message.	→	Data error?	→	1. Turn off the instrument and unplug the power cord, and restart the instrument 1 minutes later. 2.If the above method doesn't work, please contact with the dealer for a replacement.

Operational guidance

Face



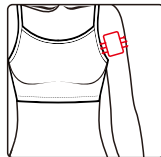
20mins



Arm




20mins



Belly



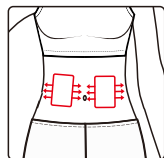
20mins



Leg



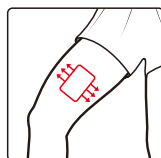20mins



Warranty & Support

We offer a one-year warranty on the main unit and a three-month warranty on components.

Product Name: LED Fat Remover
Product Model: MY-16151
Production Date: _____
Quality Inspection: _____

PASS