

# LED red light physiotherapy belt

## Product introduction:

- 1) voltage input: DC 5V 2A
- 2) charging port: Type-C
- 3) LED power: 3W
- 4) material: PU+fabric
- 5) LED quantity: 15pcs

## Matters needing attention:

- 1) After use, it should be stored in a ventilated and dry place;
- 2) Do not press the product heavily, otherwise it will be damaged;
- 3) Do not pull out the USB charging port too hard, otherwise it will break;
- 4) To clean the product, turn off the power, wipe it with a damp cloth, and dry it with a towel;
- 5) Do not soak or rinse in water, and do not use strong detergents to clean;

Picture  
breakdown:






Switch / time setting key

After the power is connected, press and hold the power button for 3 seconds to turn on and off the machine;


DC 5V power input interface (Type-C)



Gree

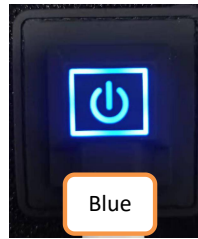

Blue

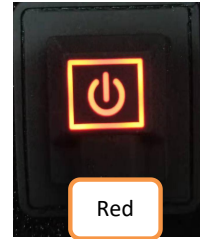

Red

Description of key distinction:

The "open key" is divided into 3 gears: green, blue, red

The three gears are controlled by timing and time distribution:

Green: 15 minutes

Blue: 30 minutes