

ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

Model: SR-AT-04



Dear Customer,thank you for buying our product,please read the manual instructions carefully and operate it correctly.The wrong usage methods would cause the instrument and its accessories damaged or faulted.After usage,please turn off the power and maintain the accessories.

ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

BASIC INFORMATION

- 1.Product Name: Ultrasonic Beauty Instrument
- 2.Therapy Probe: 3 pieces
- 3.Applicable Parts: Body、 Face skin、 Eyes skin
- 4.Main Functions: Skin Care、 Thin Body
- 5.Applicable occasions: Beauty Salons、 SPA、 Home Use

Special Announcement: In order to use the instrument safely and effectively, please read this manual carefully before using it. All beauty instruments are contraindicated for the following groups:


- Heart patients
- Patients with severe hypertension
- People with metal in their bodies
- Pregnant and lactating women、 Menstrual period
- Underage child
- A variety of infectious diseases

ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

INTRODUCTION

1.Panel

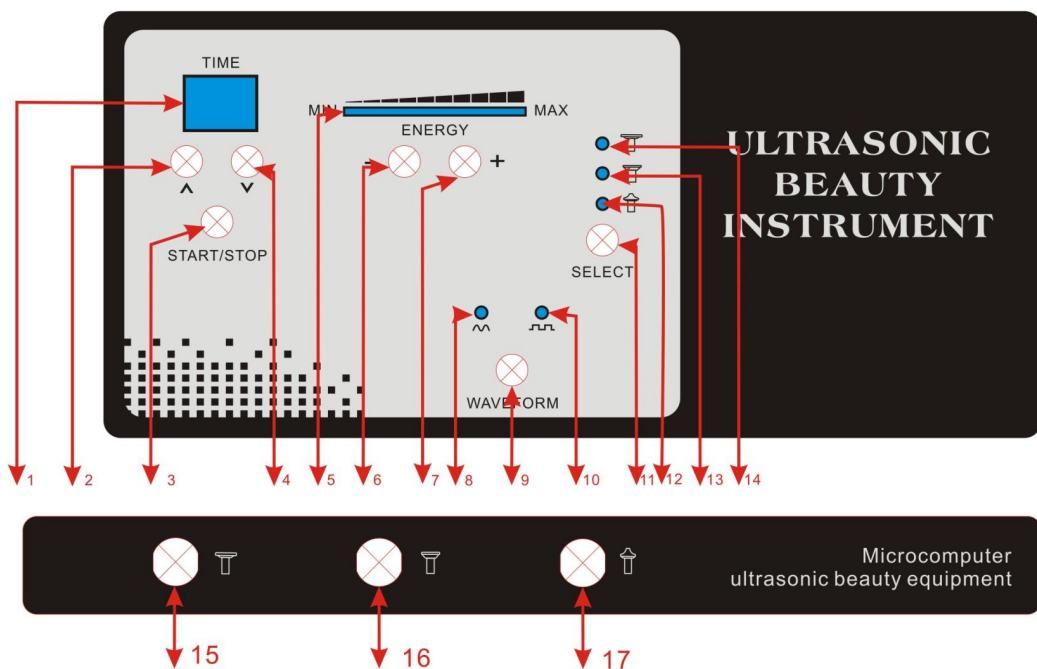

1:Timer

2&4: Time UP/DOWN

3: Power START/STOP

5: ENERGY show

6&7:ENERGY-/+

8: Continuty Wave

9: WAVEFORM Selection

10: Discontinuty Wave

11: Probes selection

12: Eyes probe

13: Face probe

ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

14: Body probe

15: Output of Body Probe

16: Output of Face Probe

17: Output of Face Probe

2.Probes

- Probe for Eyes Skin



ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

Cylindrical design for the handle, which helps to better adhere to the eyes skin. Usage with eye cream or eye essence, it can effectively eliminate dark circles, eye bags, fine lines, dry lines, firm eyes skin.

- Probe for Face Skin



To rub the essence on the skin, then massage the face skin with the probe, which can further penetrate the essence into skin so that it can be fully treated. With the usage of essence, it can effectively shrink skin

ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

pores, fade fleck, improve dark-yellow skin and remove acne, brighten the skin.

- Probe for Body



ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

To rub essence on the body skin, then massage the body with this probe, which can further penetrate the essence into the body skin so that it can be fully treated. With the usage of essence, it can effectively remove stretch marks, burn redundant fatness of arms, abdomen, hips and legs so that to tighten body skin and thin body.

STEPS FOR USAGE

- 1、 Before using the instrument, clean the treatment site of skin.
- 2、 After the skin is cleaned, use a thermal sprayer to steam the treatment site for about 10 minutes to make the pores open slightly.
- 3、 Remove facial blackheads.
- 4、 Tap the skin with a proper amount of toner.
- 5、 Gently wipe off the toner.
- 6、 Select the appropriate cream to rub evenly on the treatment site.
- 7、 Massage with the probe on the skin rubbed ointment.
- 8、 After the treatment is finished, clean the probes.

ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

PRECAUTIONS

- 1、 Wash facial skin or treatment area before ultrasound treatment.
- 2、 The plaster used should have a certain viscosity,apply enough plaster.
- 3、 The ultrasonic probe should not be close to the eyes.
- 4、 The time for using body ultrasonic probe is 30 minutes per time.
- 5、 The time for using Face probe is 20 minutes per time.
- 6、 The time for using Eye probe is 15 minutes per time.
- 7、 It is not allowed to without cosmetic operation while the ultrasonic energy is outputting,otherwise,the the probes would be burned because of overheat.

USAGE GUIDE

Ultrasonic Lead-Out Detoxification

1. Do the facial detoxification,use with detoxifying capsules or other cleansing lotions.

ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

2. Usage Frequency: 1 time each week

3. The recommended Usage Time:

- The face probe is used 15 minutes per time.
- Eye probe is used 10 minutes per time.

❖ **It is best to shorten the usage time for the sensitive skin、red blood or thin skin.**

Ultrasonic Lead-In Therapy

1.Used with beauty lotions or essence to help skin absorb

2.Usage Frequency: 2-3 times each week

3.The recommend Usage Time:

- Each time 10 minutes for the eyes
- Each time 15 minute for the face.
- Each time 20 minutes for a certain part of body.

(include arms、abdomen、hips、Legs)

❖ **Before doing any therapy on the sensitive skin,make an allergy**

ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

**test by daubing the cream after the ear.The instrument itself
does not have any allergens.**

TECHNICAL PARAMETER

Rated Voltage	AC220V/50Hz or AC110V/60Hz
Rated Power	≤50W
Working Environment	Environment Temp.: ≤-10℃~40℃
	Relative Humidity: ≤80%
Maximum Time	90 minutes
Ultrasonic Energy	1~10 levels,adjustable
Ultrasonic Frequency	1MHz
Penetration Depth	Body Probe: 0.75mm~1.25mm
	Face/Eyes Probes : 0.25mm~0.5mm
Gross.Weight	6 KG (All accessories included)
Packing Dimension	46cmx32cmx19cm

ULTRASONIC BEAUTY INSTRUMENT

USER MANUAL

PACKING LIST

ITEM	QTY	REMARKS
Controller Machine	1(pcs)	
Accessories Holder	1 (pcs)	
BODY Ultrasonic Probe	1(pcs)	
FACE Ultrasonic Probe	1(pcs)	
EYES Ultrasonic Probe	1(pcs)	
Power Cord	1(pcs)	
English Manual	1(copy)	