



SR-AF1312

Instructions

Warning

Dear users:

Thanks for choosing our company's product, in order to fully make use of the products, we sincerely suggest follow:

- ◆◆ Please read the instruction in detail and keep it well for further reading and consulting.
- ◆◆ Please follow the instruction to fit and operate the machine correctly.
- ◆◆ Please do not remove and change any accessories of the machine.
- ◆◆ It forbids others to open and remove the mainframe except the authorized people in the company.
- ◆◆ It forbids others to change and maintain the machine except the professional people.
- ◆◆ Please contact us if there is any trouble during operation.



Please read the NOTES carefully before use Machine.

1. Testing: Please turn on the machine first before connect all the handpiece to see if the machine works.
2. Please turn on the machine, then connect the sprayer
3. Please install the pure water bottle and waster water bottle in correct place according to the mark
4. Please donot connect sprayer when using hydro microdermabrasion and diamond dermabrasion.

After Sale service

◆ We offer ONE YEAR warranty. Within warranty period, all exchanges or replacement due to defects are free. (Buyer is responsible for shipping charges). If the machine is over warranty, please also do not hesitate to contact us. we will assure you of our best service all the time.

目录

| | |
|--------------------------------|----|
| I. Working Condition..... | 4 |
| II. Introduction..... | 4 |
| HYDRA/HYDRO-DERMABREASION..... | 4 |
| ULTRASNOIC..... | 4 |
| RADIO FREQUENCY..... | 4 |
| BIO MICRO CURRENT..... | 5 |
| HYDRO OXYGEN JET SPRAY..... | 5 |
| LED PDT PHOTON..... | 5 |
| III. Accessories..... | 6 |
| IV. Operation..... | 7 |
| ULTRASOUND..... | 8 |
| BIO MICRO CURRENT..... | 9 |
| SEXTUPOLAR RF..... | 9 |
| COOLING HEAD..... | 10 |
| HYDRO SPRAY..... | 11 |
| HYDRA/HYDRO-DERMABREASION..... | 12 |
| LED PHOTON MASK..... | 15 |
| V. After using..... | 15 |

I. Working Condition

Using at normal temperatures and pressure. Please do not stay in too cold, too hot, too dry and too wet.

II. Introduction

HYDRA/HYDRO-DERMABREASION

Hydro-Microdermabrasion completely changed the traditional way, which is hand-cleaning skin relying on individual's practice skills, Hydro-Microdermabrasion uses vacuum suction mode controlled by intelligent process, through the combination of products and equipment to improve skin texture.

It uses specifically designed hydropeel tips which gently exfoliate the skin using a derma planing motion. The spiral tips allow the SKIN SERUMS longer duration on the skin, whilst the spiral edges are designed to push the serums deeper into the skin – producing a plumped up effect!

Hydro-Microdermabrasion resurfacing treatment thoroughly enhances the skin using vortex technology to cleanse, exfoliate, extract, and hydrate simultaneously. It merges soothing and invigorating spa therapies with advanced medical technology to achieve instant lasting results. The procedure is smoothing, moisturizing, non-irritating, and immediately effective.

ULTRASNOIC

What is ultrasonic massage? Ultrasonic use of high frequency ultrasonic wave (1000000 / 3000000) shocks to promote human blood circulation and speed up metabolism. Application ultrasound, with all kinds of cream cosmetics,penetrate into the skin to strengthen the organization to beauty effect.

RADIO FREQUENCY

Radio frequency is also called diathermy(deep heating) which is the system to provide the treatment by generating the heat from inside the human body. Wrinkles and loose skin can take its toll on your looks. Even healthy and active people eventually show signs of aging on their faces. Since it's the most important visual feature of your body—the one that most people will recognize you by—it's important to keep your face fresh and youthful. More and more people have tried face lifts over the years to combat wrinkles and skin imperfections. Though often successful, traditional face lifts require minor surgery and prolonged recovery periods. There's increasing demand for non-surgical techniques to refresh facial appearances. That's where radio frequency face lifts come into play.

A radio frequency emitter is pressed up against the skin. These radio waves get past the outer skin layers and deliver heat energy to the muscles and tissues beneath. The

heat helps to contract these layers and build up collagen levels. The overall result tightens up the outer layers of skin and reduces wrinkling effects. Since there's a lot of heat involved, some cooling must be applied to the skin at the same time.

The radio frequency face lift was a non-surgical treatment for wrinkles and facial skin imperfections. It is a proven medical procedure that requires no scalpels or sutures. It's also a popular choice for busy professionals who care about their looks. The treatment only takes about an hour, and full recovery is expected in a few days. The results take a bit of time to show fully though. Some results are visible immediately, while the full results take a few months to develop as the deep tissue layers heal.

BIO MICRO CURRENT

By outputting simulated human bio-current, it can go through the skin into muscle cell, stipulating energy ATP contained inside the cell, enable cell to regain normal operation and funtion. The unit is application fo laest technology with preset computer program that can help for better shaping for lineament, facial sharpening, double layers chin, wrinkle, crow's-feet, pouches, black eye, etc. boosting up facial blood circulation and metastasis, reducing facial splash and pore constriction. BIO will tighten and refine the facial skin int sanguine condition, which have proofing effect to loosen skin, reaching the purpose of beauty enhancement.

HYDRO OXYGEN JET SPRAY

Nutrition or cosmetic product can be filled into the stainless steel bottle. With high pressure, the nutrition and oxygen are sprayed on the skin area which greatly enhance nutrition to be absorbed maximally, skin cleaning and skin rejuvenation.

LED PDT PHOTON

How the lights work?

LED light is the use of specific wavelengths of light which give off energy that are used in many applications,this is known as photon . LED stand for light emitting diode, diodes give off energy in the form of light,powerful,bright,efficient,and are used to kill the bacteria that cause spots

How the skin LED works

LED blue light helps to kill the bacterial that cause spots,reduce inflammation and provide a general rejuvenation effect on the skin.Spots begin when the pores become clogged or blocked and bacteria start to grow,resulting in every teenager's nightmare spots.

Blue light (415nm)

Bacteria in acne contain porphyrins, the blue light is able to kill porphyrins due to the high affinity with this wave length.The blue light furthermore has calming action which is very effective for hypersensitivity

Green light(525-550nm)

The balance of the green light can balance color pigment,reduce fine lines, nutrient aging skin,speed up the healing process of the wound,and lighten the scare

Red light(630nm)

Will stimulate the production of collagen,collagen is an essential protein used to repair damaged tissue and to replace old tissue,best for fine lines, and will also help to reduce large pores

Purple light

It is red and blue dual -band light, it is a combination of two kinds of phototherapy effect.

Yellow light

Add energy to skin cells promote gland function, assist digestion, treatment of skin diseases , enhance the immune function.

LED light (white light)

Penetrate the skin deep, speeds up the living tissue of chen xin metabolism, decompose splash improve the appearance of fine lines and sagging skin.

ICE REPAIR HANDLE

Function: Froze skin, shrink pores.

Tips: No vibration, No light. During use, there is a slight heat in the handle portion.

VANADIUM-TITANIUM

Digesting forehead wrinkles, improving the gloss of forehead skin, improving the fine wrinkles, micro spots and dryness of eye corners; Remodeling skin plump and brighten the skin; Shrinking pores; Improving sagging jowls, skin flabby, dry; Diluting the neck wrinkles.Note: The VANADIUM-TITANIUM needs to be used with a 5ml syringe.

VIBRATION MASSAGER

V face and model body;Tightening the tender skin; Relieving nerves and fatigue; Shrinking pores; Diluting edict wrinkle and eyewrinkles, can also use waist massage; Warming uterus and strengthen daily uterus care.

3.0MM RADAR SHAPING HANDLE

It helps to remove wrinkles and improve skin elasticity and to form collagen. It makes the energy at 65°C to deliver to the SMAS layer (about 3.0mm~4.5mm) in the face skin and it contracts the skin muscle layer.

Function: Shaping “V” facial lines, firming skin.

Tips: For THINNER/ YOUNGER skin types. Need to be used with gel. No vibration. No light.

4.5MM RADAR SHAPING HANDLE

It helps to remove wrinkles and improve skin elasticity and to form collagen. It makes the energy at 65°C to deliver to the SMAS layer (about 3.0mm~4.5mm) in the face skin and it contracts the skin muscle layer.

Function: Shaping “V” facial lines, firming skin.

Tips: For THICKER/ ELDER skin types. Need to be used with gel. No vibration. No light.

III. Accessories

- Hydro microdermabrasion
- BIO
- Ultrasound probe
- Sextupolar RF
- Cooling head
- Hydro spray
- LED mask
- VANADIUM-TITANIUM
- 3.0 Radar Shaping Handle
- 4.5 Radar Shaping Handle
- Vibration Massager

IV. Operation

1. Install the acrylic holder and connect all accessories with machine.
2. Fill in water to the fresh water bottle.
3. Insert the power supply like below
4. Trun on the machine.



Multi-functional Skin Care System

Bubble
Brush

Radio
Frequently

Ultrasonic

Micro
Power

Ice Seal
Imports

Radar Lines
3.0mm

Oxygen
Gun

Radar Lines
4.5mm

Spectral
Mask

Vanadium
Titanium

Bubble Brush

Attention of hydro dermabrision:

- do not use on wound
- pour out the water in the waste recycling bottle in time, to avoid damage caused by overflowing.
- different size of heads can be flush with water directly, then soaked with 75% alcohol.

Forbidden group of hydrodermabrasion

- people who has fever, infections diseases, acute diseases
- people with hemorrhagic disease, trauma, vascular rupture, skin inflammation, skin disease

Press the Bubble Brush mode and set the intensity and treatment time.


Press start.



Radio Frequency

Before Treatment, please apply the gel on the face.
Press the radio frequency mode and set the intensity and treatment time.

Press start.

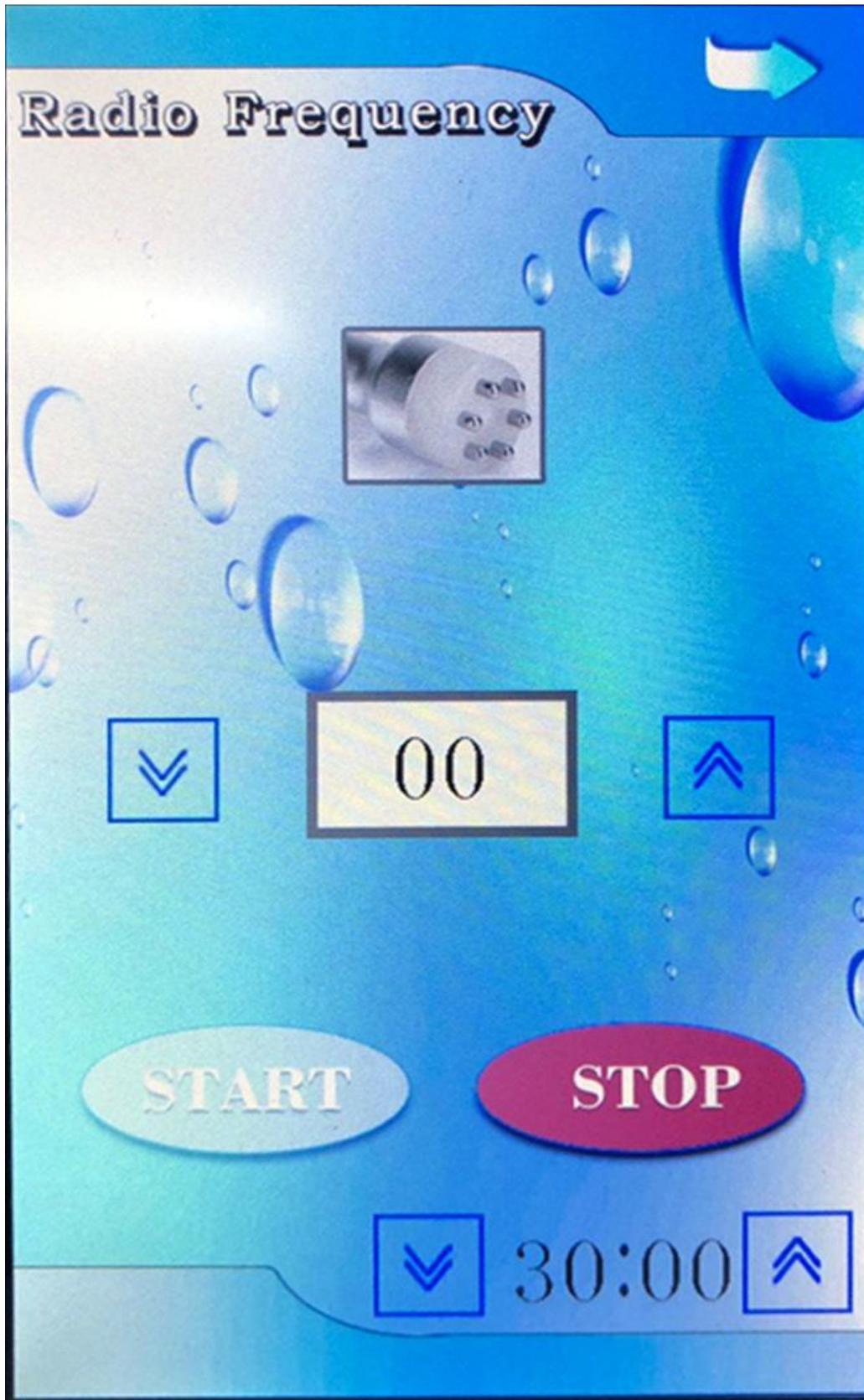

Ultrasonic

Before Treatment, please apply the gel on the face.
Press the ultrasonic mode and set the intensity and treatment time.
Press start.



Micro Power

Press the Micro mode and set the intensity and treatment time.


Press start.



Ice Seal Imports

Press the Ice Seal Imports mode and set the intensity and treatment time.


Press start.

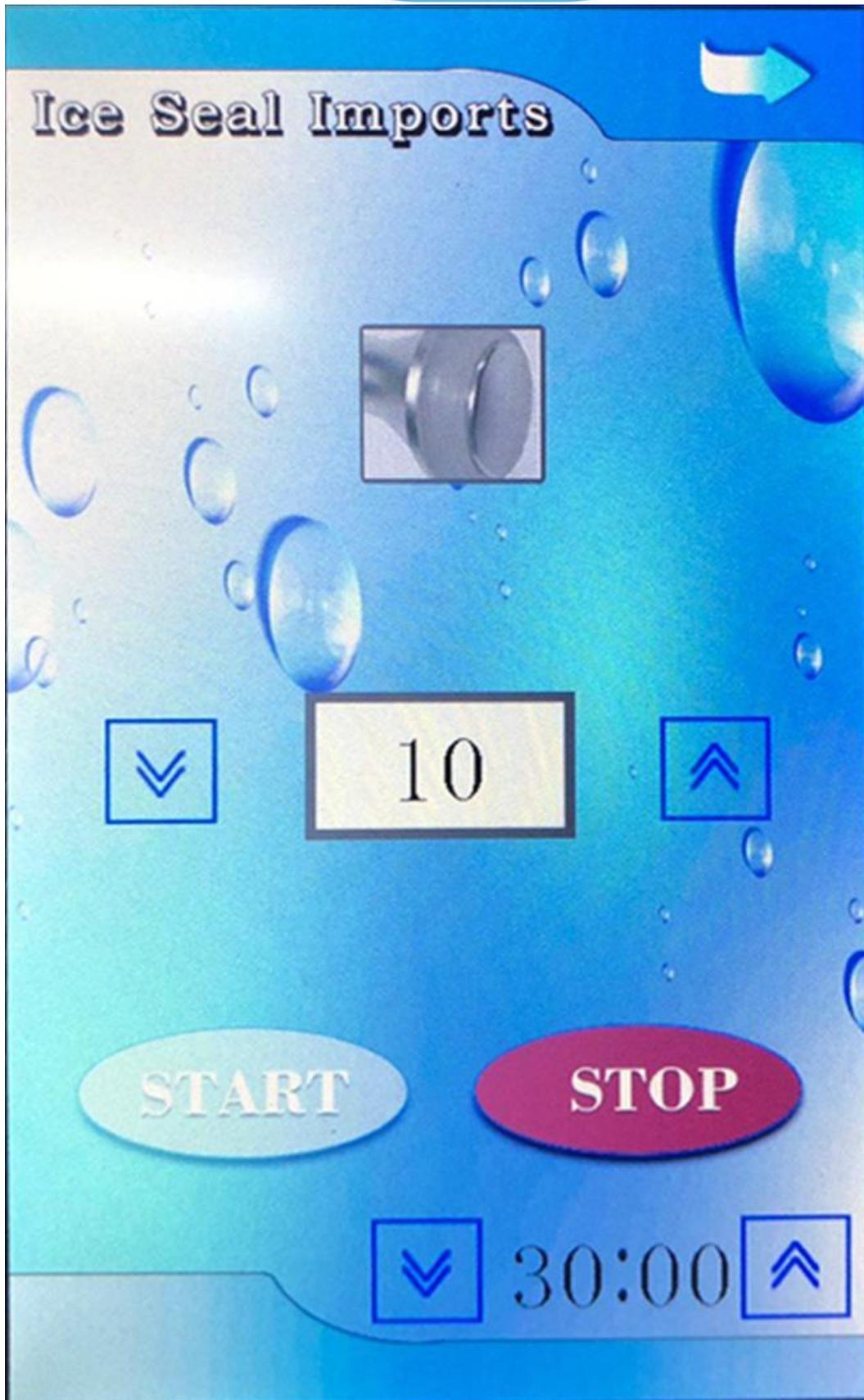

Radar Lines 3.0mm

Press the Radar handle mode and set the intensity, Output Time, Interval and treatment time.


Press start.

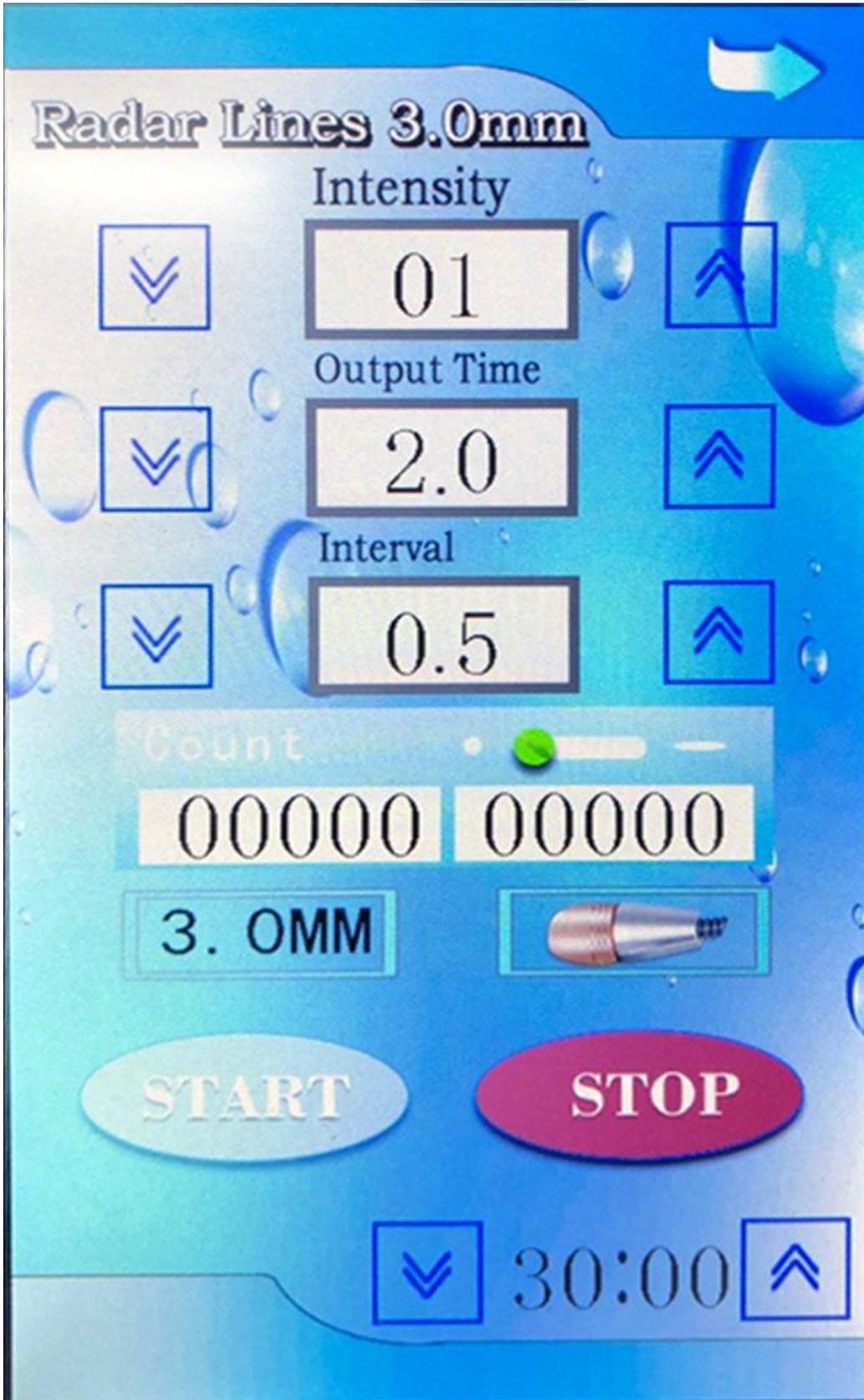

Oxygen Gun

Press the Oxygen Gun mode and set the intensity and treatment time.


Press start.

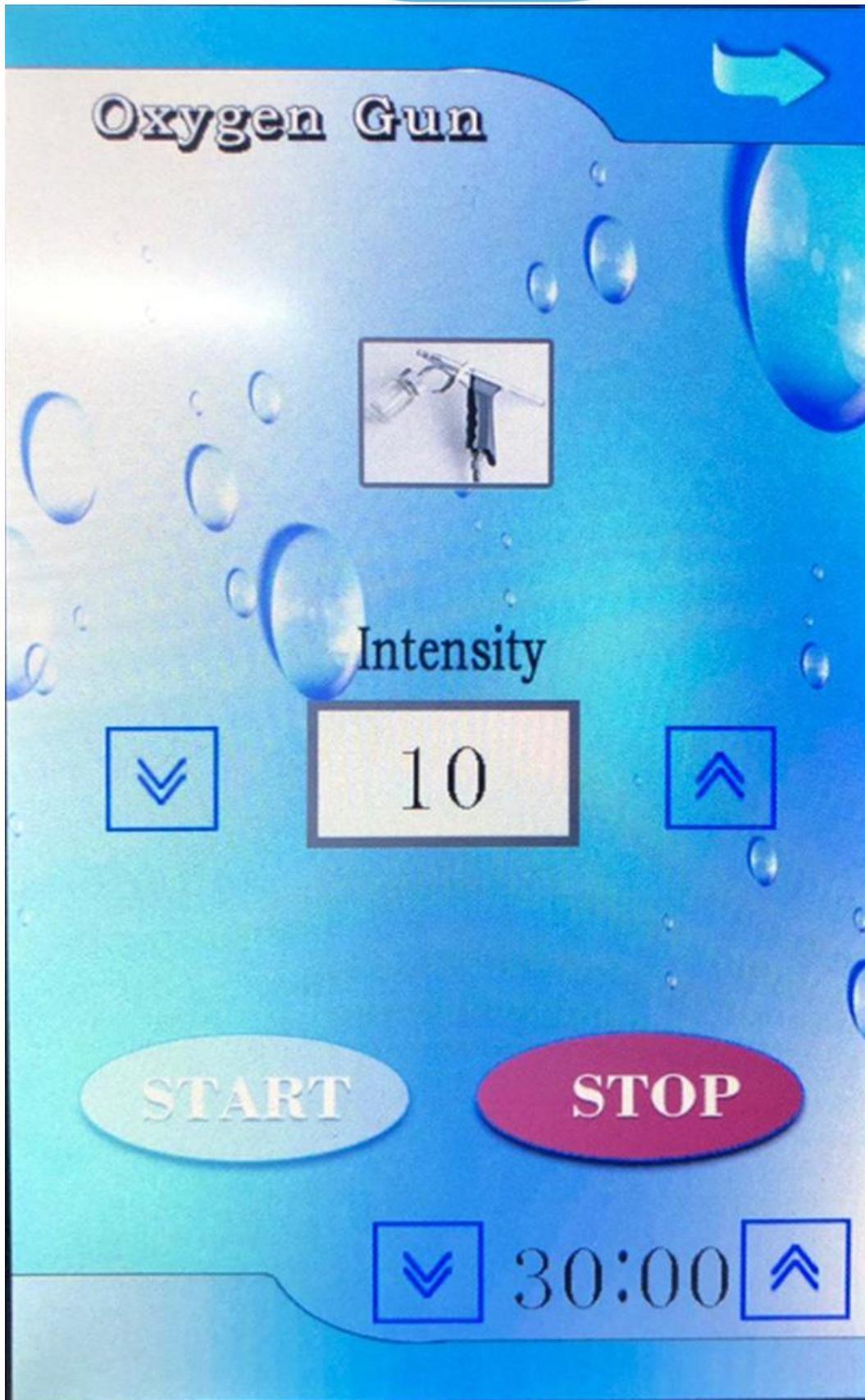

Radar Lines 4.5mm

Press the Radar handle mode and set the intensity, Output Time, Interval and treatment time.


Press start.

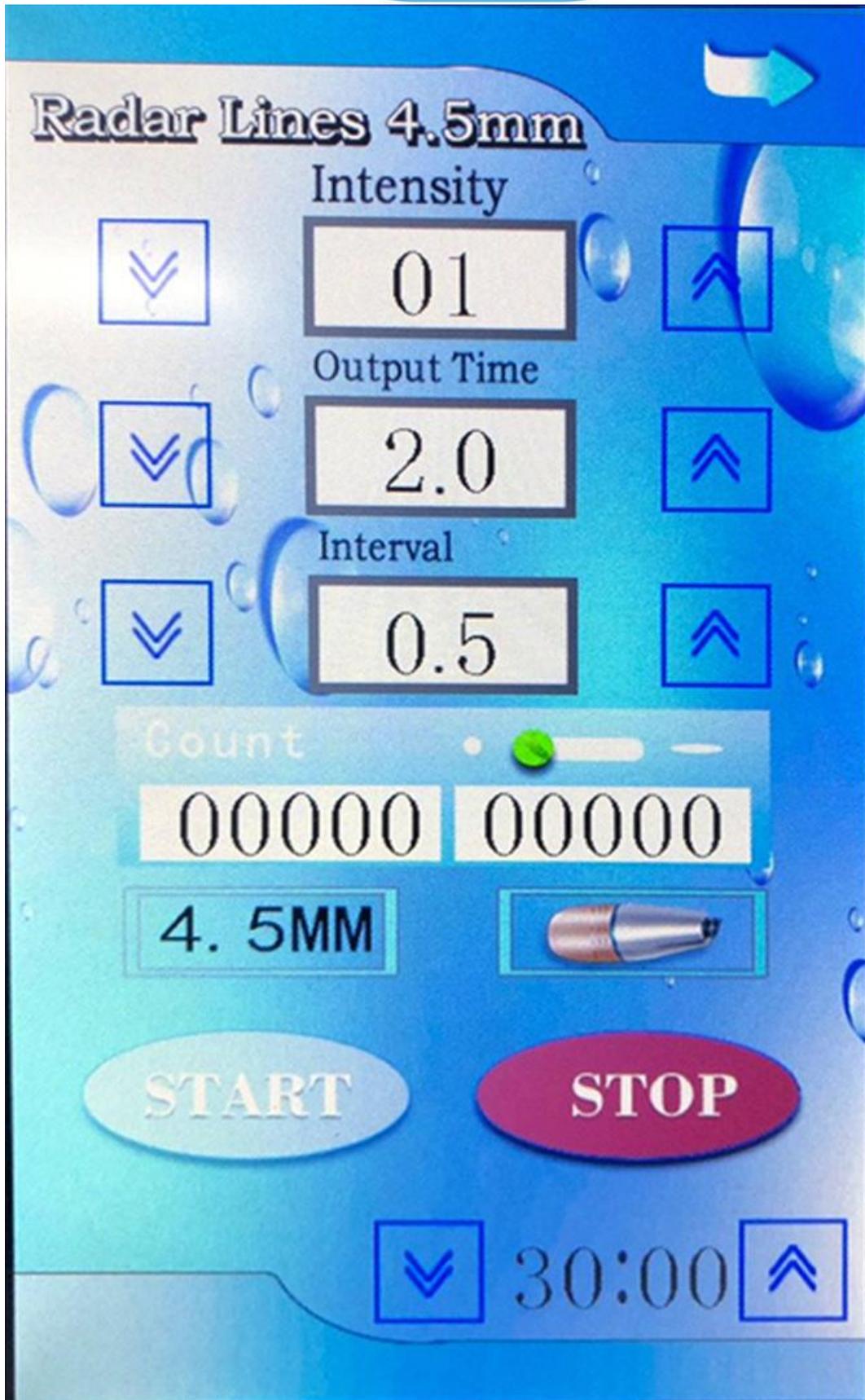

Spectral Mask

Press the Spectral Mask mode and set the intensity, Model and treatment time.

Press start.

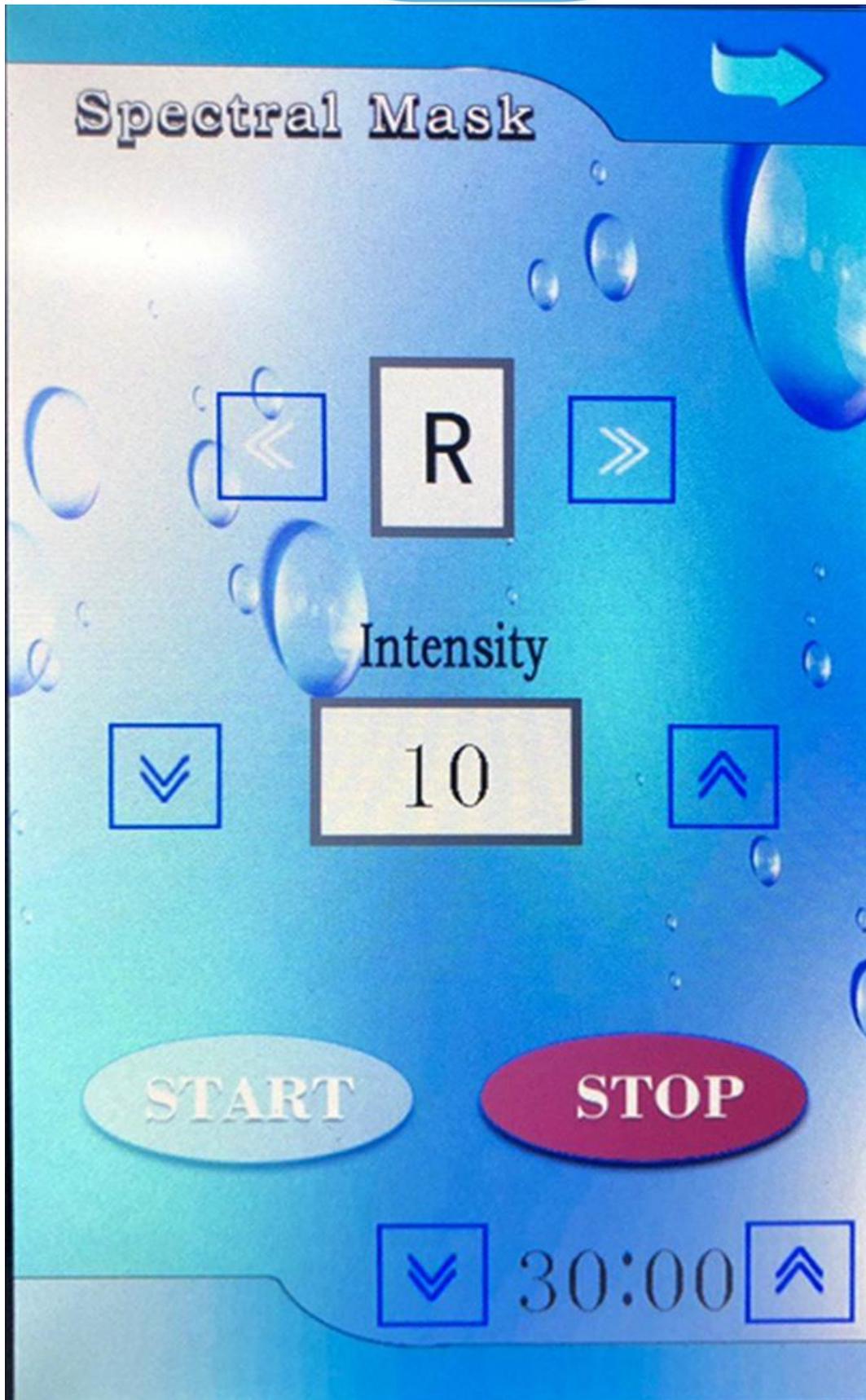

Vanadium Titanium

Press the Vanadium Titanium mode and set the intensity and treatment time.

Press start.

Vanadium Titanium

Water Output Speed



00



Water Inlet Speed



00



Propeller Energy

◀ Go Ahead

Go Back ▶

START

STOP



30:00



HYDRO SPRAY

Fill the bottle full of water.



Connect the bottle with the spray.




Set the parameter and treatment time, then press start and spray all over the face.

Please rotate the white switch to adjust the vacuum intensity.



Beautician hold the vacuum spa handle on skin for a while until the water flow into the work head, and then drop some water on skin by points, to make the skin surface wet, and then apply cleansing product on the skin, and move the vacuum work head on skin back and forth

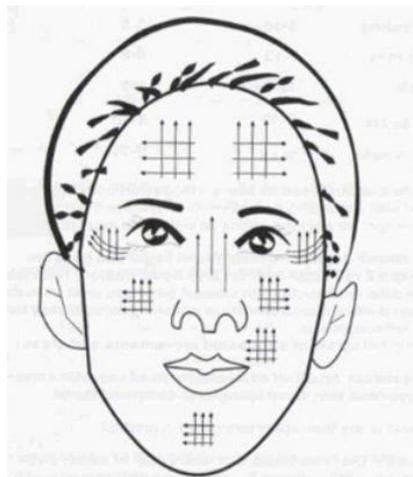

After treatment, please use one finger to press the hydro dermabrasion head, let the water flow into another bottle to wash the pipe. Open the drainage swith to let the water out.



LED PHOTON MASK

- a. Press and hold power button and you can turn on and off through it
- b. Mode button is pressed, there are setting mode and time in LCD
- c. Possible to time setting max 60 minutes using by time button
- d. Press the power button shortly and product begin to work
- e. Press and hold the power button to shut down LED output mode using

V. After using

1. Do a hydrating mask after the treatment.
2. Next treatment should be after half of month.