


# User Manual

MS-2184

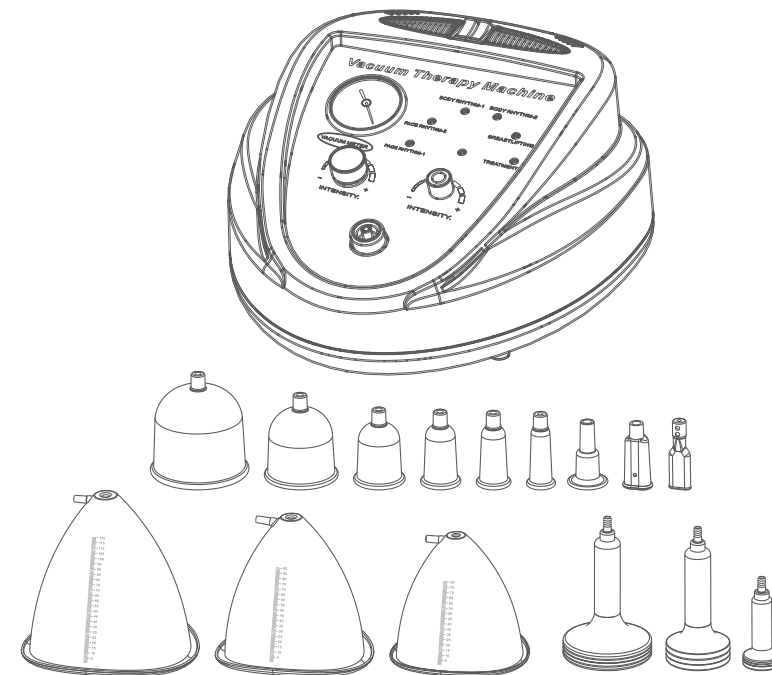

## Content

1.Safety Precautions .....	1
2.Contraindications .....	2
3.Overview .....	2
4.Control Panel Introduction .....	3
5.Technical specification .....	9
6.Use Precautions .....	9
7.Maintenance & Cleaning .....	10
8.FAQ .....	10
9.Ask & Answer .....	12
10.Operation guide .....	12
11.Warranty and support .....	15

## 1.Safety precautions

- 1.Please read this manual carefully before use and familiarize yourself with the operation of this instrument and related technology.
- 2.As a high-precision electronic device, this instrument should be stored properly. Avoid humid conditions, direct sunlight, and proximity to flammable or explosive materials to prevent damage.
- 3.Always use the designated accessory provided or recommended by the manufacturer to ensure performance throughout its service life.
- 4.Do not place the instrument in high-temperature conditions, as this may negatively impact its service life and performance.
- 5.Do not apply the instrument on eyes, thyroid, parathyroid, testicles, pregnant women's abdomen, pacemakers, etc.
- 6.Patients should consult their doctor before using this device.
- 7.Ensure that the instrument is applied on moist skin and avoid the operation on dry skin.
- 8.Power off the instrument when not in use to ensure its safety.
- 9.If taking anti-adiposity medications, discontinue use 1-2 months prior to using this device, as failure to do so may prolong the treatment course.
- 10.To maintain hygiene and extend the instrument's service life, sterilize it with saline.

## 2. Contraindications


Before using this instrument, it is essential to consult your doctor or seek advice from professionals. Details below:

1. Pregnancy or lactation.
2. Heart conditions, including the presence of a pacemaker.
3. Unhealed wounds or ongoing recovery from recent surgery.
4. Epilepsy, severe diabetes, or hyperthyroidism.
5. Malignant tumors, hemophilia, or bleeding disorders.
6. Skin conditions or infectious diseases.
7. Cesarean section within the last six months, miscarriage within the last three months, or natural birth within the last two months.
8. Recent plastic surgery (within three months), or presence of implants or metal objects in the treatment area.

## 3. Overview

With the most popular and safest technology, this vacuum beauty instrument ensures a radiant transformation for both your body and skin. It provides an ideal solution for those seeking skin and body improvements without injections, drugs, or side effects. This non-invasive device delivers visible enhancements with each use.


## 4. Control Panel Introduction

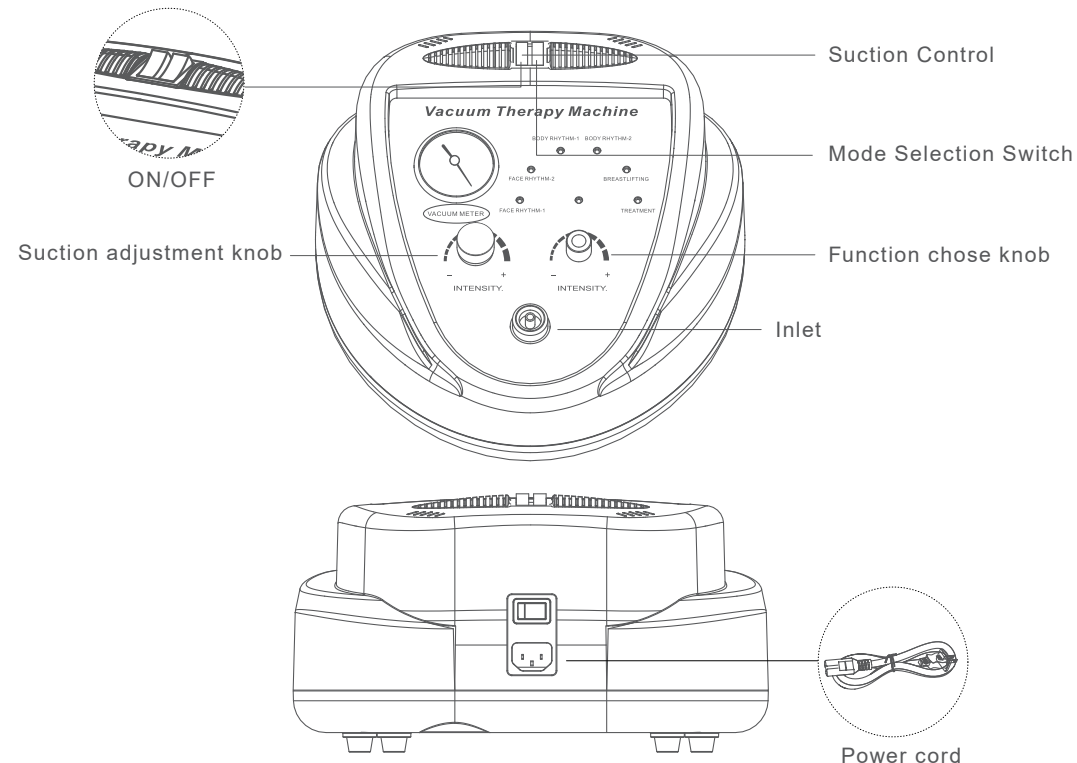


## How to replace the cotton filter




## Hose Connections

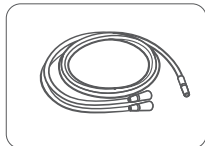

Connect the wide Y-shaped hose for the vacuum cup.

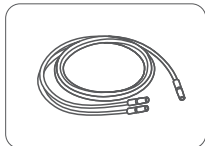

Connect the narrow Y-shaped hose for the breast cup.

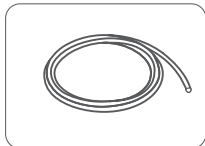

Connect the straight hose to the roller cup.

## Hose Connection Instructions

- **Wide/Narrow Y-Shaped Hose:**
  - Connect one end to the inlet on the main unit.
  - Attach the two remaining ends to the two vacuum or breast cups.
  - For single-cup use: Tie a knot on the spare end. Untie if switching to two-cup use.
- **Straight Hose:**
  - Connect one end to the inlet on the main unit.
  - Attach the other end to the roller cup.


## Mode Settings and Instructions

### 1.Power Activation

- Press the power button on the back of the main

### 2.Controls Overview

- The main unit has two top switches:unit.

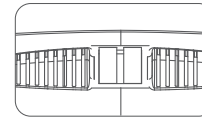

### 3.Mode Selection

- Choose the desired mode with the "Function"
  - Face Rhythm 1: Intermittent suction for facial massage and relaxation.
  - Face Rhythm 2: Relaxation mode for face lifting and tightening.
  - Body Rhythm 1: Detox mode for clearing channels.
  - Body Rhythm 2: Dredge mode to unclog meridians.
  - Breast Lifting: Mode for sculpting and firming the breast and buttocks.
  - Treatment: Massage and relaxation mode for whole-body relaxation, ideal for first-time users.

### 4.Suction Intensity Adjustment

- Use the rotary knob to adjust suction
  - Decrease Intensity: Turn left.
  - Increase Intensity: Turn right.
  - Suction-Only Mode: Turn on only the left switch, adjusting intensity as needed for individual comfort.



#### 4. How to choose the breast cup, vacuum cup, and grease cup?



3 pairs of breast cup



9 pairs of vacuum cup

Cup No. 1: For larger-sized breasts, providing firming and sculpting.

Cup No. 2: For normal-sized breasts, offering firming and enhancement.

Cup No. 3: For smaller-sized breasts, ideal for massage and relaxation.

Cup No. 4: For cupping and moving cupping therapy on large areas (e.g., thighs, buttocks, abdomen).

Cup No. 5: For cupping and moving cupping therapy on medium areas (e.g., legs, back).

Cup No. 6: For cupping and moving cupping therapy on small areas (e.g., arms, lower legs).

Cup No. 7: For moving cupping therapy on areas like the neck and arms.


Cup No. 8: For moving cupping therapy and massage on the face and neck.

Cup No. 9: For moving cupping therapy and cleansing of the face.

Cup No. 10: For moving cupping therapy and general cleansing.

Cup No. 11: For facial cleansing.

Cup No. 12: For sensitive or intimate areas.

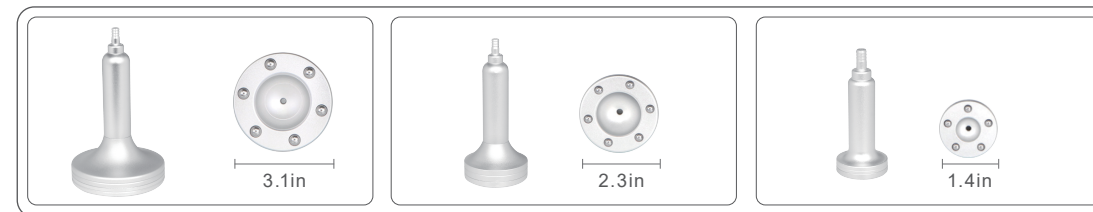

Select the appropriate roller cup from the three available sizes based on the size of the treatment area.



Large Cup: Suitable for extensive body areas, including arms, abdomen, waist, buttocks, legs, and back.



Medium Cup: Ideal for moderate areas, such as arms, neck, face, and lower legs.



Small Cup: Designed for targeted facial areas, including cheeks, around the eyes, forehead, chin, and neck.

## 5. Technical specifications

	Power supply Input: 110v-220v 50Hz/60Hz	
	Power:15W	Input Power:15W
	Dimensions: 12"×11.4"×6"	
	Net Weight: 3.80kg	Gross Weight: 4.90kg
	Power:15W	Air flow: 10LPM
	Vacuum Suction: 55cmHg	Noise level: <65db (30cm away)

## 6. Use Precautions

1. Strictly follow the instructions in the manual to use the instrument and train operator.
2. Before using the instrument, use a plug with a grounding pin and ensure the power socket is grounded.
3. Make sure the voltage of the instrument is adapted. If the local power supply voltage is unstable, we recommend adding a power-matched regulated power supply between the mains and the instrument.

- 4.Remove all the metal objects inside the body to avoid causing accidents or affecting treatment effects.
- 5.Refrain from using the instrument on an empty stomach, and initiate the treatment one hour after a meal.
- 6.Always start with the lowest suction intensity.

## 7.Maintenance & Cleaning

- 1.Disconnect the power cord when not in use and store the instrument in a dry, cool, and well-ventilated place.
- 2.To clean the vacuum cup and roller cup, wipe them with a wet wipe or saline to clean and sterilize, then dry them with a dry towel.
- 3.After cleaning the instrument, cover it with a dustcloth or plastic wrap. Before use, ensure to sanitize the handle again for hygienic reasons.
- 4.Disassembling the instrument and attempting to clean the inside is strictly prohibited. This task should only be performed by professionals.

## 8.FAQ

- 1.Are there any side effects associated with vacuum breast enlargement treatment?

A: No, there are no side effects at all. Vacuum breast enlargement therapy is entirely based on physical massage, suction, and release under negative pressure. With the integration of multiple advanced technologies, it provides users with immediate effects without any side effects or rebound effects.

- 2.When do I expect to see the result of facial treatment with the vacuum cup?

A: After treatment, you'll notice purified pores and smoother skin instantly. We recommend a weekly treatment for best results, with increasing improvements day by day.

- 3.How does negative pressure aid in detoxification during treatment?

A: Negative pressure can stimulate capillary congestion and boost cellular vitality. The alternation of negative pressure promotes the opening and closing of pores in the targeted area, improving skin breathing, increasing oxygen uptake, and aiding in the elimination of wastes. In summary, this treatment promotes Qi and circulation, dispels dampness and cold, clears meridians, detoxifies, and purges heat.

## 9. Ask & Answer

Ask	Cause	Answer
Fail to turn on the instrument.	→ The power cord is disconnected from the outlet.	→ Check if the power cord is securely connected to the outlet and ensure that the outlet is functional.
	→ The fuse is loose or burned.	→ Replace the fuse on the back of the instrument and restart it.
Weak or no suction from the vacuum cup/massage roller.	→ The hose is not securely attached.	→ Ensure the hose is tightly secured to the main unit.
	→ The cap on the inlet is loose or the rubber ring inside is broken or missed.	→ Secure the cap on the inlet tightly or replace the rubber ring inside.
	→ Check if you have set the knob.	→ Set the knob.

## 10. Operation guide


Breast enlargement:




Butt augmentation:




Face:

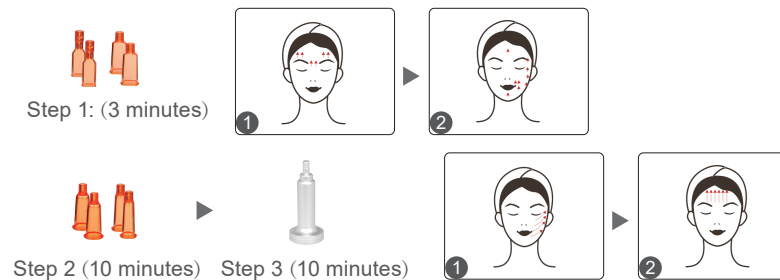


Arm:

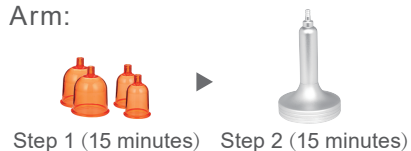




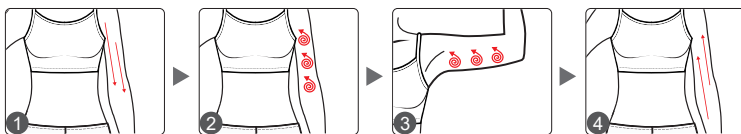


Back:

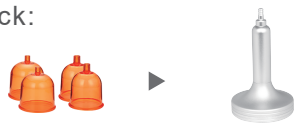


Step 1 (15 minutes) Step 2 (15 minutes)

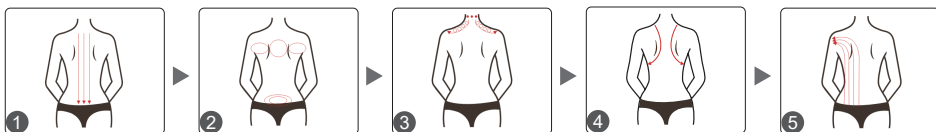

Leg:

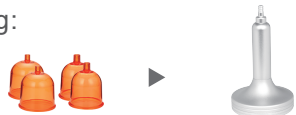

Step 1 (15 minutes) Step 2 (15 minutes)

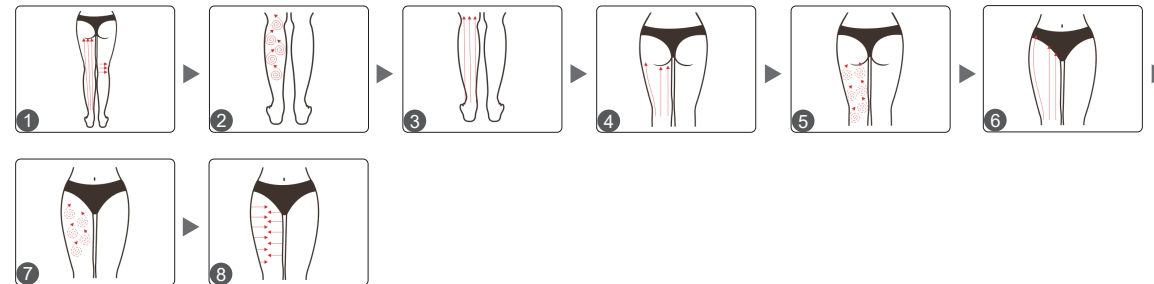

## 11. Warranty and Support

We offer a one-year warranty on the main unit and a three-month warranty on components.

Product Name: Vacuum Therapy Instrument

Product Model: MS-2184

Production Date: \_\_\_\_\_

Quality Inspection: \_\_\_\_\_

PASS